


Evolution® Duodenal

CONTROLLED-RELEASE STENT - UNCOVERED

Features




Step 1

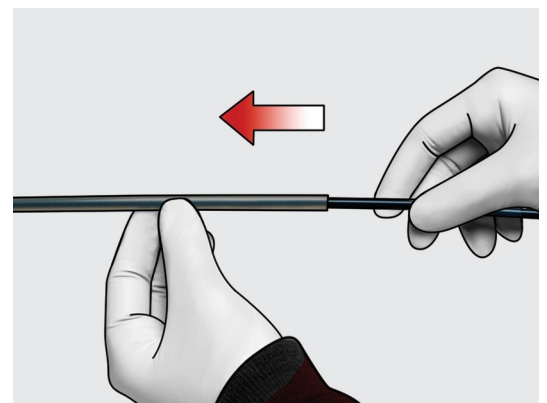


Remove the protective tubing from the stent.

Step 2

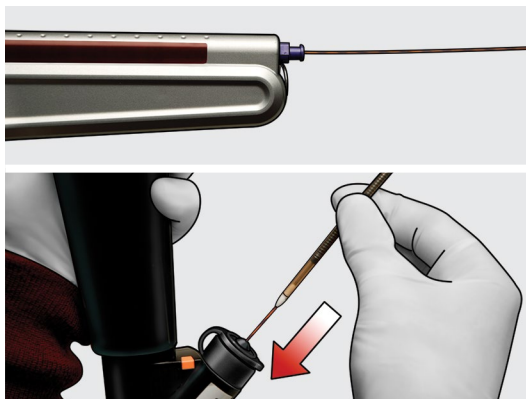

Introduce the delivery system in short increments over the wire guide, into the accessory channel, until it is endoscopically visualised exiting the scope. **NOTE:** Minimum accessory channel 3.7 mm

Step 3a

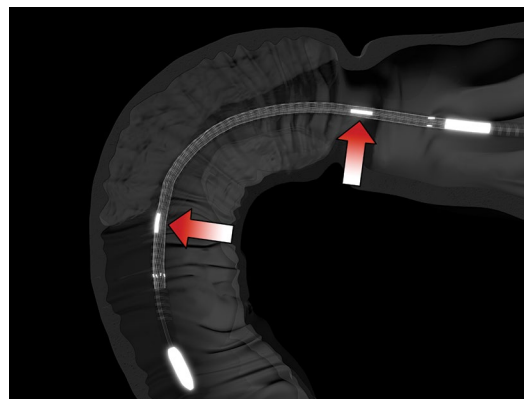


Fluoroscopically position the radiopaque markers on the inner catheter beyond the extremities of the stricture to be crossed.

Step 3b



NOTE: The yellow marker on the delivery system located at the proximal end of the stent can also serve as an endoscopic/fluoroscopic reference for the positioning of the proximal (nearest to the user) end of the stent relative to the stricture.

Step 4

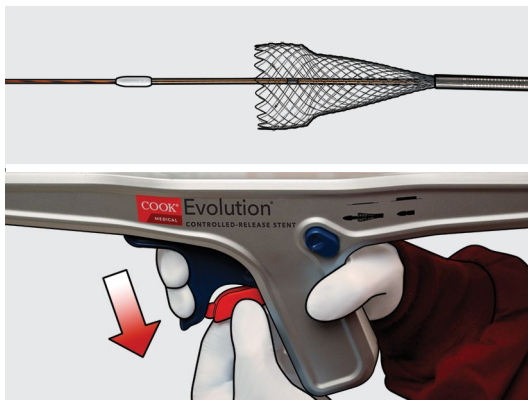

Confirm the desired stent position fluoroscopically and deploy the stent by removing the red safety guard from the handle. Continue deploying the stent by squeezing the trigger.


Step 5a



If repositioning is required during deployment, it is possible to recapture the stent.

NOTE: It is not possible to recapture the stent after passing the point-of-no-return, indicated when the stent deployment indicator on top of the introducer has passed the point-of-no-return indicator on the handle.

Step 5b

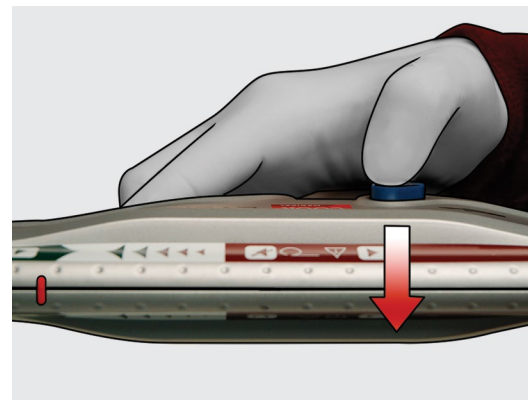

To recapture the stent, push the directional button on the side of the delivery system to opposite side.


NOTE: Hold your thumb on the button when squeezing the trigger for the first time to recapture. Continue squeezing the trigger as required to recapture the stent by the desired amount.

Step 6




To resume deployment, push the directional button to the opposite side again and hold the button for the first stroke while squeezing the trigger.

Step 7

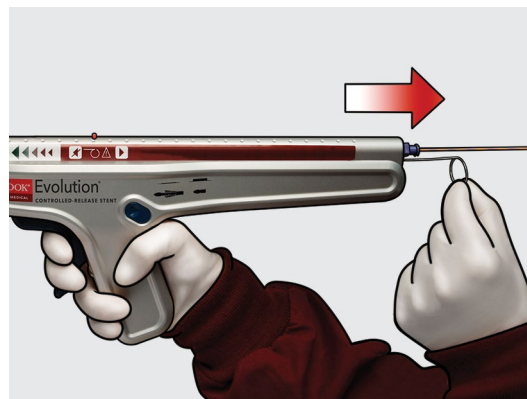

When the point-of-no-return has been passed, pull the safety wire out of the delivery handle near the wire guide port. Continue deploying the stent by squeezing the trigger.

Step 8

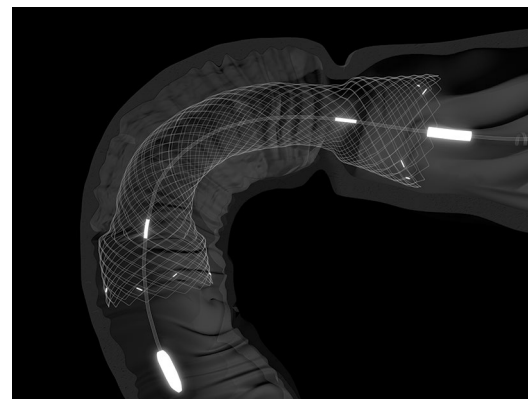

After deployment, fluoroscopically confirm the stent expansion. Once the expansion is confirmed, the introduction system can be safely removed.

Refer to current instructions for detailed system use.

Please see product risk information in the IFU at cookmedical.eu

Customer Service

EU Website: cookmedical.eu
EDI: cookmedical.eu/edi
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001702, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Iceland: +354 8007615, IS.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Norway: +47 23162968, no.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com
EDI: cookmedical.com/edi.do
Americas:
Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
E-mail: customersupport@cookmedical.com
Australia:
Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
E-mail: cau.custserv@cookmedical.com



AI, ESC, IR, OHNS, PI, RH, SUR-8.5X11