


ConquestTTC[®]

LITHOTRIPTOR CABLE

Important

Prior to mechanical lithotripsy, refer to current instructions for detailed system use. Also refer to Cook product labels to confirm compatibility of baskets and minimum accessory channel. **Wire cutters must be available for this procedure.**

Step 1

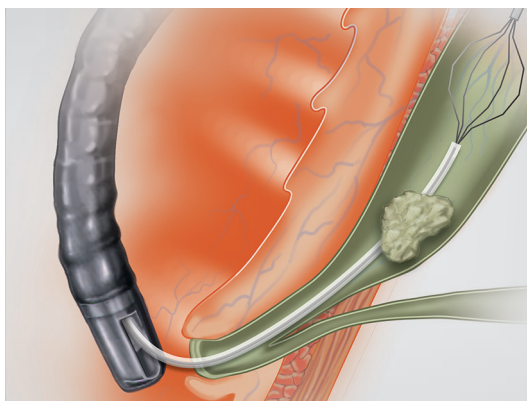


Fluoroscopically confirm stone entrapment with basket.

Step 2




Once the stone is completely captured in the basket, make a complete cut using wire cutters, at the junction of the 10 F and 7 F sheaths.

Step 3

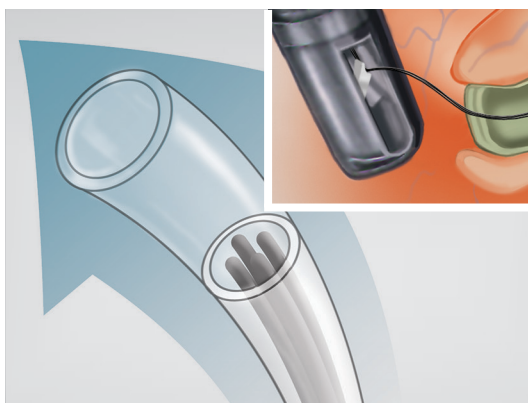


Remove the outer sheath of the basket slowly. As the distal tip of the sheath clears the endoscope elevator, gently lift the elevator to lock the basket wires in place. Then, continue removing the outer sheath. Dispose of the outer sheath.

Step 4

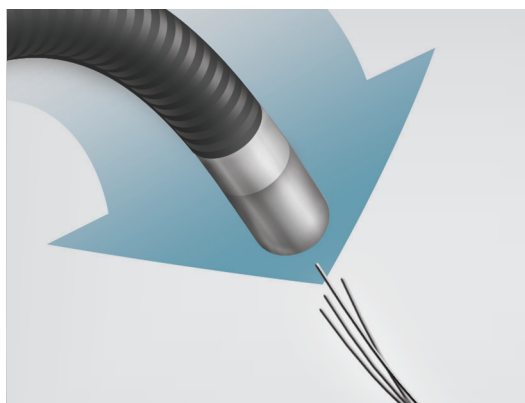


Introduce the lithotripsy cable over the basket wires and through the channel until the wires exit the proximal end of the cable. Grasp the basket wires exiting the cable.

Step 5

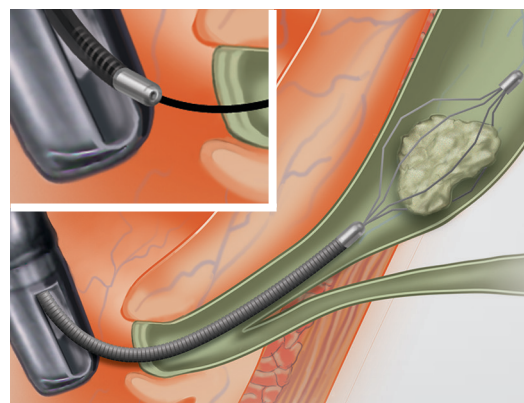



Release the elevator and continue to advance the cable until the distal tip reaches the basket within the duct.

Step 6

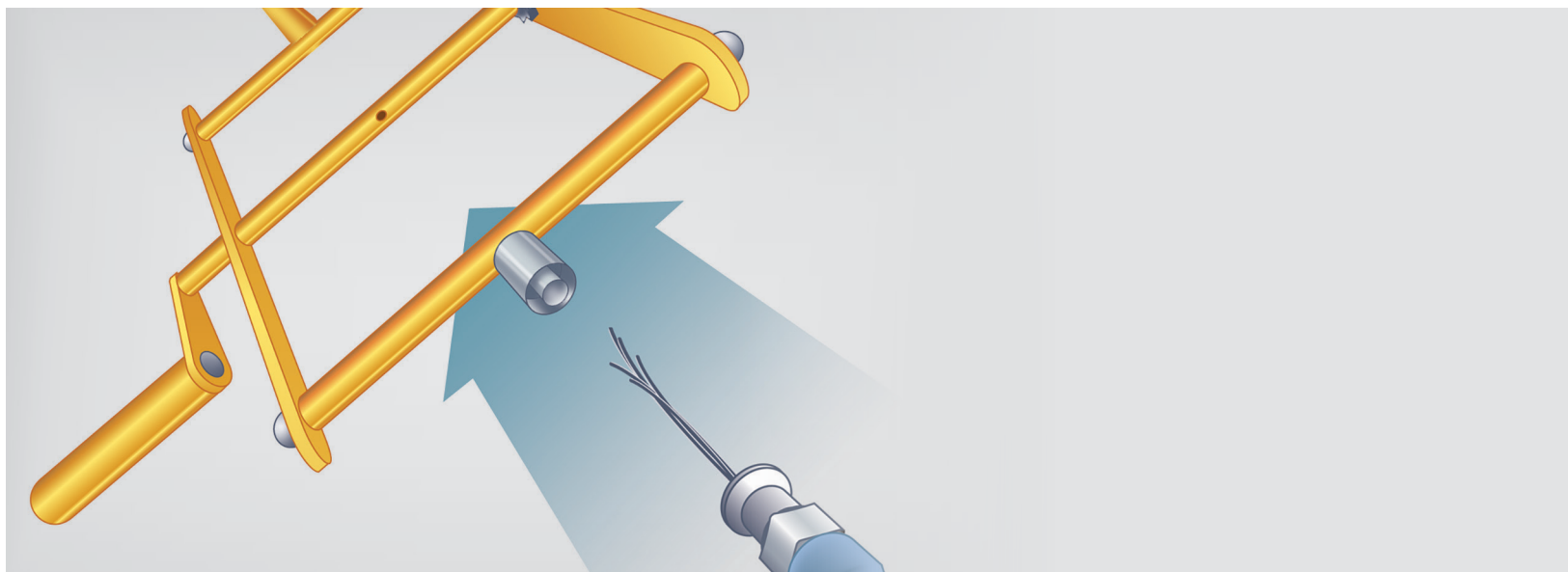

Pass the basket wires through the male Luer lock on the lithotripsy handle and through the rotating rod opening of the handle.

Step 7

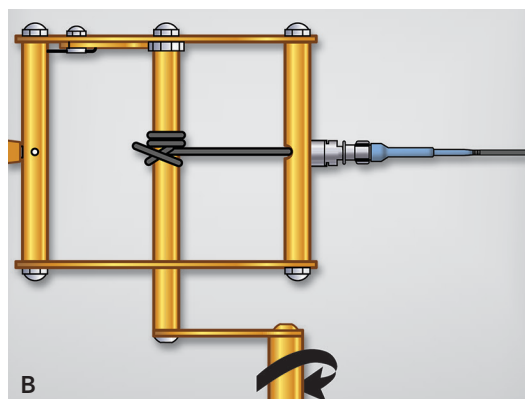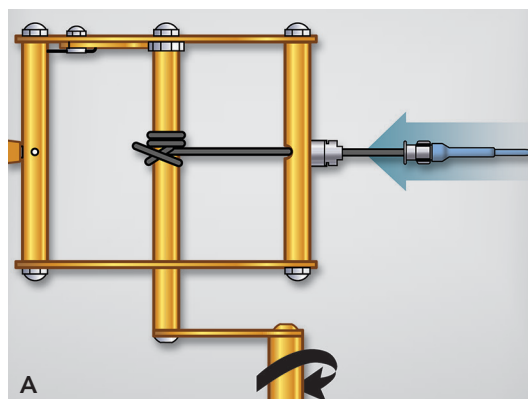


A. Connect the male Luer lock of the lithotripter handle to the female Luer lock of the cable. Wrap the basket wires around the rotating rod to secure.

B. Apply steady pressure by rotating the handle slowly. Sudden release of pressure may indicate fracture of the stone, fracture of the basket, or loss of the stone from the basket.

Warning: If basket fragmentation should occur, surgical intervention may be necessary.

Note: The handle will only move in one direction unless the ratchet is released.

Step 8

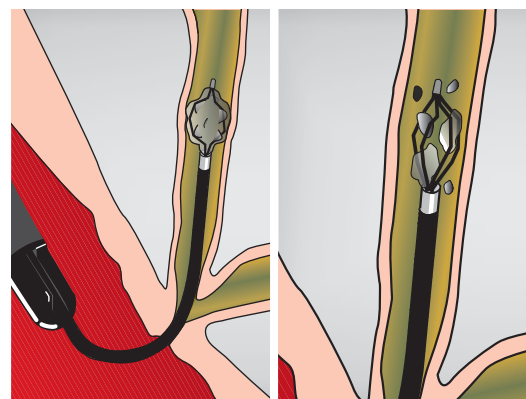

After the stone has been fractured, assess the need to remove the fractured stone from the duct. Release the ratchet on the handle and remove the basket.

Refer to current instructions for detailed system use.

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001702, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Iceland: +354 8007615, is.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com
 EDI: cookmedical.com/edi.do
Americas:
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: customersupport@cookmedical.com
Australia:
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 E-mail: cau.custserv@cookmedical.com



AI, ESC, IR, OHNS, PI, RH, SUR-8.5X11