


# Rely on CXI

Your advanced support catheter

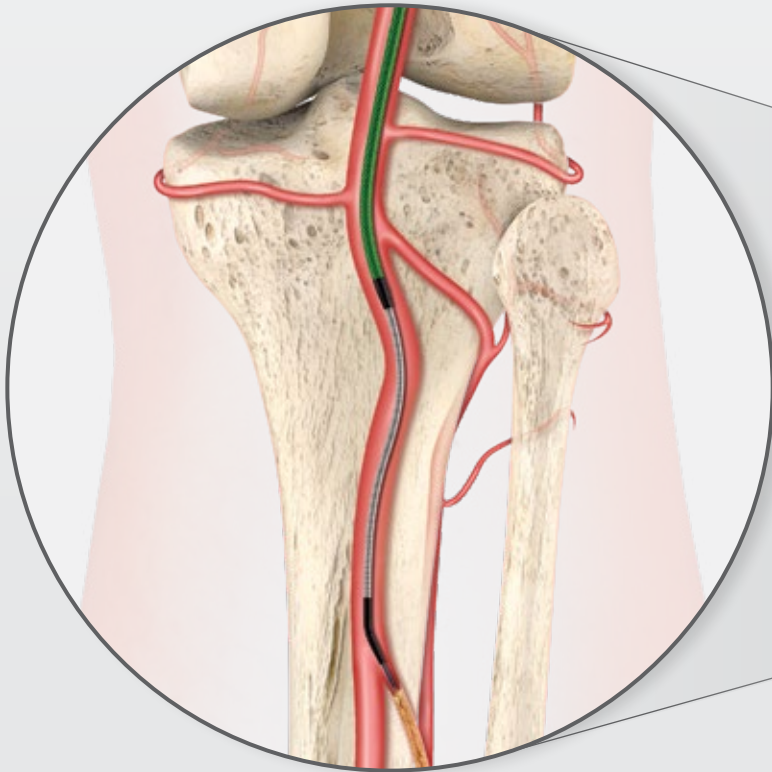

**ENHANCED  
VISIBILITY**  
.018"  
platform

**CXI**<sup>®</sup>  
SUPPORT CATHETER

**Pushability**  
CXI's stainless steel braid provides responsive torque and enhances pushability.

**Coaxial capability**  
The 2.3 Fr and 2.6 Fr sizes can fit coaxially in a 4.0 Fr size to increase the amount of support.

**Compatibility**  
CXI is compatible with our limb preservation tools: micropuncture kits, sheaths, wires, and balloons for antegrade and retrograde interventions.




**.018 Platform**  
Enhanced visibility plus additional lengths and tip configurations

An angled tip allows you to select branched vessels.

The seamless transition from the tip to the wire guide is designed to minimize vessel trauma.




Platinum-iridium marker bands are at the tip and at 5, 10, and 15 cm to enhance visibility and to allow vessel measurement.

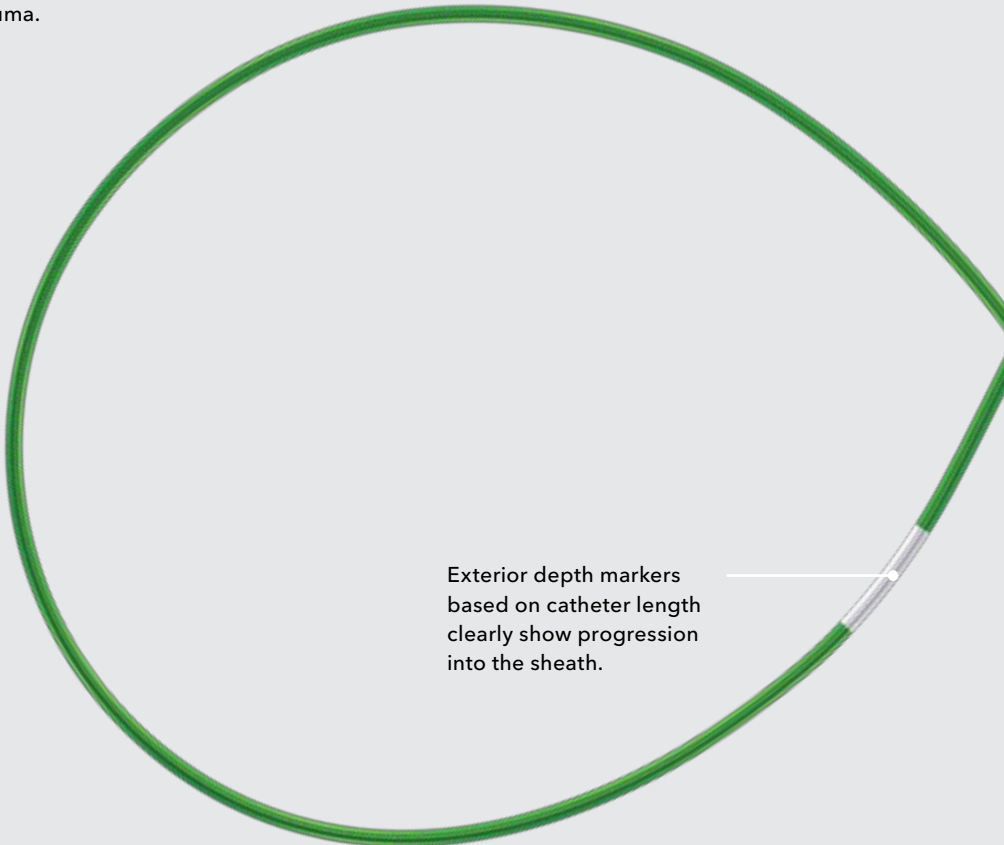

Exterior depth markers based on catheter length clearly show progression into the sheath.

Stainless steel braid provides responsive torque and enhances pushability.

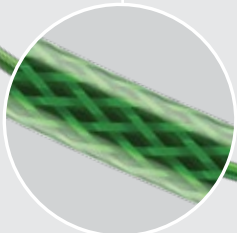
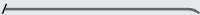



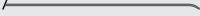






The hydrophilic coating on the distal 40 cm increases pushability and trackability.







**SUPPORT CATHETER**

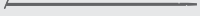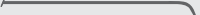Order Numbers	Reference Part Number	Fr	Accepts Wire Guide Diameter inch	Length cm	Tip Configuration
<b>With platinum-iridium marker bands</b>					
2.3 Fr					
G18369	CXI-2.3-14-65-0	2.3	.014	65	Straight 
G18371	CXI-2.3-14-90-0	2.3	.014	90	
G23160	CXI-2.3-14-135-0	2.3	.014	135	
G18373	CXI-2.3-14-150-0	2.3	.014	150	ANG 
G18370	CXI-2.3-14-65-ANG	2.3	.014	65	
G18372	CXI-2.3-14-90-ANG	2.3	.014	90	
G23162	CXI-2.3-14-135-ANG	2.3	.014	135	ANG2 
G18376	CXI-2.3-14-150-ANG	2.3	.014	150	
G23858	CXI-2.3-14-65-ANG2	2.3	.014	65	
G23859	CXI-2.3-14-90-ANG2	2.3	.014	90	ANG2 
G23860	CXI-2.3-14-135-ANG2	2.3	.014	135	
G23861	CXI-2.3-14-150-ANG2	2.3	.014	150	
<b>Now available with platinum-iridium marker bands</b>					
2.6 Fr					
G18377	CXI-2.6-18-65-0	2.6	.018	65	Straight 
G56550	CXI-2.6-18-90-0	2.6	.018	90	
G23161	CXI-2.6-18-135-0	2.6	.018	135	
G56552	CXI-2.6-18-150-0	2.6	.018	150	ANG 
G18378	CXI-2.6-18-65-ANG	2.6	.018	65	
G56551	CXI-2.6-18-90-ANG	2.6	.018	90	
G23163	CXI-2.6-18-135-ANG	2.6	.018	135	ANG2 
G56553	CXI-2.6-18-150-ANG	2.6	.018	150	
G23862	CXI-2.6-18-65-ANG2	2.6	.018	65	
G23863	CXI-2.6-18-90-ANG2	2.6	.018	90	ANG2 
G23864	CXI-2.6-18-135-ANG2	2.6	.018	135	
G23865	CXI-2.6-18-150-ANG2	2.6	.018	150	
<b>With tungsten marker bands</b>					
4.0 Fr					
G52546	CXI-4.0-35-90-P-NS-0	4.0	.035	90	Straight 
G52542	CXI-4.0-35-135-P-NS-0	4.0	.035	135	
G52544	CXI-4.0-35-150-P-NS-0	4.0	.035	150	
G52547	CXI-4.0-35-90-P-NS-DAV	4.0	.035	90	ANG 
G52543	CXI-4.0-35-135-P-NS-DAV	4.0	.035	135	
G52545	CXI-4.0-35-150-P-NS-DAV	4.0	.035	150	

Tip options

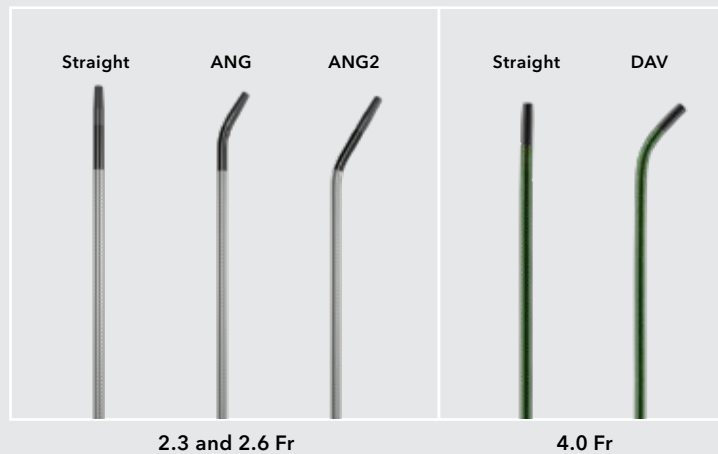

CXI distal tip under fluoroscopy

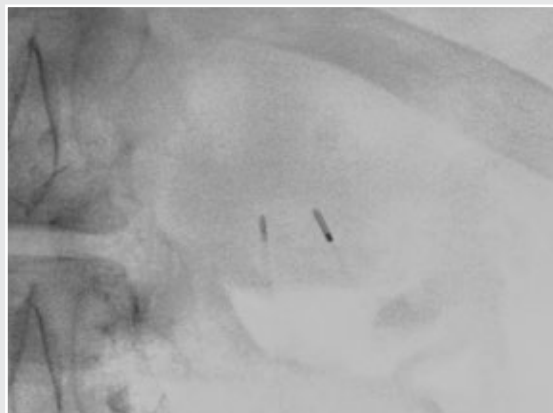

The image on the left shows the CXI with tungsten marker bands, and the image on the right shows the CXI with platinum-iridium marker bands under fluoroscopy.

Compatible products

	Antegrade Approach (Femoral Access)	Retrograde Approach (Pedal Access. Compatible with 2.9 Fr sheath)
Flexor® Ansel Guiding Sheath with High-Flex Dilator	X	
Flexor Shuttle Tibial Infrapopliteal Access System	X	
Micropuncture Pedal (Introducer) Access Set		X
Advance® 35LP Low-Profile Balloon Catheter	X	
Advance 18LP Low-Profile Balloon Catheter	X	
Advance 14LP Low-Profile Balloon Catheter	X	
Advance Micro® 14 Ultra Low-Profile Balloon Catheter		X
Roadrunner® Uniglide Hydrophilic Wire Guide	X	
Approach® CTO Microwire Guide	X	X
Approach Hydro ST Microwire Guide	X	X

Customer Service

EU Website: [cookmedical.eu](http://cookmedical.eu)  
 EDI: [cookmedical.eu/edi](http://cookmedical.eu/edi)  
 Distributors: +353 61239240, [ssc.distributors@cookmedical.com](mailto:ssc.distributors@cookmedical.com)  
 Austria: +43 179567121, [oe.orders@cookmedical.com](mailto:oe.orders@cookmedical.com)  
 Belgium: +32 27001702, [be.orders@cookmedical.com](mailto:be.orders@cookmedical.com)  
 Denmark: +45 38487607, [da.orders@cookmedical.com](mailto:da.orders@cookmedical.com)  
 Finland: +358 972519996, [fi.orders@cookmedical.com](mailto:fi.orders@cookmedical.com)  
 France: +33 171230269, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
 Germany: +49 6950072804, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
 Hungary: +36 17779199, [hu.orders@cookmedical.com](mailto:hu.orders@cookmedical.com)  
 Iceland: +354 8007615, [IS.orders@cookmedical.com](mailto:IS.orders@cookmedical.com)  
 Ireland: +353 61239252, [ie.orders@cookmedical.com](mailto:ie.orders@cookmedical.com)  
 Italy: +39 0269682853, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
 Netherlands: +31 202013367, [nl.orders@cookmedical.com](mailto:nl.orders@cookmedical.com)  
 Norway: +47 23162968, [no.orders@cookmedical.com](mailto:no.orders@cookmedical.com)  
 Spain: +34 912702691, [es.orders@cookmedical.com](mailto:es.orders@cookmedical.com)  
 Sweden: +46 858769468, [se.orders@cookmedical.com](mailto:se.orders@cookmedical.com)  
 Switzerland - French: +41 448009609, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
 Switzerland - Italian: +41 448009609, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
 Switzerland - German: +41 448009609, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
 United Kingdom: +44 2073654183, [uk.orders@cookmedical.com](mailto:uk.orders@cookmedical.com)

USA Website: [cookmedical.com](http://cookmedical.com)  
 EDI: [cookmedical.com/edi.do](http://cookmedical.com/edi.do)  
**Americas:**  
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335  
 E-mail: [orders@cookmedical.com](mailto:orders@cookmedical.com)  
**Australia:**  
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283  
 E-mail: [cau.custserv@cookmedical.com](mailto:cau.custserv@cookmedical.com)



AI, ESC, IR, OHNS, PI, RH, SUR-8.5X11