


Block with simplicity.
Direct with ease.




ENDOBONCHIAL BLOCKER SETS

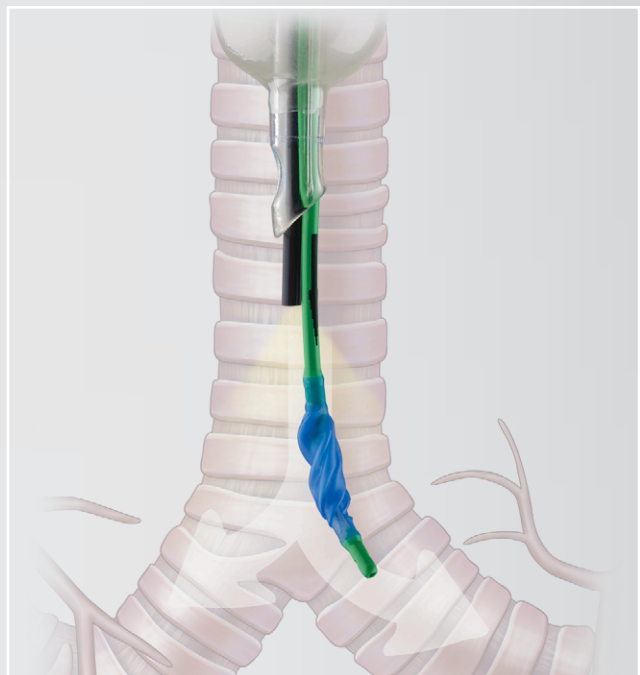

ARNDT ENDOBRONCHIAL BLOCKER

Cook's endobronchial blockers enable one-lung ventilation through a conventional single-lumen endotracheal tube (SLT). The SLT is easier to place than a double-lumen tube (DLT) and poses less risk of trauma to internal structures.¹ The loss of an airway, a potential problem during the exchange of a DLT for an SLT, is also eliminated with this method.²

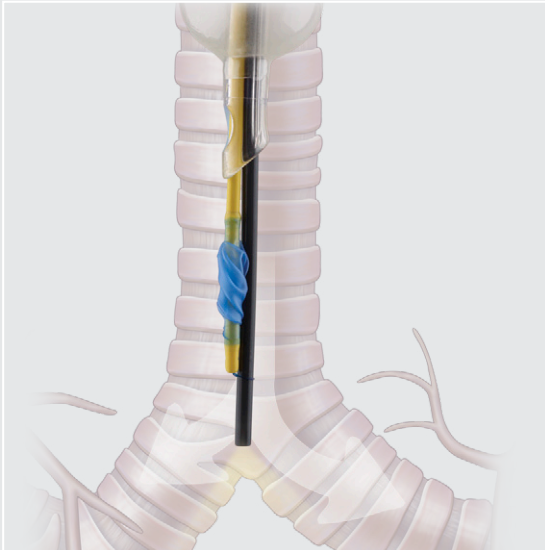

Features and benefits

- The adjustable guide loop is coupled with an appropriately sized bronchoscope to enable precise placement.
- The low-pressure, high-volume balloon creates excellent surface area contact with the inner bronchial wall while minimizing potential trauma to the bronchus.
- The distance markings indicate how far the blocker has been advanced.
- The suction adapter included with the set offers an additional method of lung deflation.


ETT and Fiberoptic Bronchoscopy (FOB) Recommendations

Order Number	Reference Part Number	Blocker Fr	ETT Diameter mm	FOB OD mm
G44109	C-AEBS-5.0-50-SPH-AS	5.0	≥ 4.5	≤ 2.8
G44114	C-AEBS-7.0-65-SPH-AS	7.0	≥ 6.0	≤ 3.6
G44120	C-AEBS-9.0-78-SPH-AS	9.0	≥ 7.5	≤ 4.2

Average Inflation Volumes

Order Number	Reference Part Number	Blocker Fr	Balloon Volume ml	Balloon Configuration
G44109	C-AEBS-5.0-50-SPH-AS	5.0	0.5-2	spherical
G44114	C-AEBS-7.0-65-SPH-AS	7.0	2-6	spherical
G44120	C-AEBS-9.0-78-SPH-AS	9.0	4-8	spherical

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

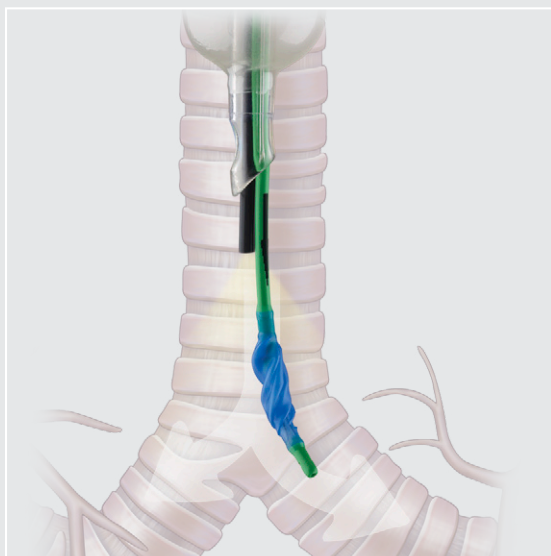

Features and benefits

- The unique deflecting tip provides a wide range of motion to precisely direct the blocker under bronchoscopic guidance.
- The high-traction silicone grip on the introducer shaft increases the operator’s control, resulting in precise device placement.
- The low-pressure, high-volume balloon creates excellent surface area contact with the inner bronchial wall while minimizing potential trauma to the bronchus.
- Distance markings increase visualization by indicating how far the blocker has been advanced.
- The suction adapter included with the set offers an additional method of lung deflation.

ETT and Fiberoptic Bronchoscopy (FOB) Recommendations

Order Number	Reference Part Number	Blocker Fr	ETT Diameter mm	FOB OD mm
G44122	C-AEBTS-9.0-65-SPH-AS	9.0	≥ 7.0	≤ 4.2

Average Inflation Volume

Order Number	Reference Part Number	Blocker Fr	Balloon Volume mL	Balloon Configuration
G44122	C-AEBTS-9.0-65-SPH-AS	9.0	6-9	spherical

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

References

1. Knoll H, Ziegeler S, Schreiber JU, et al. Airway injuries after one-lung ventilation: a comparison between double-lumen tube and endobronchial blocker: a randomized, prospective, controlled trial. *Anesthesiology*. 2006;105(3):471-477.
2. Arndt GA, DeLessio ST, Kranner PW, et al. One-lung ventilation when intubation is difficult—presentation of a new endobronchial blocker. *Acta Anaesthesiol Scand*. 1999;43(3):356-358.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001633, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Norway: +47 23162968, no.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
E-mail: orders@cookmedical.com

Australia:
Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
E-mail: cau.custserv@cookmedical.com

