



Secure safe passage.

A dynamic duo for protected ureteral access

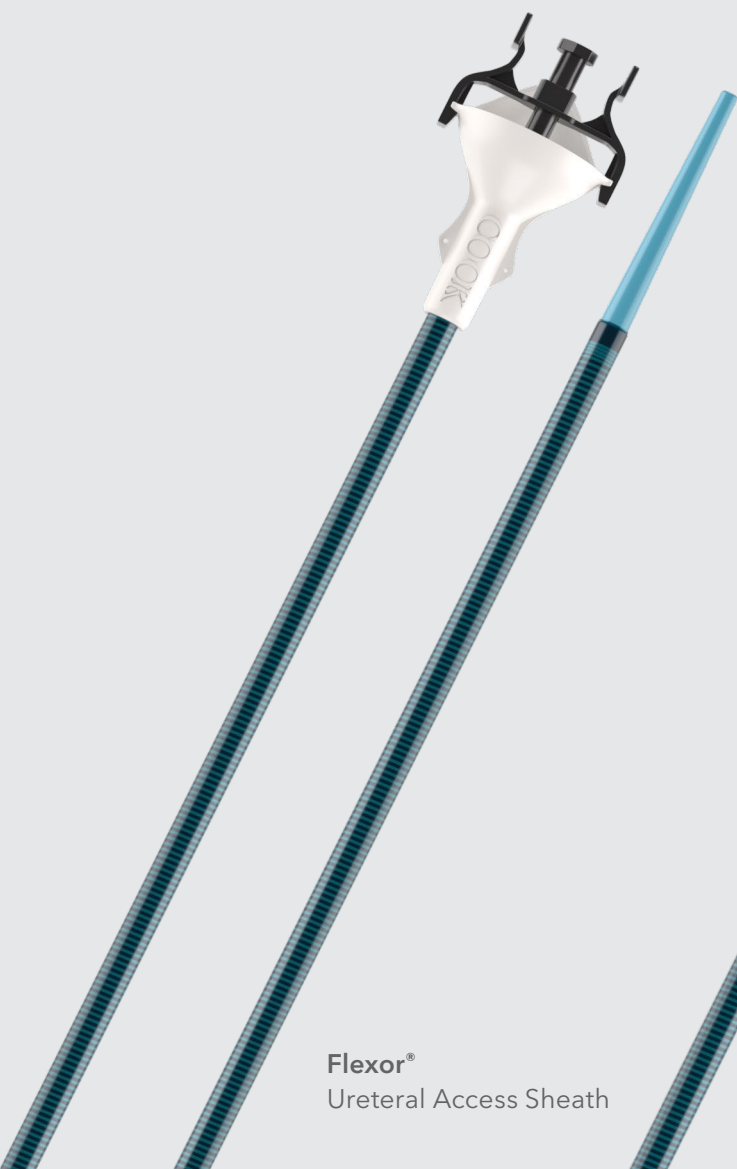


Flexor®
Ureteral Access Sheath

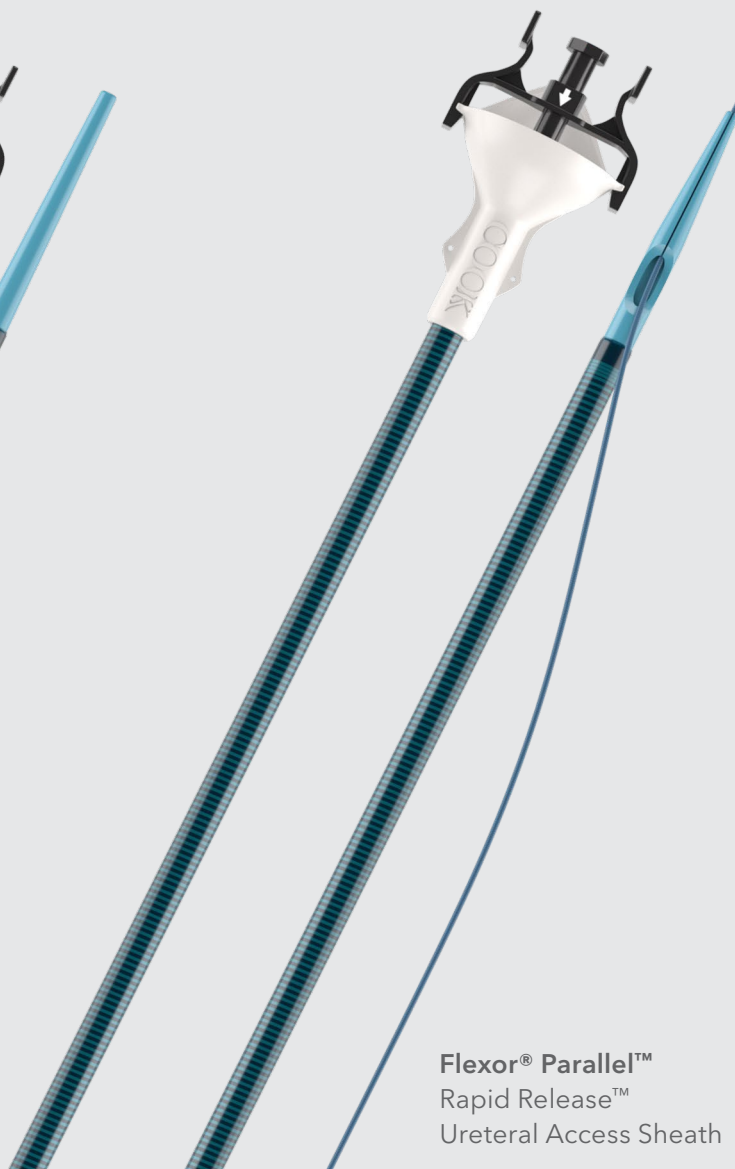

Flexor® Parallel™
Rapid Release™
Ureteral Access Sheath

We were one of the first to make a ureteral access sheath with a reinforced coil core—to resist kinking and compression—and with a hydrophilic coating for ease of placement.

Flexor

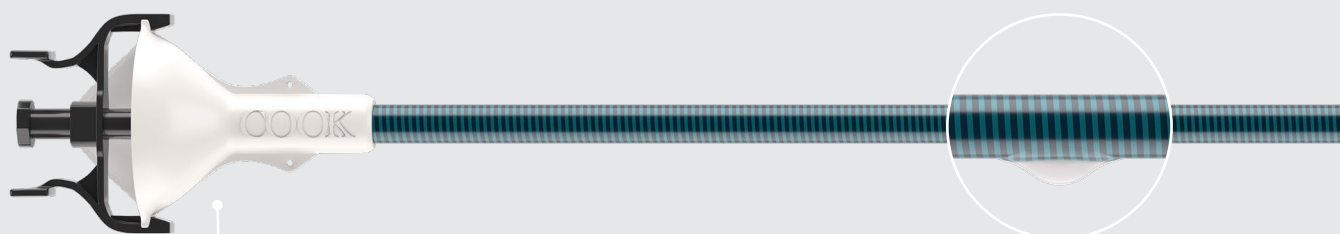

Dilator hub


The sheath's locking mechanism on the dilator hub secures the dilator to the sheath for simultaneous advancement of the sheath and dilator.

Funnel

The ergonomic funnel acts as a handle during insertion. The large trough facilitates instrument introduction.

Hydrophilic coating

The external surface of the sheath and the dilator tip are AQ® hydrophilic coated to create a low-friction surface to ease the insertion process.

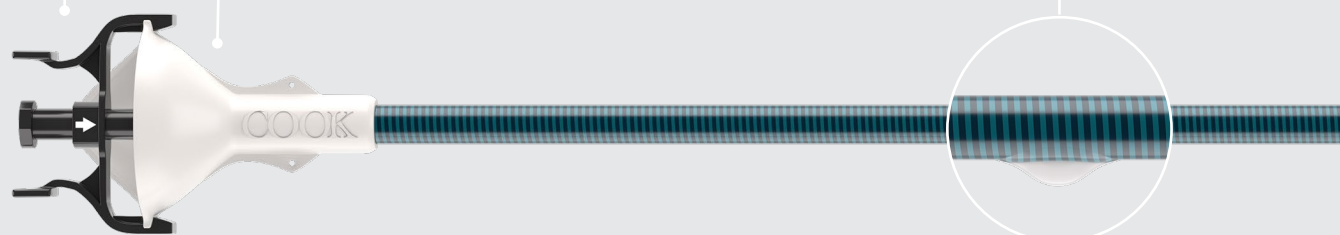

Flexor Parallel


The original

The Flexor Ureteral Access Sheath is used to establish a conduit during endoscopic urological procedures, facilitating the passage of endoscopes and other instruments into the urinary tract. The sheath's working channel protects the ureter during repeated instrument exchanges and protects delicate instruments and smaller flexible scopes from damage. When compared to several competitor sheaths, the Flexor sheath was shown to be less likely to buckle while advancing up the ureter.¹

The upgrade

We doubled down on protected ureteral access when we introduced the Flexor Parallel Rapid Release Ureteral Access Sheath. The Flexor Parallel sheath gives you two options for placement—the traditional Flexor placement or the rapid-release technique. The latter technique is designed to save time in the operating theatre and save product by eliminating the need for a second wire guide.

1. Pedro RN, Hendlin K, Durfee WK, et al. Physical characteristics of next-generation ureteral access sheaths: buckling and kinking. *Urology*. 2007;70(3):440-442.

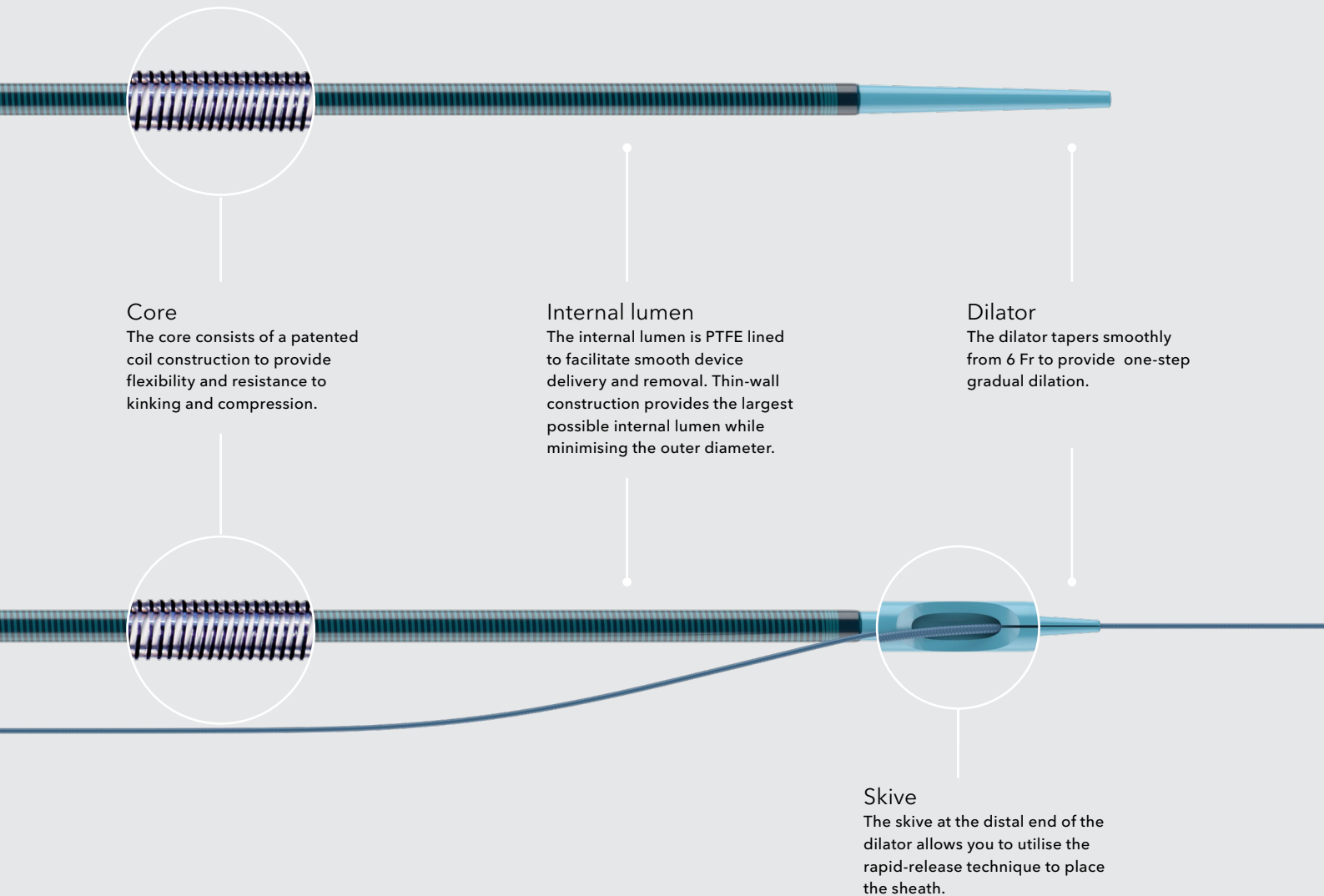

Flexor®

URETERAL ACCESS SHEATH

Order Number	Reference Part Number	ID Fr	Length cm
G19177	FUS-095020	9.5	20
G19176	FUS-095028	9.5	28
G19174	FUS-095035	9.5	35
G46167	FUS-095045	9.5	45
G19175	FUS-095055	9.5	55
G34877	FUS-107020	10.7	20
G34878	FUS-107028	10.7	28
G34879	FUS-107035	10.7	35
G34880	FUS-107045	10.7	45
G34881	FUS-107055	10.7	55
G19173	FUS-120020	12.0	20
G19172	FUS-120028	12.0	28

Order Number	Reference Part Number	ID Fr	Length cm
G19168	FUS-120035	12.0	35
G46168	FUS-120045	12.0	45
G19169	FUS-120055	12.0	55
G30391	FUS-140020	14.0	20
G30393	FUS-140028	14.0	28
G30389	FUS-140035	14.0	35
G46169	FUS-140045	14.0	45
G30390	FUS-140055	14.0	55
Paediatric			
G32853	FUS-095013	9.5	13
G34876	FUS-107013	10.7	13
G32852	FUS-120013	12.0	13

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Flexor® Parallel™

RAPID RELEASE™ URETERAL ACCESS SHEATH

Order Number	Reference Part Number	ID Fr	Length cm
G56428	FUS-095020-P	9.5	20
G56429	FUS-095028-P	9.5	28
G56430	FUS-095035-P	9.5	35
G56431	FUS-095045-P	9.5	45
G56432	FUS-095055-P	9.5	55
G34871	FUS-107020-P	10.7	20
G34872	FUS-107028-P	10.7	28
G34873	FUS-107035-P	10.7	35
G34874	FUS-107045-P	10.7	45
G34875	FUS-107055-P	10.7	55

Order Number	Reference Part Number	ID Fr	Length cm
G56433	FUS-120020-P	12.0	20
G56434	FUS-120028-P	12.0	28
G56435	FUS-120035-P	12.0	35
G56436	FUS-120045-P	12.0	45
G56438	FUS-120055-P	12.0	55
G56439	FUS-140020-P	14.0	20
G56440	FUS-140028-P	14.0	28
G56441	FUS-140035-P	14.0	35
G56442	FUS-140045-P	14.0	45
G56443	FUS-140055-P	14.0	55

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001702, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Poland: +48 223060159, pl.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com
 EDI: cookmedical.com/edi.do
Americas:
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: customersupport@cookmedical.com
Australia:
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 E-mail: cau.custserv@cookmedical.com



CC, URO-A4