


ERCP is all about smooth, easy access.



D.A.S.H.[™]
SPHINCTEROTOME WITH DOMETIP[®]



Image courtesy of Dr. Douglas Howell,
Maine Medical Center, Portland, ME.

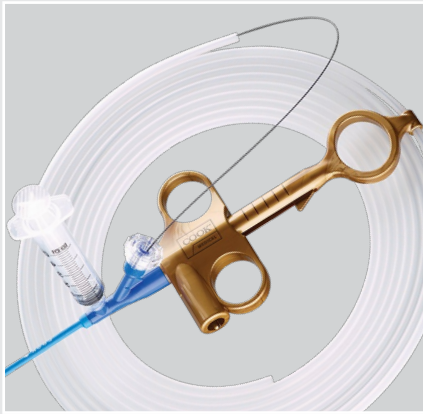

D.A.S.H.[™]

SPHINCTEROTOME WITH DOMETIP[®]

Used for cannulation of the ductal system and for sphincterotomy.

- DomeTip designed with anatomy in mind for access
- Custom 3D forming wire enables reliable cutting wire orientation
- Integrated hub provides triple-lumen functionality on a double-lumen platform, giving you the ability to inject over the wire guide
- Pre-attached racetrack for procedural convenience
- Preloaded options for convenience with a 0.021, 0.025, or 0.035 inch diameter wire guide
- Gold, silver, and bronze depth reference markings endoscopically identify relative position of the 25 mm long monofilament cutting wire within the duct

Order Number	Reference Part Number	Catheter Fr	Tip Shape	Wire Guide Diameter inch	Wire Guide Length cm	Minimum Accessory Channel Diameter mm
G26850	DASH-21**	5.4	DomeTip	0.021	–	2.8
G25224	DASH-1*	6.3	DomeTip	0.025	–	2.8
G54193	DASH-35*	6.6	DomeTip	0.035	–	2.8

Preloaded with Acrobat[®] 2 Calibrated Tip Wire Guide

G57298	DASH-AWG2-25-450	6.3	DomeTip	0.025	450	2.8
G57299	DASH-AWG2-35-260	6.6	DomeTip	0.035	260	2.8
G57300	DASH-AWG2-35-450	6.6	DomeTip	0.035	450	2.8

Preloaded with Tracer Metro[®] Direct Wire Guide

G25790	DASH-21-480	5.4	DomeTip	0.021	480	2.8
G25225	DASH-260	6.3	DomeTip	0.025	260	2.8
G25227	DASH-480	6.3	DomeTip	0.025	480	2.8
G25666	DASH-35-480	6.6	DomeTip	0.035	480	2.8

Minimum accessory channel diameter 2.8 mm. *Recommend Acrobat[®] 2 Calibrated Tip Wire Guide, sold separately. **Recommend Tracer Metro[®] Direct Wire Guide, sold separately. Active cord available separately.

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Delivery for details.

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001633, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: orders@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283

E-mail: cau.custserv@cookmedical.com-A4



AI, ESC, IR, OHNS, PI, RH, SUR-A4