
Get glowing results.



GIO-TIPII® double lumen ercp catheter

The Glo-Tip II catheter balances two significant factors in cannulating procedures – access and visualization. The catheter's DomeTip® permits smooth, easy cannulation and a radiopaque band inside the wire lumen emits a distinctive "glow" under fluoroscopy. And, for convenience, its double lumen design permits contrast injection through a dedicated lumen while the catheter is placed over a wire guide. The Glo-Tip II: part of Cook Medical's full line of innovative endoscopic accessories designed to meet your patients' needs.



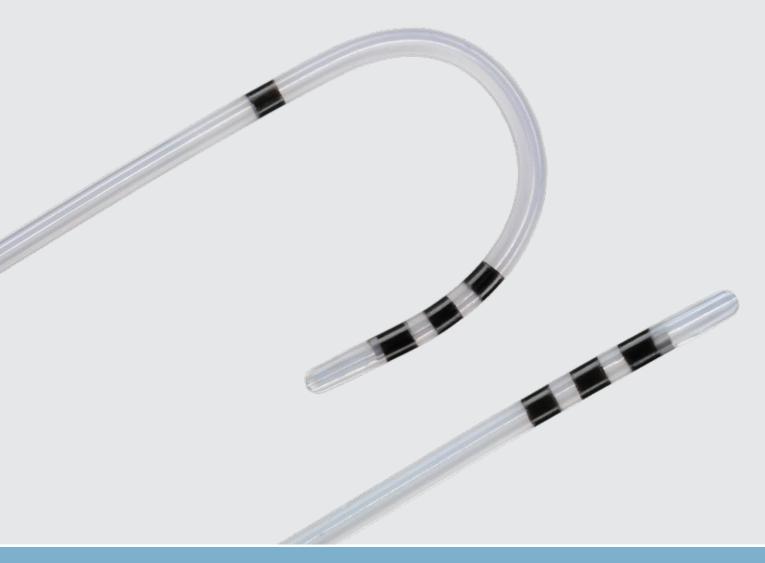

Glo-TipII[®] DOUBLE LUMEN ERCP CATHETER

This device is used for endoscopic cannulation of the ductal system.

Supplied sterile and is disposable - intended for single use only.

- DomeTip® designed with anatomy in mind for smooth, easy access
- Radiopaque tip provides excellent fluoroscopic visualization
- Double lumen design permits contrast injection through a dedicated lumen while catheter is placed over a wire guide
- Mark V System[®] with endoscopic markers at 5, 10, and 15 cm from tip to visualize depth of insertion

Order Number	Reference Part Number	Catheter Size Fr	Catheter Length cm	Tip	Remarks	Color Connection
G22226	GT-2-T	6	200	Dometip	Standard	Purple
G24694	GT-2-T-RB	6	200	Dometip	Radiopaque Bands	Purple
G24506	GT-2-T-ANG	6	200	Dometip	Angled	Purple

Minimum accessory channel 2.8 mm. Stylet length 100 cm. Recommend Acrobat® 2 Calibrated Tip Wire Guide, sold separately.

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service

EU Website: cookmedical.eu

EDI: cookmedical.eu/edi

Distributors: +353 61239240, ssc.distributors@cookmedical.com Austria: +43 179567121, oe.orders@cookmedical.com Belgium: +32 27001633, be.orders@cookmedical.com Denmark: +45 38487607, da.orders@cookmedical.com Finland: +358 972519996, fi.orders@cookmedical.com France: +33 171230269, fr.orders@cookmedical.com Germany: +49 6950072804, de.orders@cookmedical.com Hungary: +36 17779199, hu.orders@cookmedical.com Ireland: +353 61239252, ie.orders@cookmedical.com Italy: +39 0269682853, it.orders@cookmedical.com Netherlands: +31 202013367, nl.orders@cookmedical.com Norway: +47 23162968, no.orders@cookmedical.com Spain: +34 912702691, es.orders@cookmedical.com Sweden: +46 858769468, se.orders@cookmedical.com Switzerland - French: +41 448009609, fr.orders@cookmedical.com Switzerland - Italian: +41 448009609, it.orders@cookmedical.com Switzerland - German: +41 448009609, de.orders@cookmedical.com United Kingdom: +44 2073654183, uk.orders@cookmedical.com

COOK• MEDICAL USA Website: cookmedical.com EDI: cookmedical.com/edi.do Americas: Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335 E-mail: orders@cookmedical.com Australia: Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283 E-mail: cau.custserv@cookmedical.com-A4

AI, ESC, IR, OHNS, PI, RH, SUR-A4