


Seeing isn't just believing,
it's everything.


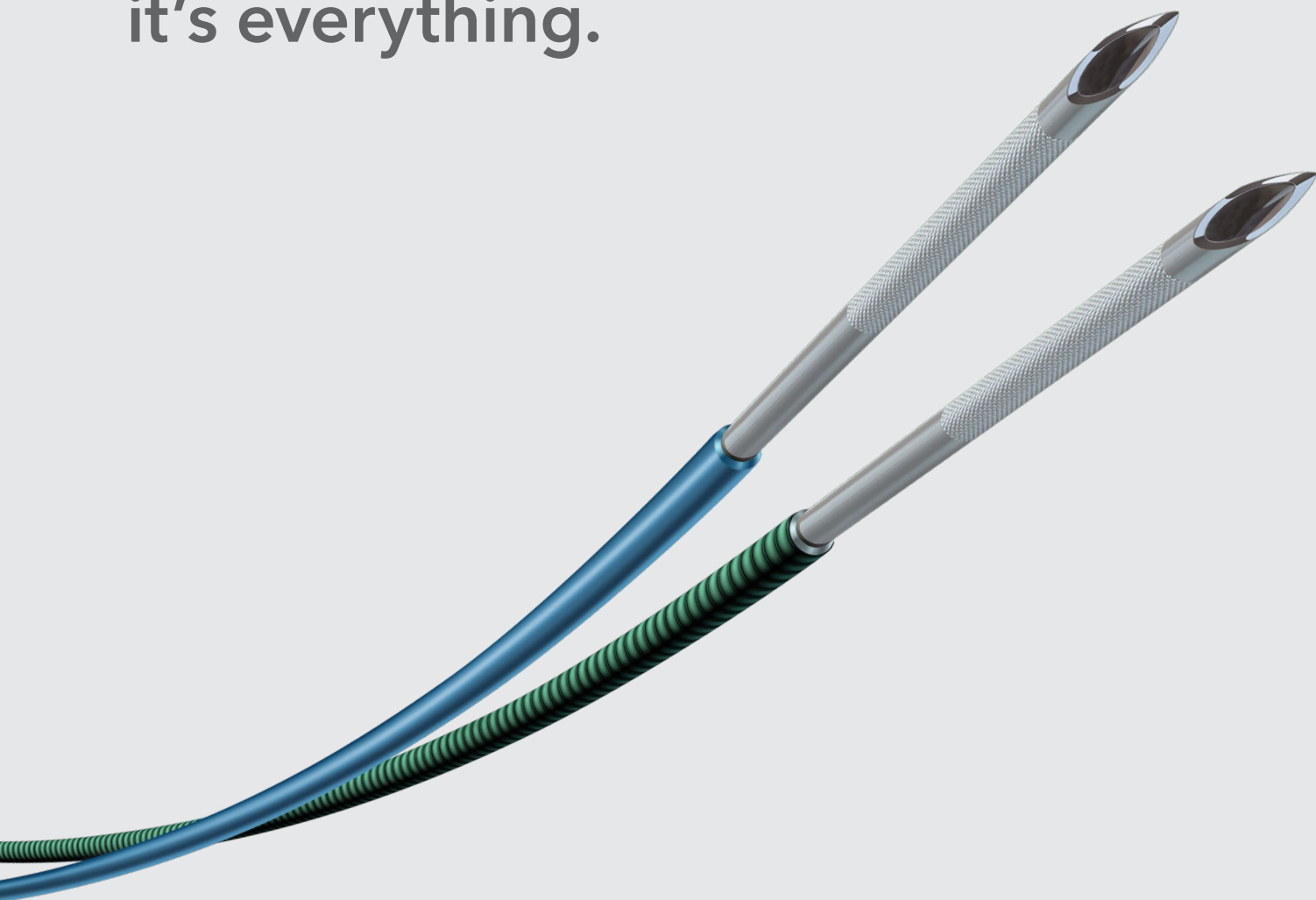


EchoTip[®] Ultra

ENDOSCOPIC ULTRASOUND NEEDLE

3 platforms – all with enhanced needle visualization

No matter what needle size your ultrasound procedure requires, a quality image is key to successful outcomes. That's why Cook offers the highly visible dimpling technology on our EchoTip Ultra needles.

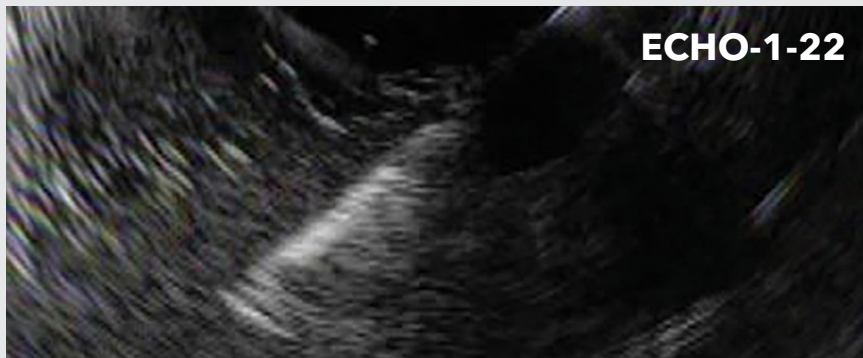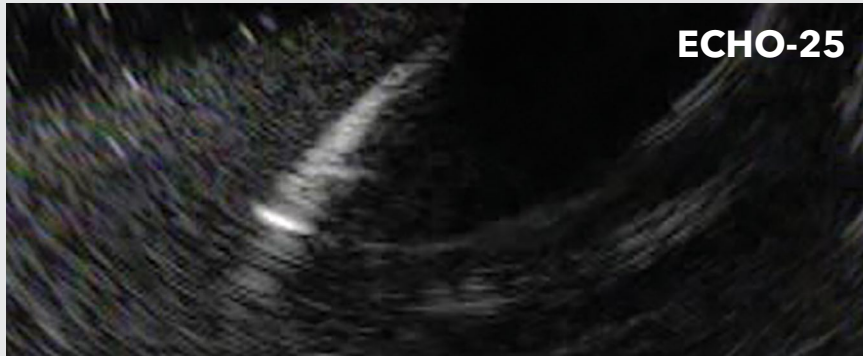


Balancing needle visualization, flexibility and ergonomics.


In today's complex, challenging EUS procedures, needle visualization, flexibility and precise control are more important than ever. We designed the EchoTip® Ultra to give you all three of these attributes.

Thanks to its highly visible needle and flexible coiled sheath*, the EchoTip® Ultra helps you reach difficult positions for sampling lesions.

Even when the device is torqued, the ergonomic handle gives you the stability you need for precise needle control.

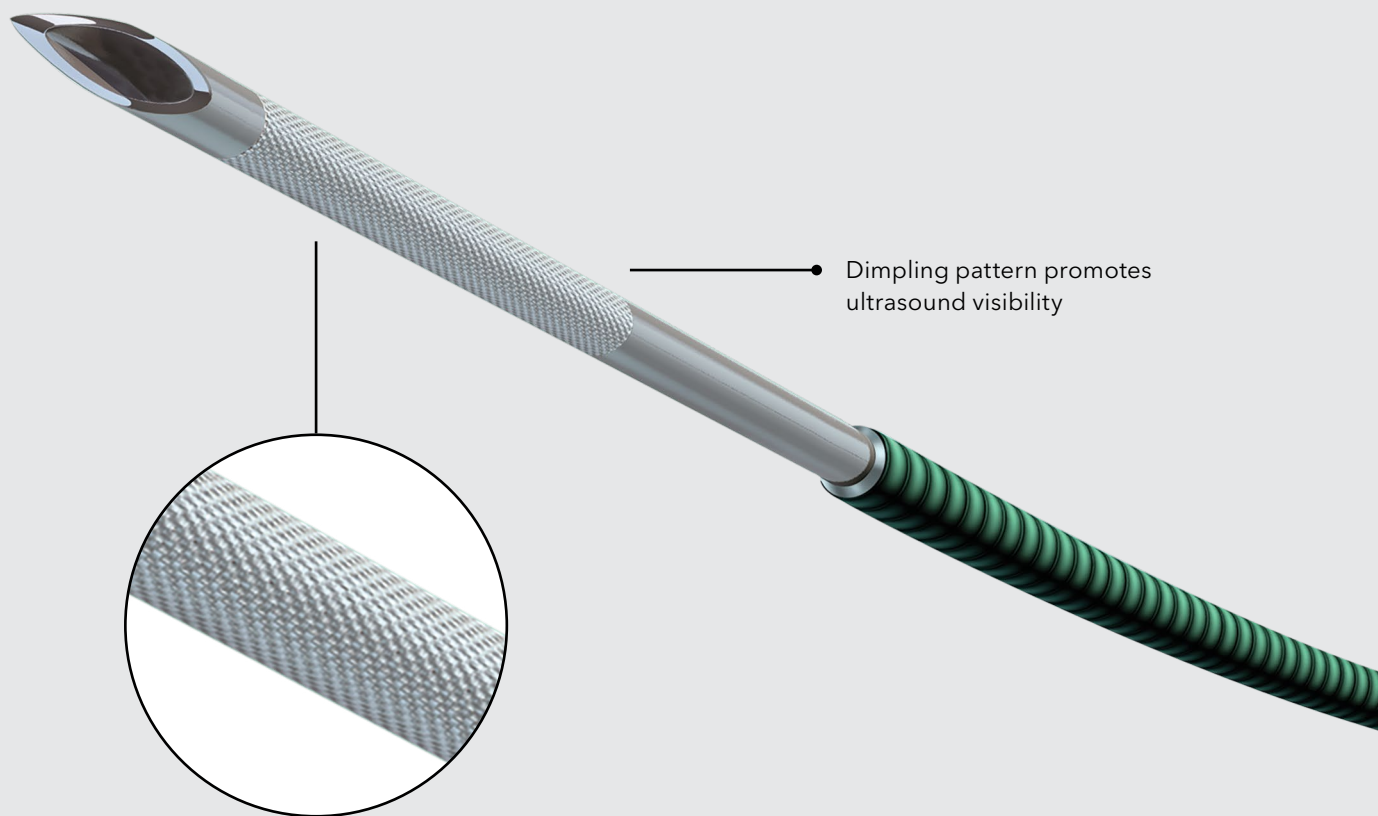

*Featured on the EchoTip Ultra 22 ga (ECHO-3-22) depicted above.

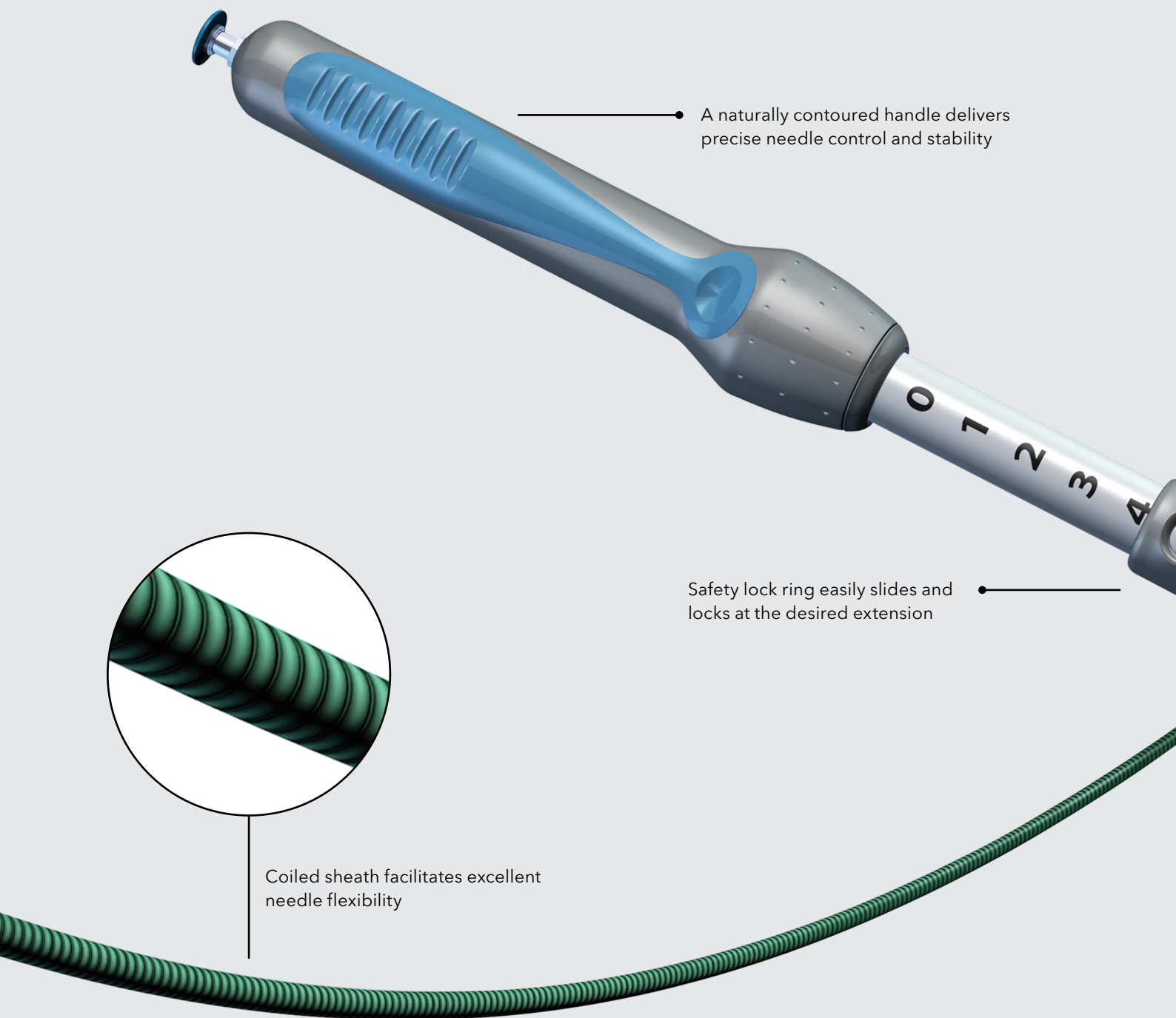


A naturally contoured handle delivers precise needle control and stability

Safety lock ring easily slides and locks at the desired extension

Coiled sheath facilitates excellent needle flexibility

Reaching challenging targets.

Even when you are targeting small lesions in difficult locations, EchoTip Ultra's flexible sheath and naturally contoured handle give you the control and stability you need. And that can help you confidently perform any EUS procedure referred to you.

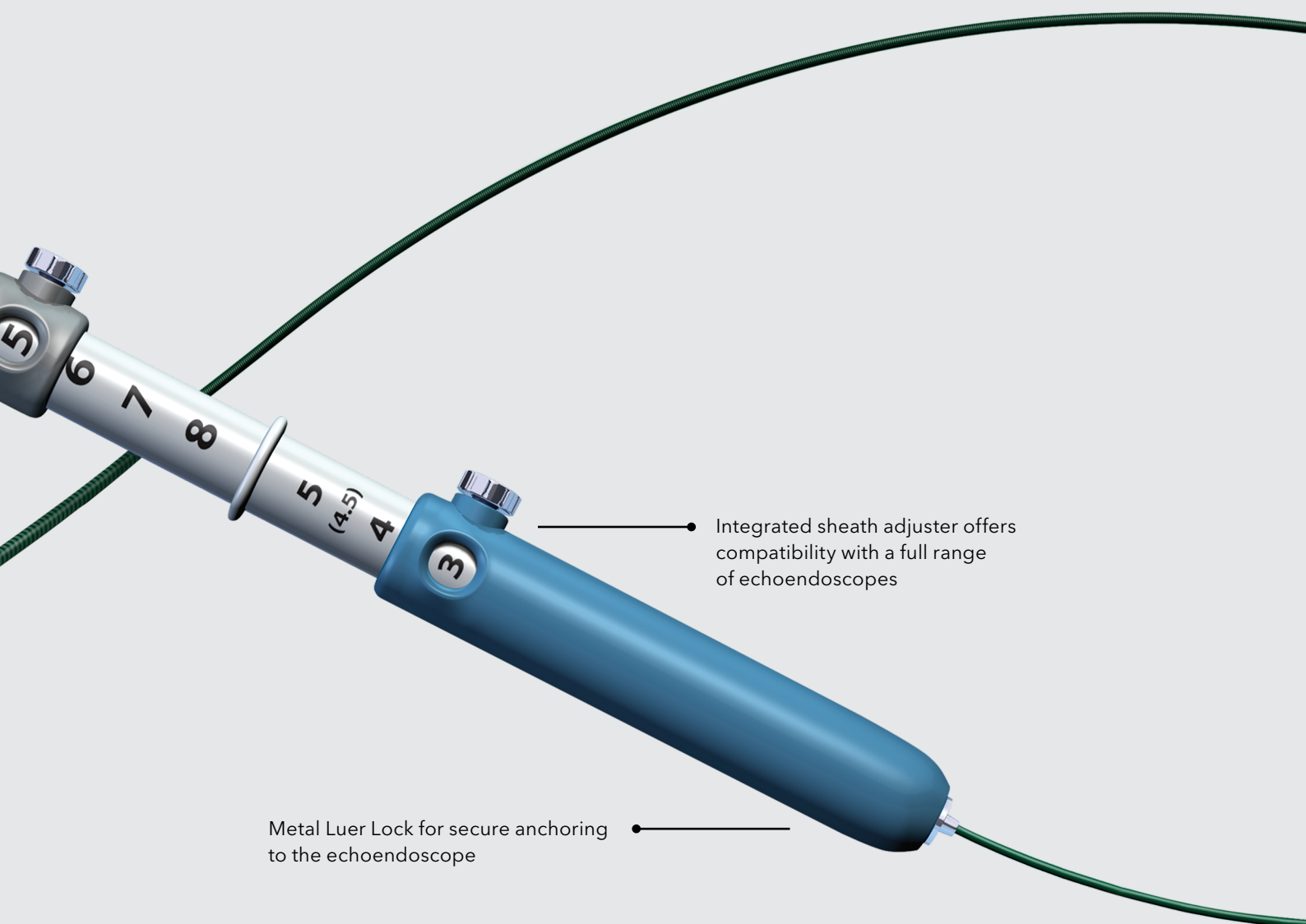

Integrated sheath adjuster offers compatibility with a full range of echoendoscopes

Metal Luer Lock for secure anchoring to the echoendoscope

Product Specifications

Used to sample targeted submucosal gastrointestinal lesions through the accessory channel of an ultrasound endoscope.

Supplied sterile and is disposable - intended for single use only.

Order Number	Reference Part Number	Needle Size gauge	Sheath Size Fr	Adjustable Needle Extension cm	Stylet Tip	Minimum Accessory Channel mm
G31520	ECHO-19	19	5.2-4.3	0-8	ball tip	2.8
G31518	ECHO-1-22	22	5.2	0-8	ball tip	2.8
G31521	ECHO-3-22*	22	5.4	0-8	ball tip	2.8
G31519	ECHO-25	25	5.2	0-8	ball tip	2.8

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details regarding: stock/non-stock status (expect longer lead time if non-stock), pricing and/or availability.

**This device is indicated for adult use only.*

Customer Service

EU Website: cookmedical.eu
EDI: cookmedical.eu/edi
Distributors: +353 61239240, ssc.distributors@cookmedical.com
Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001702, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Iceland: +354 800 7615, is.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Norway: +47 23162968, no.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: orders@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283

E-mail: cau.custserv@cookmedical.com



AI-ESC-IR-OHNS-PI-RH-SUR-A4