


Directly measure portal
pressure under EUS.

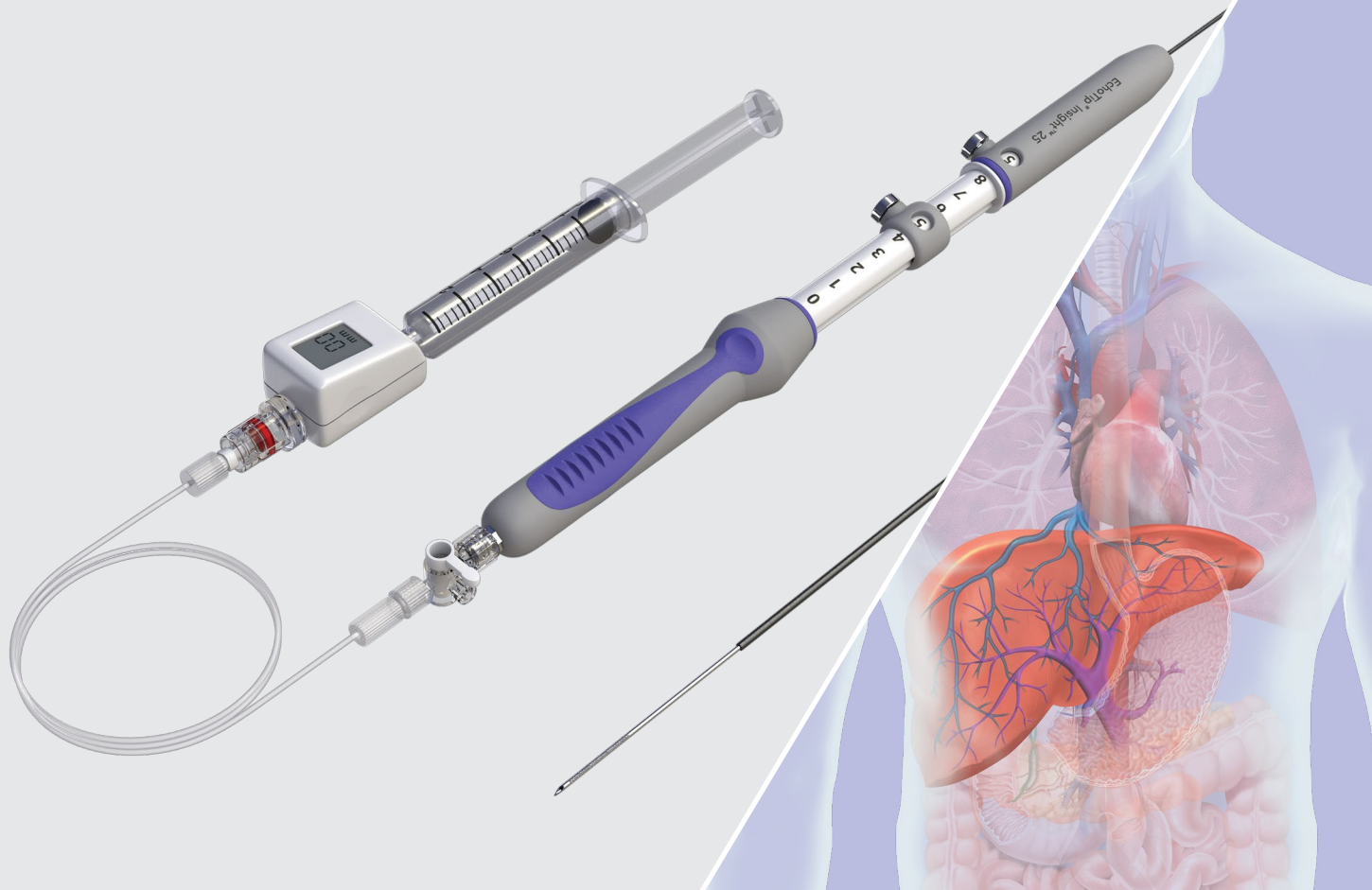

EchoTip[®] Insight[™]

PORTOSYSTEMIC PRESSURE GRADIENT MEASUREMENT SYSTEM

Accessible portal pressure measurement for endohepatology


The EchoTip Insight system allows endoscopists to directly measure portal pressure under endoscopic ultrasound (EUS). A new, minimally invasive procedure which now allows for comprehensive liver assessments through endoscopy procedures.


Direct Measurement

With the EchoTip Insight system, you can now directly measure the pressure in the portal and hepatic veins. Under EUS guidance the average of three measurements in each vein is used to calculate the gradient.

Accurate Measurement

The transducer allows you to accurately measure the portal pressure gradient so you can determine the appropriate treatment pathway for your patient.

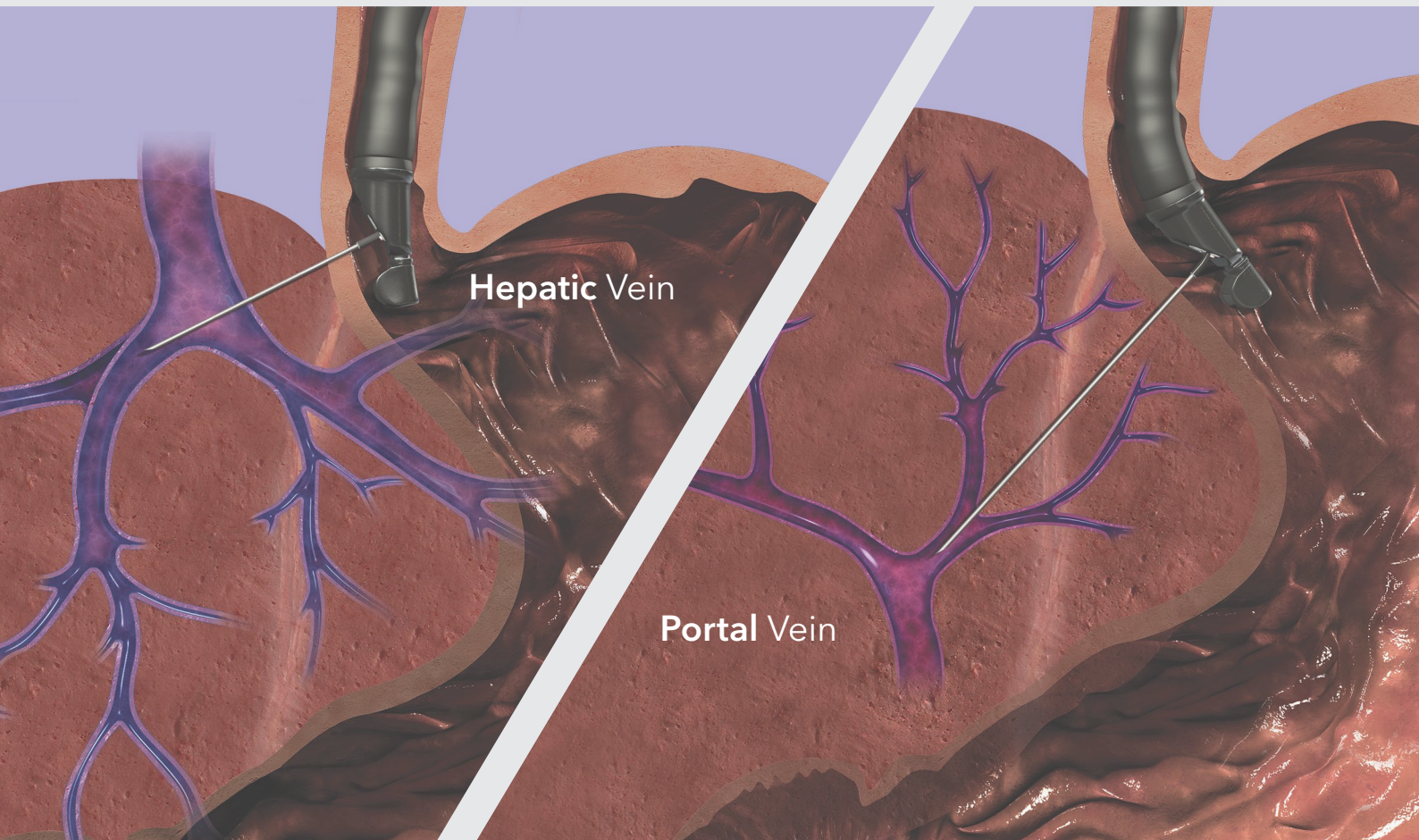

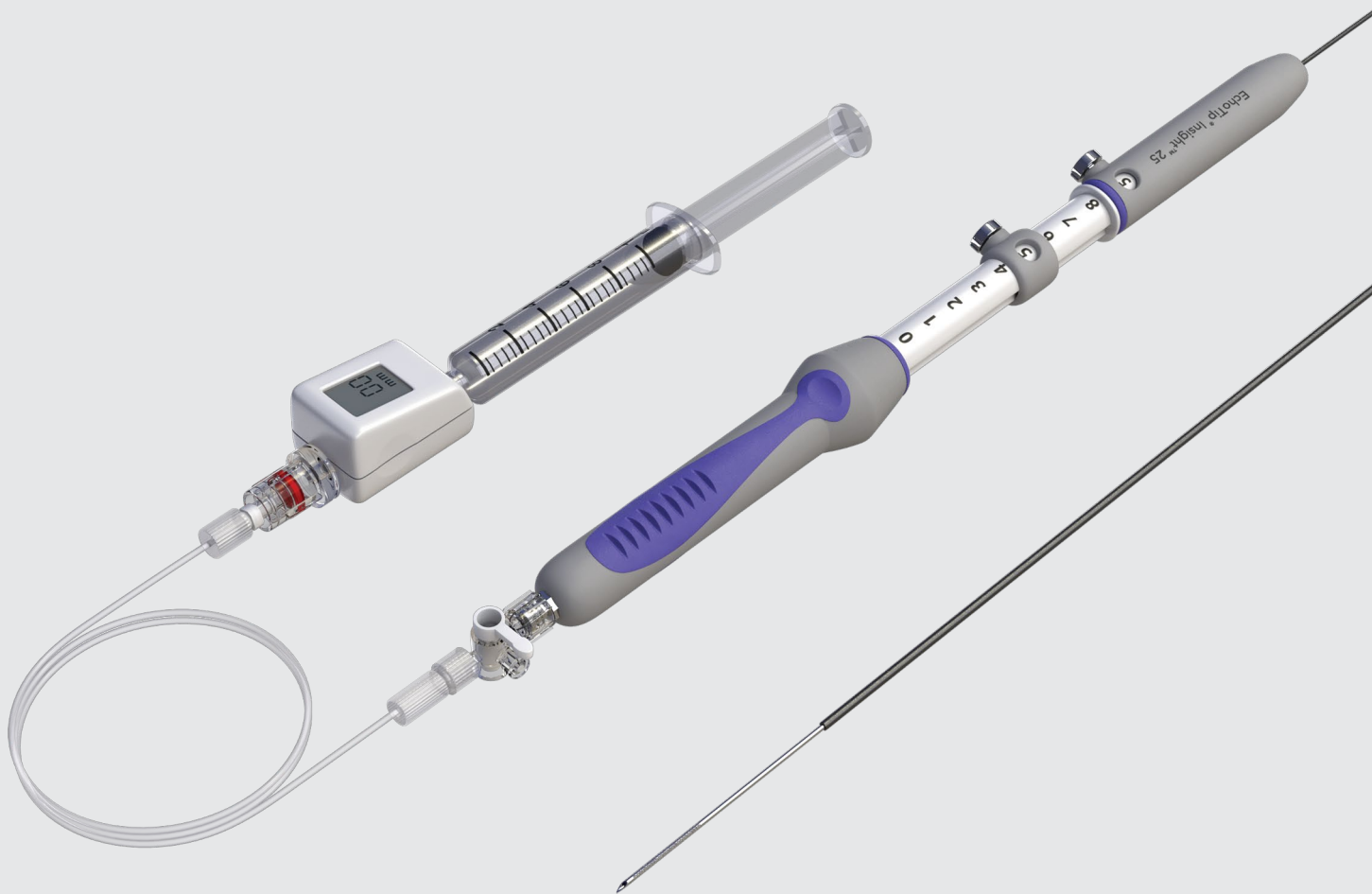

Needle

- Echogenic tip - so you can ensure that the needle tip is within the target
- Small gauge needle (25 gauge) - to facilitate access of the portal and hepatic veins

Transducer

- Digital display of the pressure reading
- Self-calibrating transducer
- Compact size

EchoTip[®] Insight[™]

PORTOSYSTEMIC PRESSURE GRADIENT MEASUREMENT SYSTEM

Product Specifications

Used to directly measure pressures in the hepatic and portal venous vasculatures and is used in conjunction with an ultrasound endoscope.

Order Number	Reference Part Number	Needle Size Gauge	Sheath Size Fr	Adjustable Needle Extension cm	Connecting Tube Length cm	Pressure Range mm Hg	Minimum Accessory Channel mm
G46560	ECHO-PPG	25	5.2	8	90	-199 to +999	2.8

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Vista® | Education and Training

To learn more about EchoTip Insight please contact your local Cook representative or visit www.cookmedical.com.

If you would like to attend one of our Vista Education and Training workshops, please find information on our Vista website: vista.cookmedical.com.

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001702, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Iceland: +354 800 7615, is.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com
 EDI: cookmedical.com/edi.do
Americas:
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: orders@cookmedical.com
Australia:
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 E-mail: cau.custserv@cookmedical.com



AI-ESC-IR-OHNS-PI-RH-SUR-A4