


Directly measure portal
pressure under EUS.



EchoTip[®] Insight[™]
PORTOSYSTEMIC PRESSURE GRADIENT

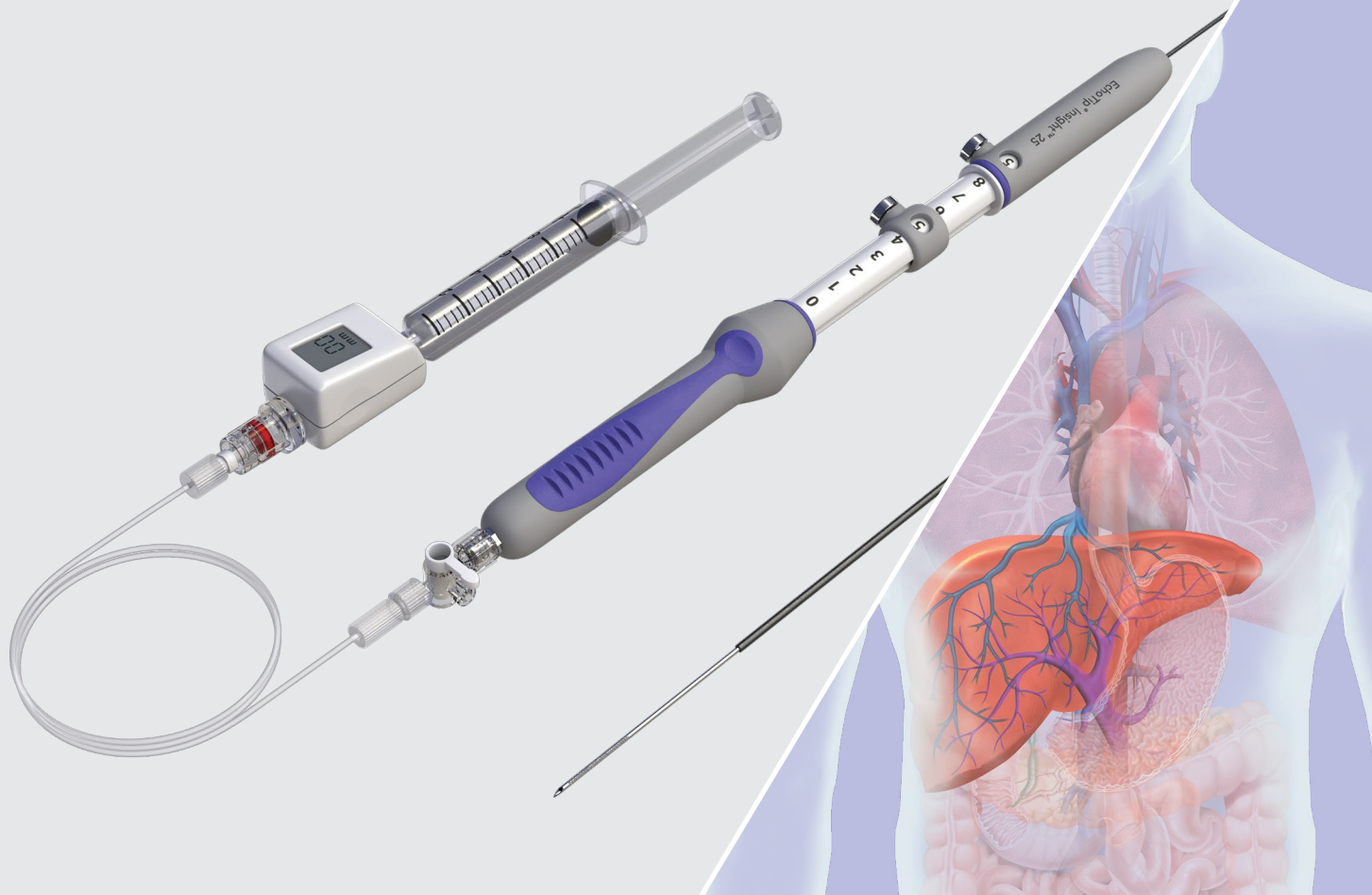


EchoTip[®] Insight[™]

PORTOSYSTEMIC PRESSURE GRADIENT

The EchoTip Insight Needle and Connecting Tube are packaged together with the Compass CT transducer in an outer box. Once assembled, these devices are used with an ultrasound endoscope in patients with portal hypertension.


Compass CT Pressure Transducer:

- Digital display of the pressure reading
- Self-calibrating transducer
- Compact size

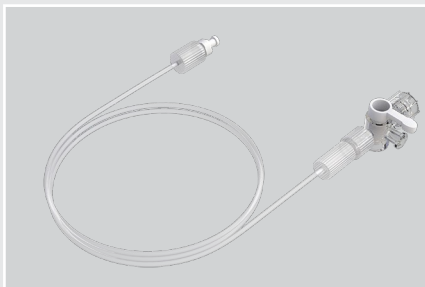

Connecting Tube:

- 90 cm long connecting tube allowing the transducer to be positioned level with phlebostatic axis

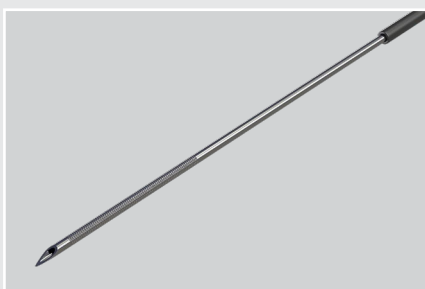

Endoscopic Ultrasound Needle:

- Echogenic tip - so you can ensure that the needle tip is within the target
- Small gauge needle (25 gauge) - to facilitate access of the portal and hepatic veins

Order Number	Reference Part Number	Needle Size Gauge	Sheath Size Fr	Adjustable Needle Extension cm	Connecting Tube Length cm	Pressure Range mm Hg	Minimum Accessory Channel mm
G58000	ECHO-PPG-I	25	5.2	8	90	-199 to +999	2.8

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details. Please see product risk information in the IFU at cookmedical.eu

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001633, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: orders@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283

E-mail: cau.custserv@cookmedical.com



AI, ESC, IR, OHNS, PI, RH, SUR-A4