


Choose strength.

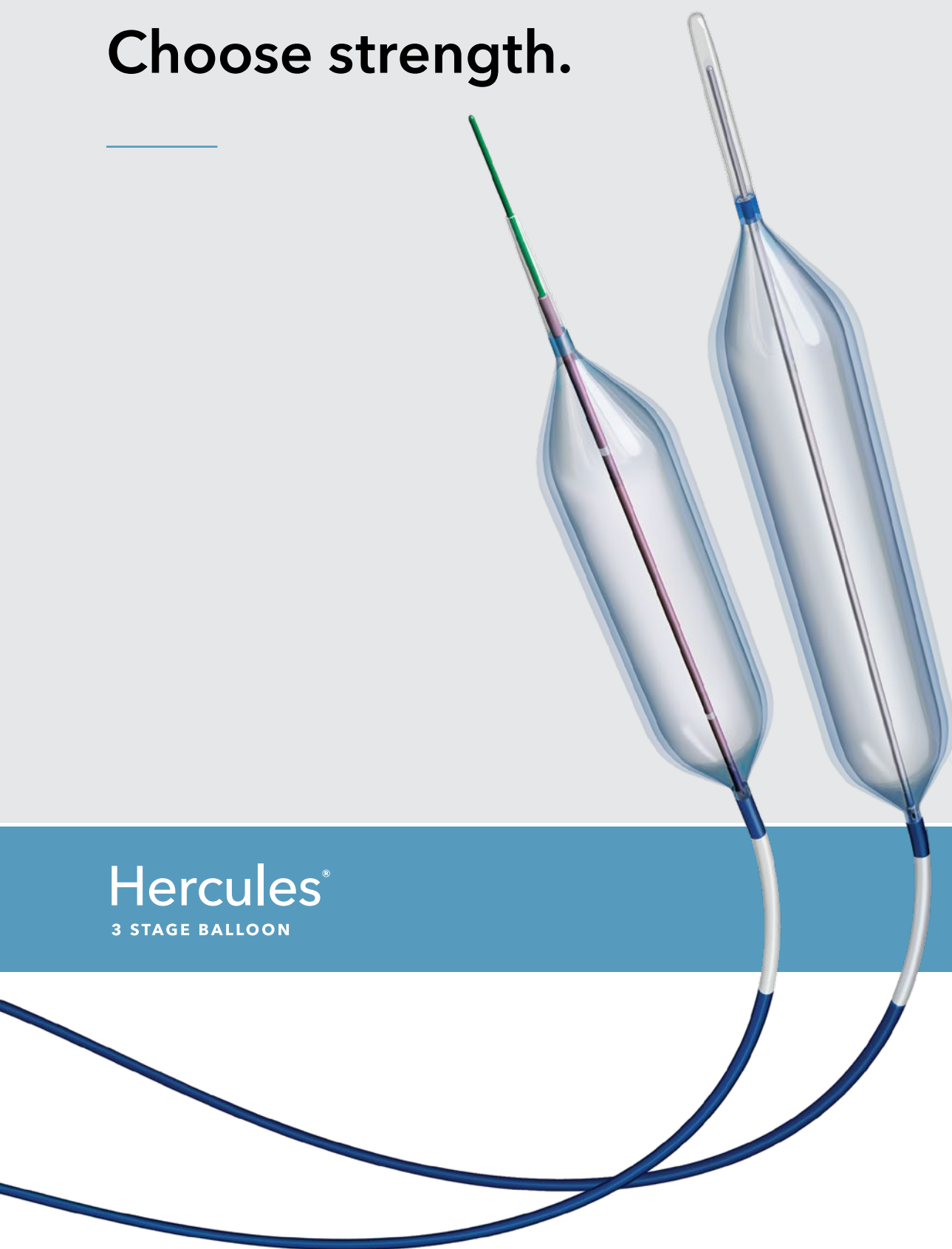

Hercules[®]
3 STAGE BALLOON

Hercules®

3 STAGE BALLOON

Hercules balloons are made of PET—a durable material that has the strength to stand up to a GI stricture. With PET, Hercules balloons not only exert radial force on the stricture but also maintain their position during the dilation procedure.

The flexible, atraumatic tip provides access through tight strictures.


A clear balloon allows direct visualisation of the dilation process.

Two radiopaque markers enhance the visibility of the balloon positioning.

Rapid balloon deflation is achieved with the express-evacuation catheter.


Rapid balloon deflation is achieved with the express-evacuation catheter and rapid deflation sleeve.

The kink-resistant nitinol catheter provides manoeuvrability.

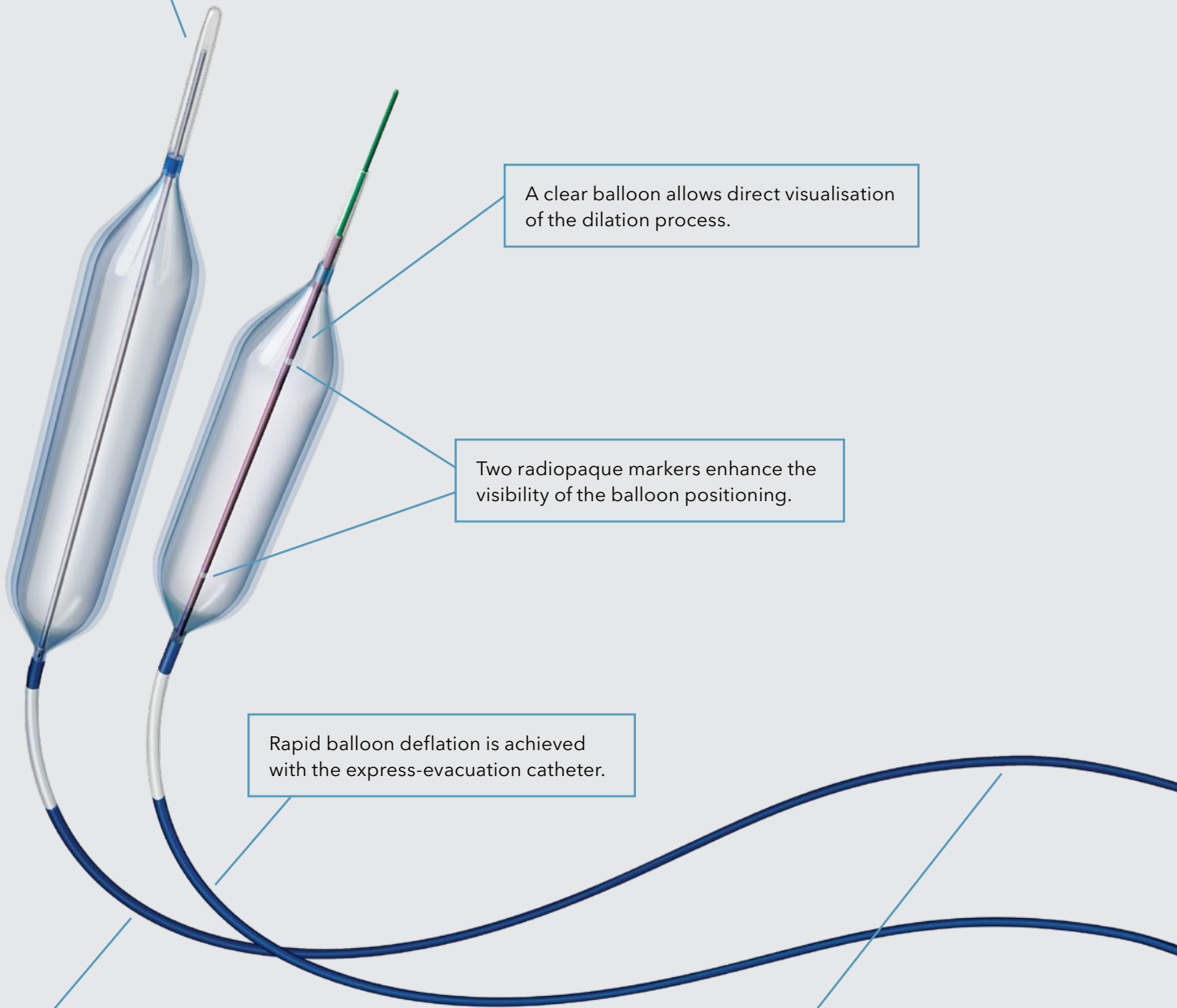


Clinical strength

Hercules balloons have the dilating strength to relieve your patient's stricture in the GI tract, including in the oesophagus, pylorus, duodenum, and colon, and may support longer-term luminal patency.

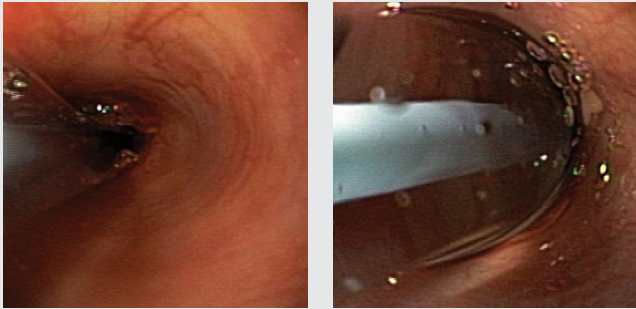


Rapid deflation

Hercules is engineered to facilitate rapid deflation after a dilation procedure. By applying negative pressure, you can deflate the balloon quickly with the express-evacuation catheter* and rapid deflation sleeve**.

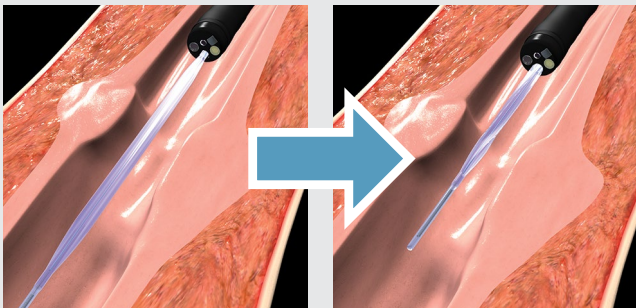

The size and inflation information on the catheter tag allow for easy identification.

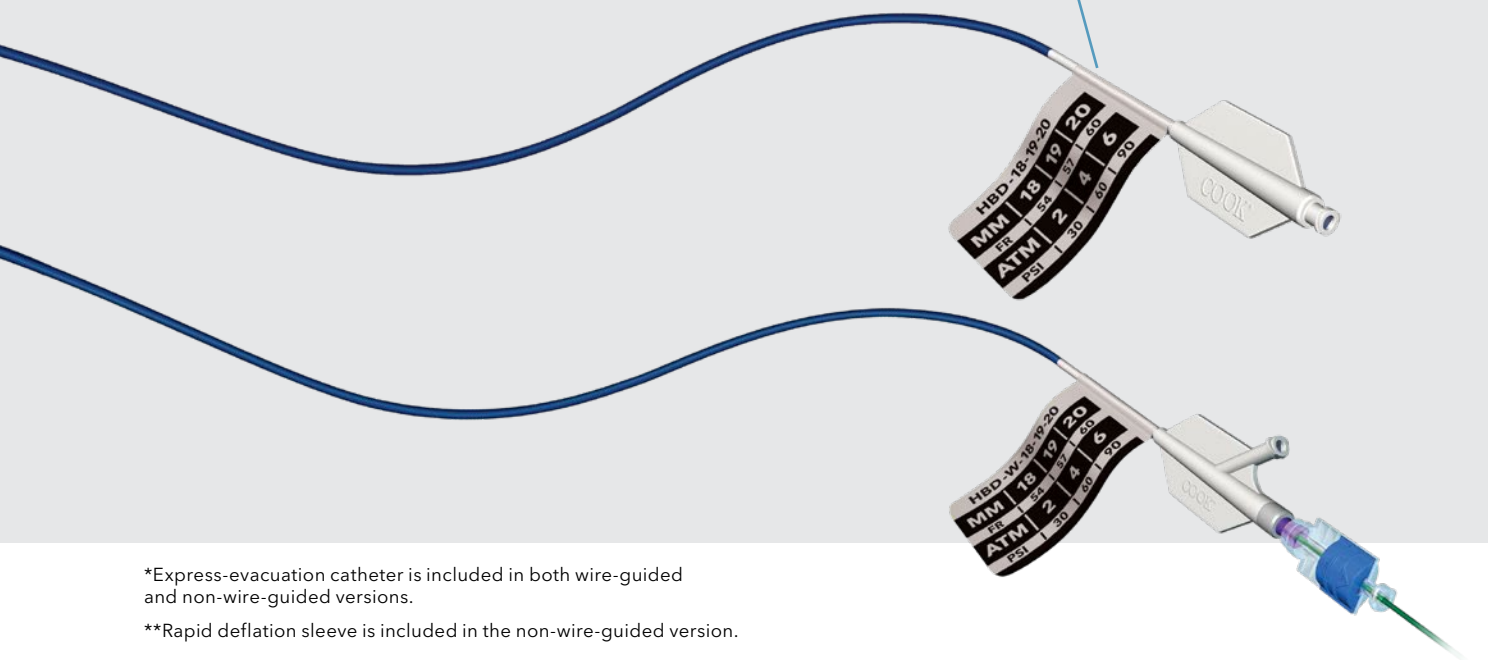

*Express-evacuation catheter is included in both wire-guided and non-wire-guided versions.

**Rapid deflation sleeve is included in the non-wire-guided version.

Hercules® 3 Stage Balloon

Used to endoscopically dilate strictures of the oesophagus.

Order Number	Reference Part Number	Catheter Fr	Working Length cm	Balloon Diameter mm	Balloon Fr	Balloon Length cm	Pressure atm	Pressure psi	Minimum Accessory Channel Diameter mm
G31928	HBD-8-9-10	7	174	8-9-10	24-27-30	8	4-6-8	60-90-120	2.8
G31925	HBD-10-11-12	7	174	10-11-12	30-33-36	8	2-4-6	30-60-90	2.8
G31926	HBD-12-13.5-15	7	174	12-13.5-15	36-40.5-45	8	2-4-6	30-60-90	2.8
G48732	HBD-15-16.5-18	7	174	15-16.5-18	45-49.5-54	8	2-4-6	30-60-90	2.8
G31927	HBD-18-19-20	7	174	18-19-20	54-57-60	8	2-4-6	30-60-90	2.8

Hercules® 3 Stage Wire Guided Balloon

Used to endoscopically dilate strictures of the gastrointestinal tract including the oesophagus, pylorus, duodenum, and colon.

Order Number	Reference Part Number	Catheter Fr	Working Length cm	Balloon Diameter mm	Balloon Fr	Balloon Length cm	Pressure atm	Pressure psi	Minimum Accessory Channel Diameter mm
G51799	HBD-W-8-9-10	8	240	8-9-10	24-27-30	5.5	4-6-8	60-90-120	2.8
G51795	HBD-W-10-11-12	8	240	10-11-12	30-33-36	5.5	2-4-6	30-60-90	2.8
G51796	HBD-W-12-13.5-15	8	240	12-13.5-15	36-40.5-45	5.5	2-4-6	30-60-90	2.8
G51797	HBD-W-15-16.5-18	8	240	15-16.5-18	45-49.5-54	5.5	2-4-6	30-60-90	2.8
G51798	HBD-W-18-19-20	8	240	18-19-20	54-57-60	5.5	2-4-6	30-60-90	2.8

Wire guide length is 300 cm.

Cook Sphere® Inflation Device

Used to inflate the balloon with fluid, to monitor the pressure in the balloon during the procedure, and to deflate the balloon after dilation.

Order Number	Reference Part Number	Fluid Capacity mL	Maximum Pressure atm
G34903	CID-60-15	60	15

Please follow any step-by-step instructions carefully, including the risk, technical information, Device Name and intended use provided in the IFU at cookmedical.eu. Note that some products or part numbers may not be available in all markets; consult your Cook representative for details.

Customer Service

EU website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, emea.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001702, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 97251996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 21612723770, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Iceland: +354 8007615, is.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA website: cookmedical.com
 EDI: cookmedical.com/edi.do

Americas:
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 Email: customersupport@cookmedical.com

Australia:
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 Email: cau.custserv@cookmedical.com



AI, ESC, IR, PI-A4