

Give bowel control back
to your patients.



Chait


PERCUTANEOUS CECOSTOMY CATHETER



Chait

PERCUTANEOUS CECOSTOMY CATHETER

- Placing the catheter is a minimally invasive, reversible procedure.
- The catheter facilitates antegrade colonic cleansing for patients who have fecal incontinence.
- The unique, low-profile trapdoor fitting is designed to lie flat against the skin.
- The catheter enables patients to manage their fecal incontinence themselves.
- The catheter is 10.2 Fr, features three sideports, and is available in three lengths.

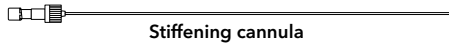

Stiffening cannula



Chait cecostomy catheter
Polyurethane



Chait access adapter with connecting tube

Order Number	Reference Part Number	Tract Length cm	Extended Catheter Length cm	Catheter Fr
Chait Percutaneous Cecostomy Catheter				
G11632	TDCS-100	0-6	20	10.2
G36392	TDCS-100-M	3-9	24	10.2
G12685	TDCS-100-L	6-14	28	10.2
Chait Access Adapter with Connecting Tube				
G11630	CTU10.0-18-CAA	10.2	-	10.0

Accessories

Dawson-Mueller Multipurpose Drainage Catheter

The drainage catheter is placed during initial cecostomy as a temporary catheter, allowing the tract to mature. After approximately six weeks, this catheter is exchanged for the Chait Percutaneous Cecostomy Catheter.

- The Mac-Loc® locking loop helps prevent catheter displacement.
- Soft, smooth Ultrathane® material helps maintain patient comfort.
- The hydrophilic coating eases insertion.
- The catheter is 8.5 Fr and 25 cm long and is packaged with a trocar stylet, a flexible stiffener, fixation device, and an introduction stiffening cannula with loading obturator.


Enterostomy Suture Anchor Set

Components

- Introducer needle, 18 gage, 7 cm long
- Loading cartridge containing two suture anchors
- Amplatz Extra-Stiff Wire Guide, 0.035 inch (0.89 mm) in diameter, PTFE coated, 80 cm long
- Coons dilator, 8 Fr, radiopaque, 20 cm long


Order Reference Number Part Number

Dawson-Mueller Multipurpose Drainage Catheter

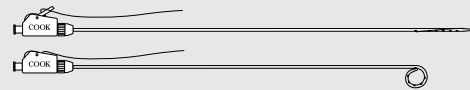


G09705 ULT8.5-38-25-P-5S-CLDM-HC

Enterostomy Suture Anchor Set

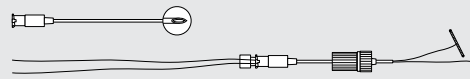

G12230 GIAS-100-CHAIT

Amplatz Extra-Stiff Wire Guide

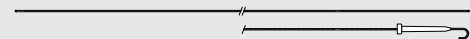

G04003 THSF-35-80-AES

Coons Dilator



G03927 JCD8.0-38-20-COONS

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001702, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Iceland: +354 800 7615, is.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: customersupport@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 E-mail: cau.custserv@cookmedical.com



AI-ESC-IR-OHNS-PI-RH-SUR-A4