

Retrieve with ease



Over the wire (OTW)
Compatible with 0.035 inch
diameter wire guides



Four-loop configuration
Interconnected
nitinol loops




Platinum markers
Enhanced visibility



Radiopaque band
Highly visible



Flexor® technology
Kink resistant

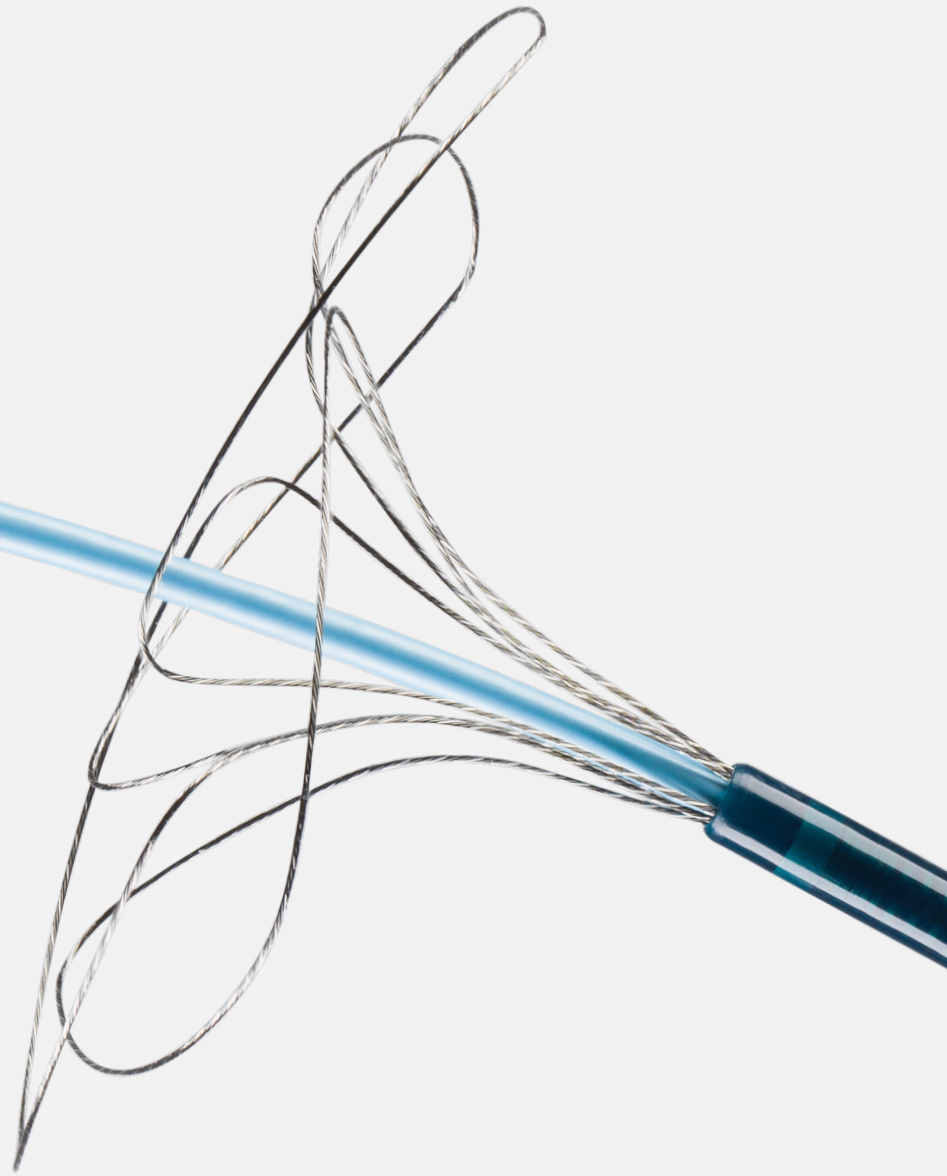

Indy OTW®
VASCULAR RETRIEVER



Indy OTW Vascular Retriever

Cook Medical offers the only four-loop snare on the market. Its unique, interconnected four-loop snare design with over-the-wire capability allows greater access to vessel coverage, reducing the need for device exchanges and saving you valuable time.

Other unique features include kink-resistant Flexor technology and visibility features throughout the entire device.

Over the wire (OTW)

Compatible with 0.035 inch diameter wire guides



Product features

Over the wire (OTW)

The Indy OTW snare is easier to deliver and has enhanced stability compared to non-OTW snares, reducing the need for additional device exchanges and saving time.

Four-loop configuration

The four-loop Indy OTW snare fully expands to 40 mm diameter, facilitating vessel coverage and ensuring a confident retrieval.

Platinum markers


Nitinol loops with platinum markers retain memory shape, delivering dependable performance, enhanced visibility, and reduced need for multiple device exchanges.

Radiopaque band

The radiopaque band markers indicate the precise location of the distal tip for visibility.

Flexor technology


The Flexor sheath features kink-resistant Flexor technology, providing maximum trackability across tortuous anatomy.

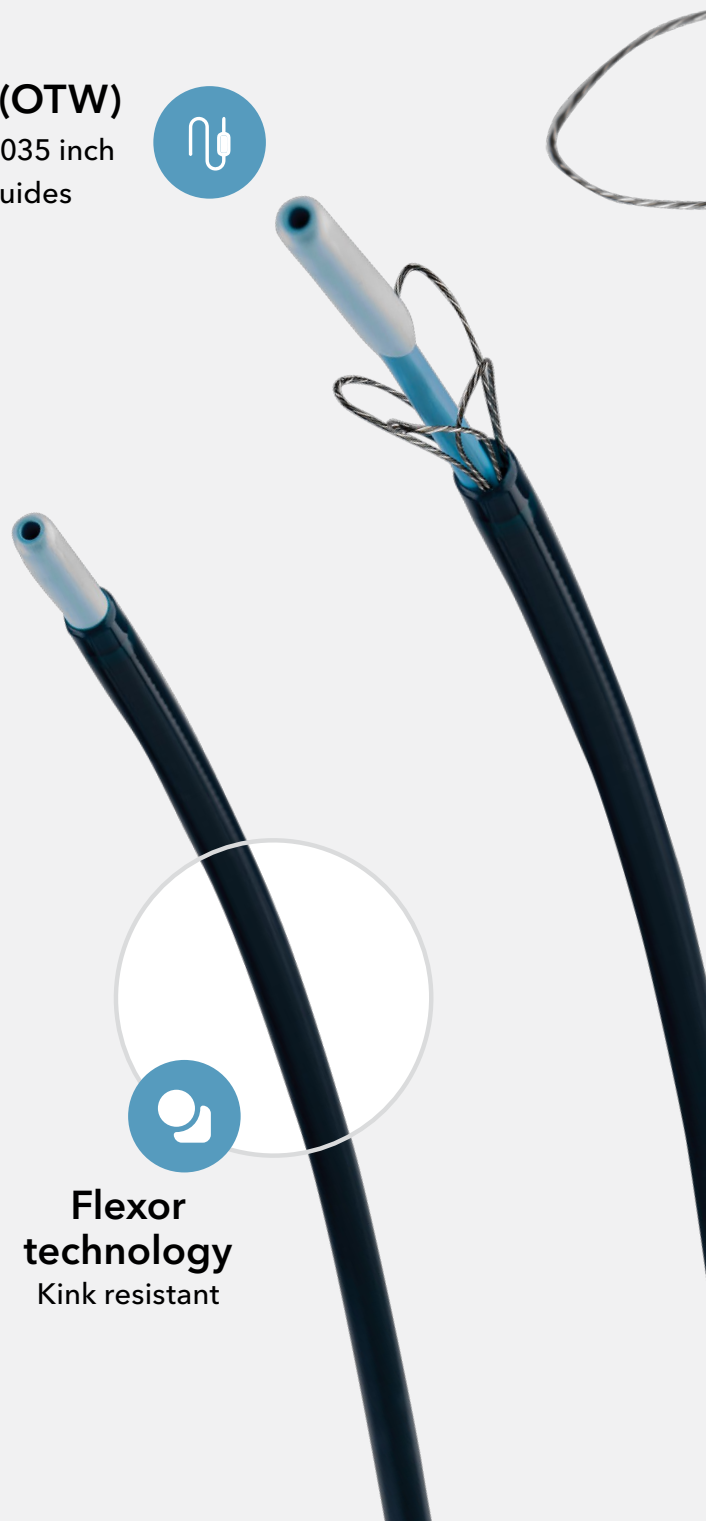

Flexor technology
Kink resistant

Four-loop configuration

Interconnected nitinol loops

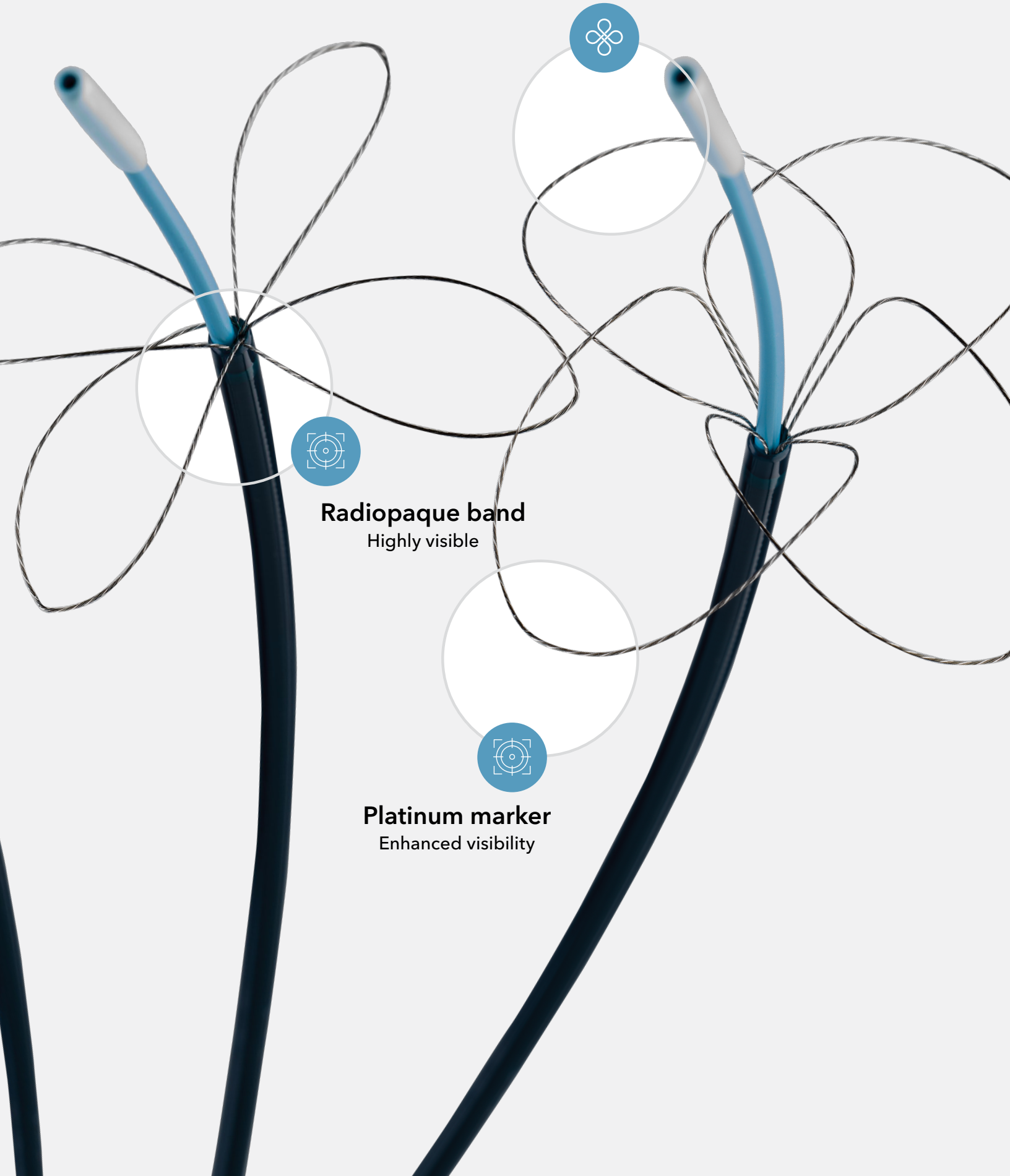


Radiopaque band

Highly visible



Platinum marker

Enhanced visibility

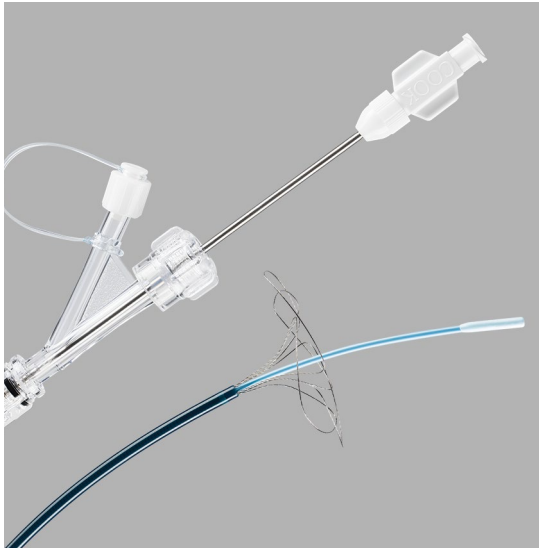

Indy OTW[®] Vascular Retriever

With its distal tip and radiopaque tubing, the Indy OTW snare provides greater access, allowing you to target and contain a foreign body and withdraw it to a peripheral vascular location.

Scan for more information.



Order Number	Reference Part Number (RPN)	Sheath Outer Diameter Fr	Accepts Wire Guide Diameter inch	Sheath cm	Basket Diameter mm
G51835	INDY-8.0-35-55-40	8	0.035	55	40
G51836	INDY-8.0-35-100-40	8	0.035	100	40

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001702, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Iceland: +354 8007615, IS.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com
 EDI: cookmedical.com/edi.do
Americas:
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: customersupport@cookmedical.com
Australia:
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 E-mail: cau.custserv@cookmedical.com



AI, ESC, IR, OHNS, PI, RH, SUR-8.5X11