

**Brilliant visibility.
Lubricious navigation.**



Slip-Cath® Beacon® Tip
HYDROPHILIC ANGIOGRAPHIC CATHETER



Learn more
about Slip-Cath.





Visibility

Brilliant under fluoroscopy

Visibility

Under fluoroscopy, the easy-to-see tip lights the way through every turn in the body's anatomy.

Trackability

The hydrophilic coating and tapered tip facilitate vessel selection, ranging from simple vasculature to more challenging anatomy. The increased lubricity allows the user to navigate the vasculature system smoothly.

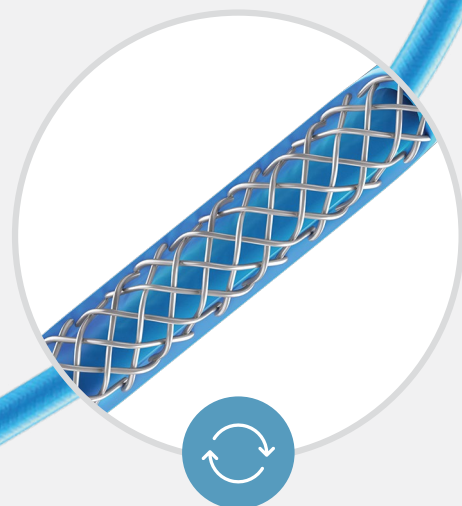
Torqueability

A nylon catheter construction with stainless steel braiding through the shaft provides control, pushability, and torque response.



Trackability

Navigate tortuous anatomy



Torqueability

Kink resistant

Flow rates

Fr	Catheter Length cm	MAXIMUM FLOW RATE		MAXIMUM RECOMMENDED PRESSURE	
		Flow Rate mL/sec	kg/cm ²	psi	
5.0	40	13	82	1200	
5.0	65	11	82	1200	
5.0	80	9	82	1200	
5.0	90	9	82	1200	
5.0	100	9	82	1200	
5.0	125	7	82	1200	
5.5	75	9	82	1200	
5.5	125	7	82	1200	
6.5	125	7	82	1200	

Packaging

The foil packaging protects device integrity.

- Updated packaging protects catheter tips from separating and is safer for patients.
- The foil pouch serves as an additional protective layer against environmental factors.
- A plastic back slit protects against kinking during transportation.

Three-year
shelf life















Protection against
light and humidity damage

Slip-Cath Beacon Tip Hydrophilic Angiographic Catheter



Order Number	Reference Part Number (RPN)	Fr	Accepts Wire Guide Diameter inch	Length cm	Tip Configuration	
G23688	SCBR5.0-38-65-P-NS-DAV	5.0	0.038	65		DAV
G11996	SCBR5.0-38-100-P-NS-DAV	5.0	0.038	100		
G13413	SCBR5.0-38-125-P-NS-DAV	5.0	0.038	125		
G12716	SCBR5.0-38-65-P-NS-SIM1	5.0	0.038	65		SIM1
G12000	SCBR5.0-38-100-P-NS-SIM1	5.0	0.038	100		
G32369	SCBR5.0-38-125-P-NS-SIM1	5.0	0.038	125		
G23691	SCBR5.0-38-65-P-NS-SIM2	5.0	0.038	65		SIM2
G11992	SCBR5.0-38-100-P-NS-SIM2	5.0	0.038	100		
G13200	SCBR5.0-38-125-P-NS-SIM2	5.0	0.038	125		
G12007	SCBR5.0-38-100-P-NS-VERT	5.0	0.038	100		VERT
G12168	SCBR5.0-38-125-P-NS-VERT	5.0	0.038	125		
G13319	SCBR5.0-38-65-P-NS-H1	5.0	0.038	65		H1
G11991	SCBR5.0-38-100-P-NS-H1	5.0	0.038	100		
G23696	SCBR5.0-38-125-P-NS-H1	5.0	0.038	125		
G23693	SCBR5.0-38-100-P-NS-HN4	5.0	0.038	100		HN4
G11997	SCBR5.0-38-100-P-NS-JB1	5.0	0.038	100		JB1
G13530	SCBR5.0-38-125-P-NS-JB1	5.0	0.038	125		
G11998	SCBR5.0-38-100-P-NS-JB2	5.0	0.038	100		JB2
G11999	SCBR5.0-38-100-P-NS-MAN	5.0	0.038	100		MAN
G12001	SCBR5.0-38-100-P-NS-SIM3	5.0	0.038	100		SIM3
G13163	SCBR5.0-38-100-P-NS-VTK	5.0	0.038	100		VTK
G12082	SCBR5.0-38-125-P-NS-VTK	5.0	0.038	125		
G36214	SCBR5.0-38-80-P-NS-RH	5.0	0.038	80		RH

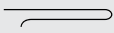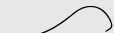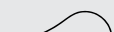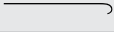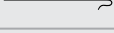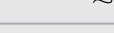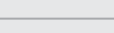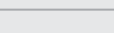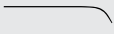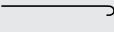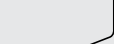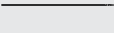Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Distribution for details.

Order Number	Reference Part Number (RPN)	Fr	Accepts Wire Guide Diameter inch	Length cm	Tip Configuration	
G59854	SCBR5.0-35-65-P-NS-RBT	5.0	0.035	65		RBT
G59855	SCBR5.0-35-90-P-NS-RBT	5.0	0.035	90		
G59856	SCBR5.0-38-65-P-NS-RBT	5.0	0.038	65		
G59857	SCBR5.0-38-90-P-NS-RBT	5.0	0.038	90		
G12003	SCBR5.0-38-65-P-NS-C1	5.0	0.038	65		C1
G11994	SCBR5.0-38-65-P-NS-C2	5.0	0.038	65		C2
G12006	SCBR5.0-38-80-P-NS-C2	5.0	0.038	80		
G12008	SCBR5.0-38-100-P-NS-C2	5.0	0.038	100		
G12004	SCBR5.0-38-65-P-NS-C3	5.0	0.038	65		C3
G13272	SCBR5.0-38-65-P-NS-RC2	5.0	0.038	65		RC2
G13092	SCBR5.0-38-65-P-NS-RDC	5.0	0.038	65		RDC
G11995	SCBR5.0-38-65-P-NS-RIM	5.0	0.038	65		RIM
G13431	SCBR5.0-38-80-P-NS-VS	5.0	0.038	80		VS
G13432	SCBR5.0-38-80-P-NS-VS1	5.0	0.038	80		VS1
G13433	SCBR5.0-38-80-P-NS-VS2	5.0	0.038	80		VS2
G13434	SCBR5.0-38-80-P-NS-VS3	5.0	0.038	80		VS3
G12798	SCBR5.0-38-65-P-NS-LEV1	5.0	0.038	65		LEV1
G12156	SCBR5.0-38-65-P-NS-KMP	5.0	0.038	65		KMP
G28535	SCBR5.0-38-100-P-NS-KMP	5.0	0.038	100		
G11993	SCBR5.0-38-40-P-NS-KMP	5.0	0.038	40		
G23689	SCBR5.0-38-65-P-NS-MPA	5.0	0.038	65		MPA
G23692	SCBR5.0-38-80-P-NS-MPA	5.0	0.038	80		
G12002	SCBR5.0-38-100-P-NS-MPA	5.0	0.038	100		
G12179	SCBR5.0-38-65-P-NS-TEGT	5.0	0.038	65		TEGT
G23694	SCBR5.0-38-100-P-NS-TEGT	5.0	0.038	100		
G32367	SCBR5.0-38-125-P-NS-TEGT	5.0	0.038	125		
G36320	SCBR5.0-38-65-P-NS-TC	5.0	0.038	65		TC
G36322	SCBR5.0-38-65-P-NS-TC-BNK	5.0	0.038	65		TC-BNK
G43120	SCBR5.0-38-100-P-NS-RVC	5.0	0.038	100		RVC
G12009	SCBR5.0-38-65-P-NS-0	5.0	0.038	65		0
G12005	SCBR5.0-38-100-P-NS-0	5.0	0.038	100		
G32366	SCBR5.0-38-125-P-NS-0	5.0	0.038	125		

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Distribution for details.

The ultimate coaxial system

Pairing Slip-Cath Shuttle with Flexor®—the versatility you can rely on

Slip-Cath Selective Catheter

- **Visibility**

Under fluoroscopy, the easy-to-see tip lights the way through every turn in the body's anatomy.

- **Trackability**

The hydrophilic coating and tapered tip facilitate vessel selection, ranging from simple vasculature to more challenging anatomy. The increased lubricity allows the user to navigate the vasculature system smoothly.

- **Torqueability**

A nylon catheter construction with stainless steel braiding through the shaft provides control, pushability, and torque response.

Flexor Shuttle Select® Guiding Sheath

- **Kink resistance**

Patented coil-reinforced construction provides optimal flexibility and maximum resistance to kinking or compression.

- **Trackability**

Select configurations with hydrophilic coating to reduce friction and enhance trackability.

- **Low friction**

Large, low-friction, PTFE-coated lumen facilitates effortless delivery.

- **Visibility**

Radiopaque band incorporated within the sheath material identifies location of the distal tip for positioning accuracy.

Shuttle Select System







The Shuttle Select can be paired coaxially with a selective catheter to simplify access in challenging arch anatomy while minimizing device exchanges. Flexor Shuttle Select balances pushability and trackability.

Flexor Shuttle Select Guiding Sheath

Order Number	Reference Part Number (RPN)	Introducer Size Fr*	Introducer Lumen Diameter inch/mm	Accepts Wire Guide Diameter inch
G31133	KSAW-5.0-38-90-RB-SHTL-FLEX-HC	5.0	0.074/1.9	0.038
G31138	KSAW-6.0-38-90-RB-SHTL-FLEX-HC	6.0	0.087/2.2	0.038

*The maximum diameter of the instrument or catheter to be introduced should be determined to ensure its passage through the sheath. Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Distribution for details.

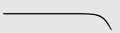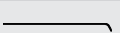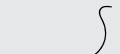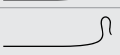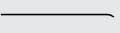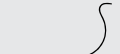
Slip-Cath Beacon Tip Hydrophilic Angiographic Catheter

Order Number	Reference Part Number (RPN)	Catheter Outer Diameter Fr	Accepts Wire Guide Diameter inch	Length cm	Tip Configuration
5.0 Fr Sheath Compatible					
G23697	SCBR5.5-35-125-P-NS-MPA-SHTL	5.5	0.035	125	 MPA
G23699	SCBR5.5-35-125-P-NS-SIM1-SHTL	5.5	0.035	125	 SIM1
G36360	SCBR5.5-35-125-P-NS-SIM2-SHTL	5.5	0.035	125	 SIM2
G23698	SCBR5.5-35-125-P-NS-VERT-SHTL	5.5	0.035	125	 VERT
G31129	SCBR5.5-35-125-P-NS-H1-SHTL	5.5	0.035	125	 H1
G31130	SCBR5.5-35-125-P-NS-JB1-SHTL	5.5	0.035	125	 JB1
G36361	SCBR5.5-35-125-P-NS-JB2-SHTL	5.5	0.035	125	 JB2
G36359	SCBR5.5-35-125-P-NS-VTK-SHTL	5.5	0.035	125	 VTK
G48152	SCBR5.5-35-75-P-NS-ANG-SHTL	5.5	0.035	75	 ANG
6.0 Fr Sheath Compatible					
G23701	SCBR6.5-35-125-P-NS-SIM1-SHTL	6.5	0.035	125	 SIM1
G36363	SCBR6.5-35-125-P-NS-SIM2-SHTL	6.5	0.035	125	 SIM2
G23700	SCBR6.5-35-125-P-NS-VERT-SHTL	6.5	0.035	125	 VERT
G31134	SCBR6.5-35-125-P-NS-H1-SHTL	6.5	0.035	125	 H1
G31135	SCBR6.5-35-125-P-NS-JB1-SHTL	6.5	0.035	125	 JB1
G36364	SCBR6.5-35-125-P-NS-JB2-SHTL	6.5	0.035	125	 JB2
G36362	SCBR6.5-35-125-P-NS-VTK-SHTL	6.5	0.035	125	 VTK


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Distribution for details.

Common catheter curves

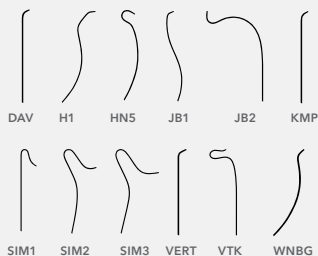CEREBRAL




VISCERAL/AORTIC

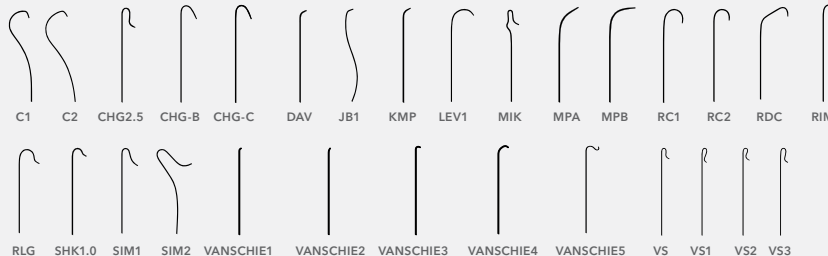


BIFURCATION/PELVIC



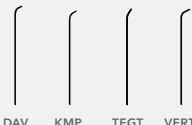MULTIPUPOSE



Vessel (average diameter)

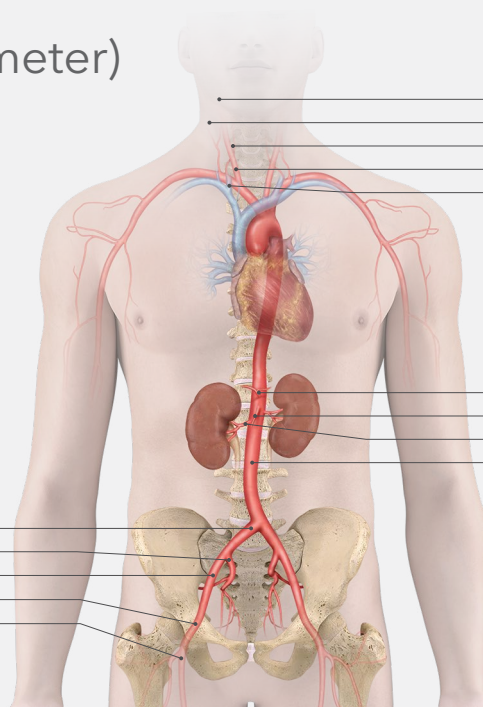

CEREBRAL

- External carotid (5–6 mm)
- Vertebral (3–5 mm)
- Internal carotid (5–6 mm)
- Common carotid (5–8 mm)
- Brachiocephalic (innominate) (8–10 mm)

VISCERAL/AORTIC

- Celiac trunk (5–7 mm)
- Superior Mesenteric Artery (6–8 mm)
- Renal (5–7 mm)
- Descending aorta (16–18 mm)

BIFURCATION/PELVIC

- Common iliac (7–10 mm)
- Internal iliac (5–6 mm)
- External iliac (6–8 mm)
- Common femoral (6–8 mm)
- Superficial femoral (4–6 mm)

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001702, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Iceland: +354 8007615, is.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: customersupport@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 E-mail: cau.custserv@cookmedical.com



AI, ESC, IR, OHNS, PI, RH, SUR-8.5X11