


Opportunità di accesso^{1,2}

Progettato per supportare la tecnologia dei dispositivi utilizzati per l'accesso e il crossing delle ostruzioni venose^{1,2,3}

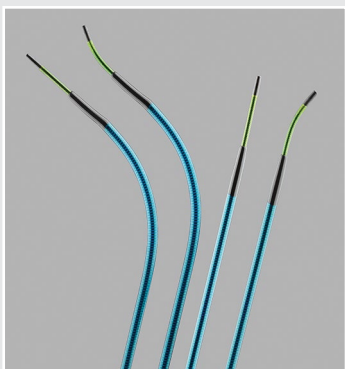

Quattro configurazioni e due lunghezze disponibili

TriForce[®]

PERIPHERAL CROSSING SET

Specifico per l'anatomia e la patologia


Progettato per supportare la tecnologia dei dispositivi utilizzati per l'accesso e il crossing delle ostruzioni venose.^{1,2,3}

Forza di torsione reattiva

Unisce la forza colonnare intrinseca dell'introduttore Flexor[®] e l'intreccio in acciaio inossidabile del catetere per crossing CXI[®], consentendo l'avanzamento continuo.

Capacità di spinta e trackability avanzate

Punta con rivestimento idrofilo e carica di tungsteno per migliorare la capacità di crossing.

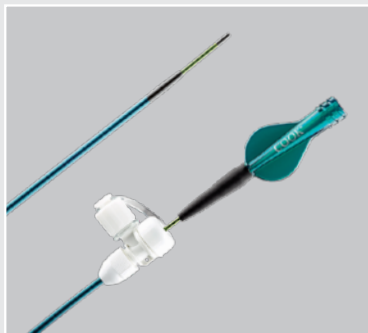

TriForce®

PERIPHERAL CROSSING SET

- Il set per crossing periferico TriForce include una guaina esterna Flexor da 5,0 Fr, un catetere di supporto CXI da 4,0 Fr e una guaina Peel-Away®.
- La guaina esterna Flexor è rastremata fino a un diametro interno di 4,0 Fr in corrispondenza della punta, per garantire una transizione fluida fino al catetere di supporto CXI.
- La guaina esterna Flexor è caratterizzata da punta radiopaca e rivestimento idrofilo.
- TriForce È compatibile con guide da 0,035 inch di diametro.

Codice di ordinazione	Codice prodotto	Diametro Fr/lunghezza del catetere cm	Lunghezza dell'introduttore cm	Punta dell'introduttore Fr	Diametro delle guide compatibili inch
Set con introduttore a punta diritta e catetere a punta diritta					
G56412	KCXS-5.0-35-65-RB-0/0-HC	4,0/65	55	4,0	0,035
G56416	KCXS-5.0-35-100-RB-0/0-HC	4,0/100	90	4,0	0,035
Set con introduttore a punta diritta e catetere a punta curva					
G56413	KCXS-5.0-35-65-RB-0/DAV-HC	4,0/65	55	4,0	0,035
G56417	KCXS-5.0-35-100-RB-0/DAV-HC	4,0/100	90	4,0	0,035
Set con introduttore a punta curva e catetere a punta diritta					
G56414	KCXS-5.0-35-65-RB-MPB/0-HC	4,0/65	55	4,0	0,035
G56418	KCXS-5.0-35-100-RB-MPB/0-HC	4,0/100	90	4,0	0,035
Set con introduttore a punta curva e catetere a punta curva					
G56415	KCXS-5.0-35-65-RB-MPB/DAV-HC	4,0/65	55	4,0	0,035
G56419	KCXS-5.0-35-100-RB-MPB/DAV-HC	4,0/100	90	4,0	0,035

Alcuni prodotti o codici prodotto potrebbero non essere disponibili in tutti i mercati. Per maggiori informazioni, rivolgersi al rappresentante Cook di zona o al team Customer Support & Distribution.

Consultare le informazioni sui rischi del prodotto nelle Istruzioni per l'uso (IFU) all'indirizzo cookmedical.eu.

- Hofmann L. What's new and what's needed in chronic venous occlusion therapy. Sito web di Endovascular Today. <https://evtoday.com/articles/2015-oct/whats-new-and-whats-needed-in-chronic-venous-occlusion-therapy>. Data di pubblicazione: 2015. Data di consultazione: 20 dicembre 2021.
- Caso clinico - Erbacher G. Case study: Crossing a difficult venous occlusion with the TriForce® Peripheral Crossing Set, Oklahoma State University Medical Center, PI-D64062-EN. https://www.cookmedical.com/data/resources/PI-D52143-EN-F_M3_1576703674956.pdf
- Lu T, Schwein A, Bechara C, Bismuth J. Improved percutaneous Central venous recanalization: rates with novel telescoping catheter system. JVS 2016;63(6) (suppl):132S-133S. [https://www.jvascsurg.org/article/S0741-5214\(16\)00539-5/fulltext](https://www.jvascsurg.org/article/S0741-5214(16)00539-5/fulltext)

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001702, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Iceland: +354 800 7615, is.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: customersupport@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283
 E-mail: cau.custserv@cookmedical.com



AI-ESC-IR-OHNS-PI-RH-SUR-A4