

Simplifying your source for complex crossings

BeBack
crossing catheter



The **BeBack crossing catheter** is a low-profile catheter that has complex crossing and targeted re-entry capabilities. It features a length-adjustable, 360° nitinol needle. The special radiopaque marker indicates in which direction the curved needle protrudes from the center of the catheter.




Available in 2.9 F and 4 F
for antegrade, retrograde, and
up-and-over approaches



The effective low-profile
catheter for complex crossing
and re-entry



A directional needle
at the tip of the catheter
has adjustable lengths.

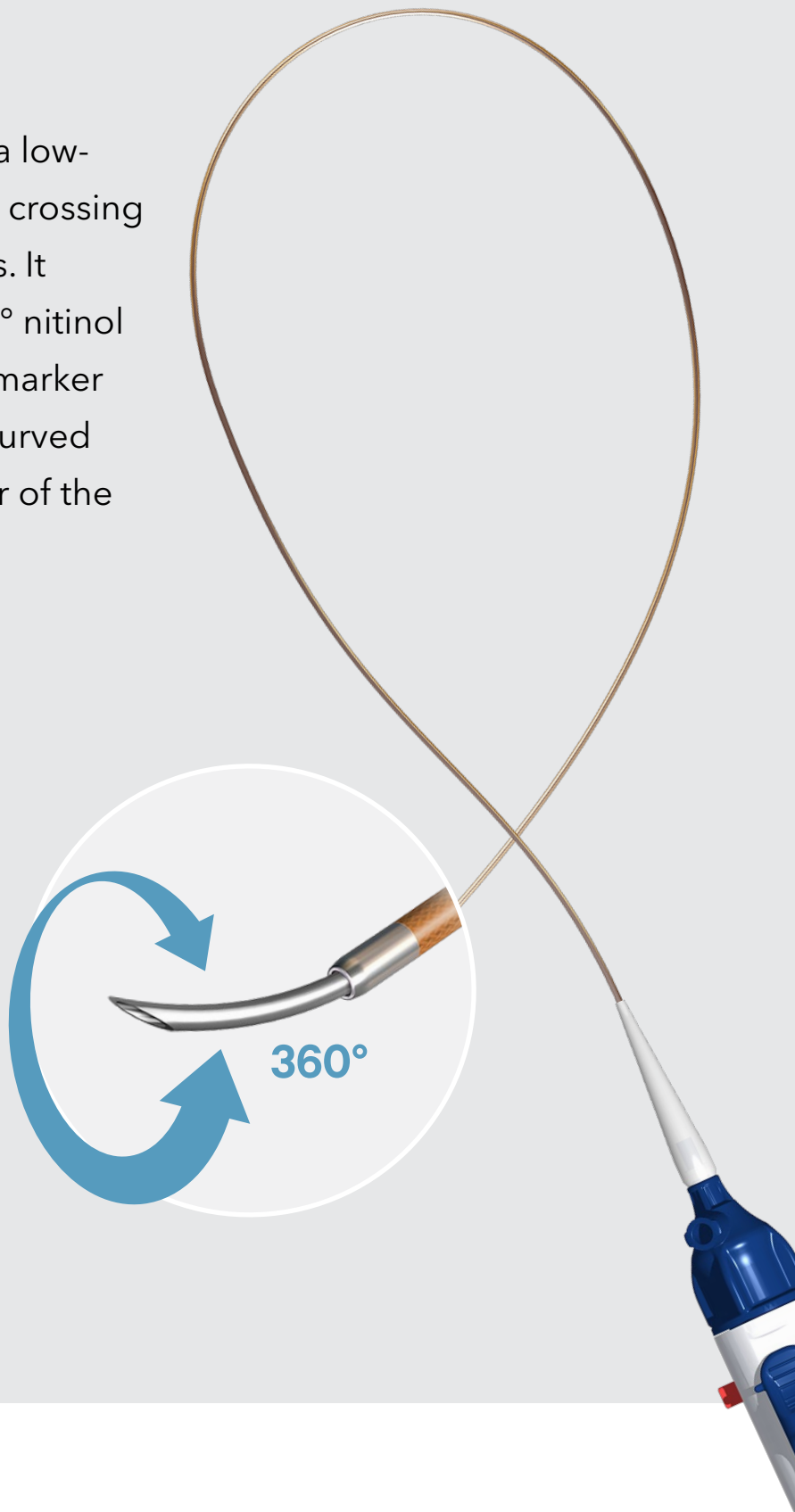

Cross chronic total occlusions (CTOs) with success


92.3%
technical success*¹

92.3%
procedural success*¹


The BeBack catheter's "adept ability to **navigate and recanalize** occluded segments provides a robust alternative, especially when traditional wire guide and catheter techniques falter."¹



Chronic total occlusion crossing

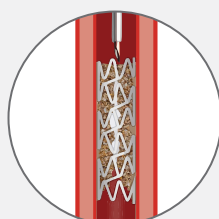

In-stent restenosis

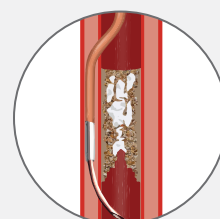


Subintimal crossing and re-entry

*In 72/78 CTO procedures where traditional wire guide and catheter techniques had failed.

Technical success was defined as "CTO recanalization by occlusion crossing or luminal re-entry via the application of BeBack and demonstration of arterial patency with completion angiogram."¹ **Procedural success** was defined as "arterial recanalization with < 30% residual stenosis and improvement in ABI by at least 0.2 after 24 hours."¹

1. Cohen I, Raskin D, Yiftach B, et al. Performance of BeBack Crossing Catheter in chronic total occlusion in patients with chronic limb-threatening ischemia: a single-center experience. J Endovasc Ther. 2023. Doi: 10.1177/15266028231204264 (E-pub. ahead of print)

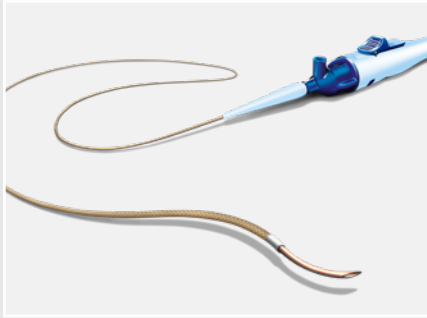

Indications

The purpose and intended use of the BeBack crossing catheter is to be used in conjunction with a steerable wire guide to access discrete regions of the peripheral vasculature. It may be used to facilitate placement and exchange of wire guides.

Product specifications

	BeBack crossing catheter 2.9 F	BeBack crossing catheter 4 F
Introducer sheath compatibility	2.9 F 5 F for up-and-over approach	4 F 6 F for up-and-over approach
Wire guide	0.014 inch	0.018 inch
Catheter length	80 cm / 120 cm	80 cm / 120 cm
Material needle	Nickel-titanium alloy (nitinol)	Nickel-titanium alloy (nitinol)
Profile needle	23.8 gage	21.5 gage
Rotation	360°	360°
Marker	Radiopaque C-shaped marker; the opening of the C indicates needle direction	

Ordering information

Cook Order Number	Cook RPN	Previous Bentley Product #	Size Fr	Length cm	Wire Guide Diameter inch
G60417	GB-600-US-018P-80L	GB 600-US 018P 80L	4	80	0.018
G60418	GB-600-US-018P-120L	GB 600-US 018P 120L	4	120	0.018
G60419	GB-603-US-014P-80L	GB 603-US 014P 80L	2.9	80	0.014
G60420	GB-603-US-014P-120L	GB 603-US 014P 120L	2.9	120	0.014

No information provided in this brochure is intended to be professional medical advice. The Instructions for Use (IFU) must be consulted before use of any product. Before purchasing this product, please review whether the product and the indication are approved in your country.

Please contact your Cook Medical representative for more information.

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001702, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Iceland: +354 8007615, IS.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: customersupport@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283

E-mail: cau.custserv@cookmedical.com



AI, ESC, IR, OHNS, PI, RH, SUR-8.5X11