


Challenge lesions with focal force

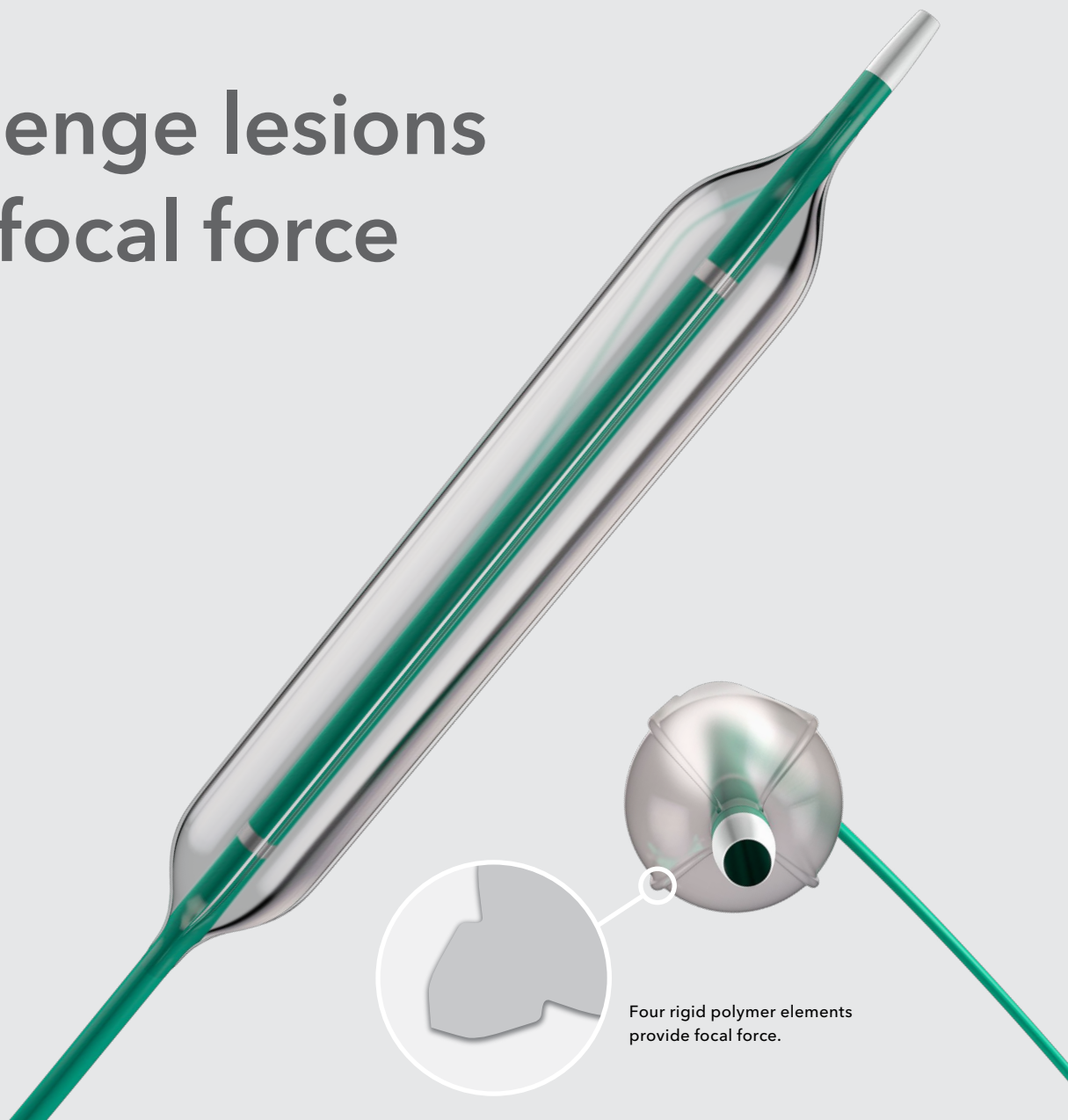

Advance® Enforcer® 35

FOCAL-FORCE PTA BALLOON CATHETER

Focal force


The balloon construction uniquely incorporates four elements that transfer focal force to the lesion. You can now dilate lesions using less pressure.

Enhanced flexibility

The single extrusion of the Advance Enforcer 35—with no attached components—enhances flexibility.

Performance

Enforcer 35 facilitates procedural success.¹

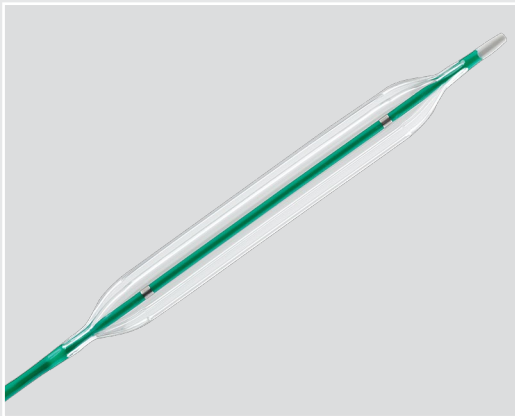

Advance® Enforcer® 35

FOCAL-FORCE PTA BALLOON CATHETER

Intended for percutaneous transluminal angioplasty of lesions in peripheral arteries including iliac, renal, popliteal, infrapopliteal, femoral, and iliofemoral, as well as obstructive lesions of native or synthetic arteriovenous dialysis fistulas.



Single extrusion element

Order Number	Reference Part Number	Accepts Wire Guide Diameter inch	Inflated Balloon Diameter mm	Inflated Length cm	Compatible Sheath Fr	Nominal Pressure atm	Rated Burst Pressure atm
80 cm Catheter Length							
G35252	ASB5-35-80-6-4	0.035	6	4	6.0	10	16
G35253	ASB5-35-80-8-4	0.035	8	4	7.0	10	12
G35254	ASB5-35-80-10-4	0.035	10	4	7.0	8	11
G35255	ASB5-35-80-12-4	0.035	12	4	8.0	6	10
135 cm Catheter Length							
G35257	ASB5-35-135-6-4	0.035	6	4	6.0	10	16
G35258	ASB5-35-135-8-4	0.035	8	4	7.0	10	12
G35259	ASB5-35-135-10-4	0.035	10	4	7.0	8	11
G35261	ASB5-35-135-12-4	0.035	12	4	8.0	6	10

Please follow any step-by-step instructions carefully, including the risk and technical information provided in the IFU. Note that some products or part numbers may not be available in all markets; consult your Cook representative for details.

¹ Holden A, Hill AA, Buckley B, et al. Clinical evaluation of the Cook Advance Enforcer 35 Focal-Force PTA Balloon Catheter for treatment of hemodialysis fistula stenoses: a feasibility study. *J Vasc Interv Radiol.* 2019;30(1):61-68. <https://doi.org/10.1016/j.jvir.2018.08.015>.

Customer Service

EU Website: cookmedical.eu
 EDI: cookmedical.eu/edi
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001633, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: orders@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283

E-mail: cau.custserv@cookmedical.com-A4



AI, ESC, IR, OHNS, PI, RH, SUR-A4