

# Enhance accuracy in oocyte collection.

RETRIEVE

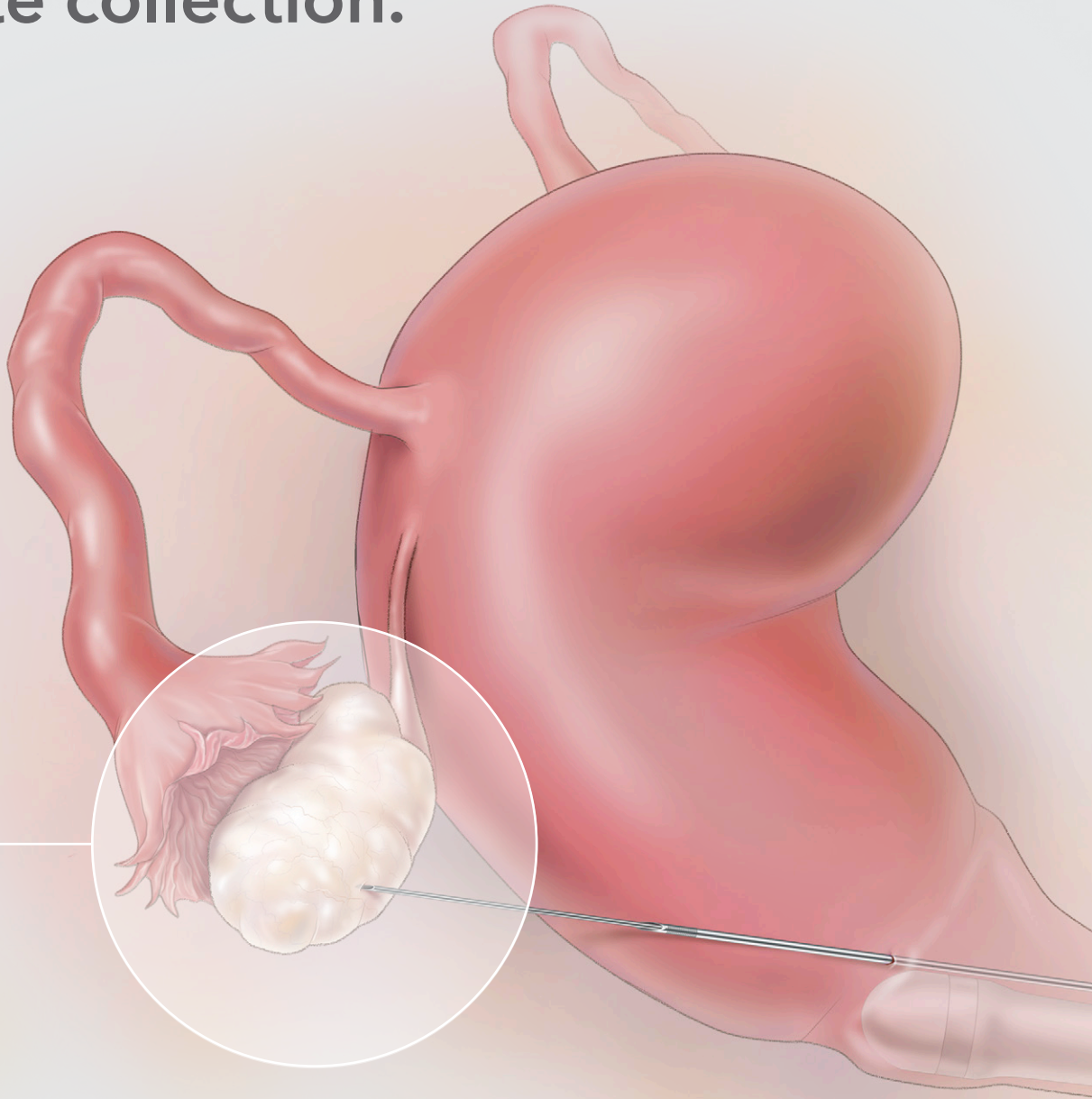


Illustration by Lisa Clark

## SMALL GAGE OVUM ASPIRATION NEEDLES




[www.cookmedical.com](http://www.cookmedical.com)

## SMALL GAGE OVUM ASPIRATION NEEDLES


**Collecting immature oocytes requires precision and control.**

The Small Gage Ova-Stiff™ Ovum Aspiration Needle from Cook Medical offers a finely tuned bevel tip specifically designed for effectiveness in small antral follicles. Highly visible under ultrasonic imaging, this tip enables accuracy, and the stiff cannula gives the clinician further control of placement. As a result, more viable oocytes can be collected.

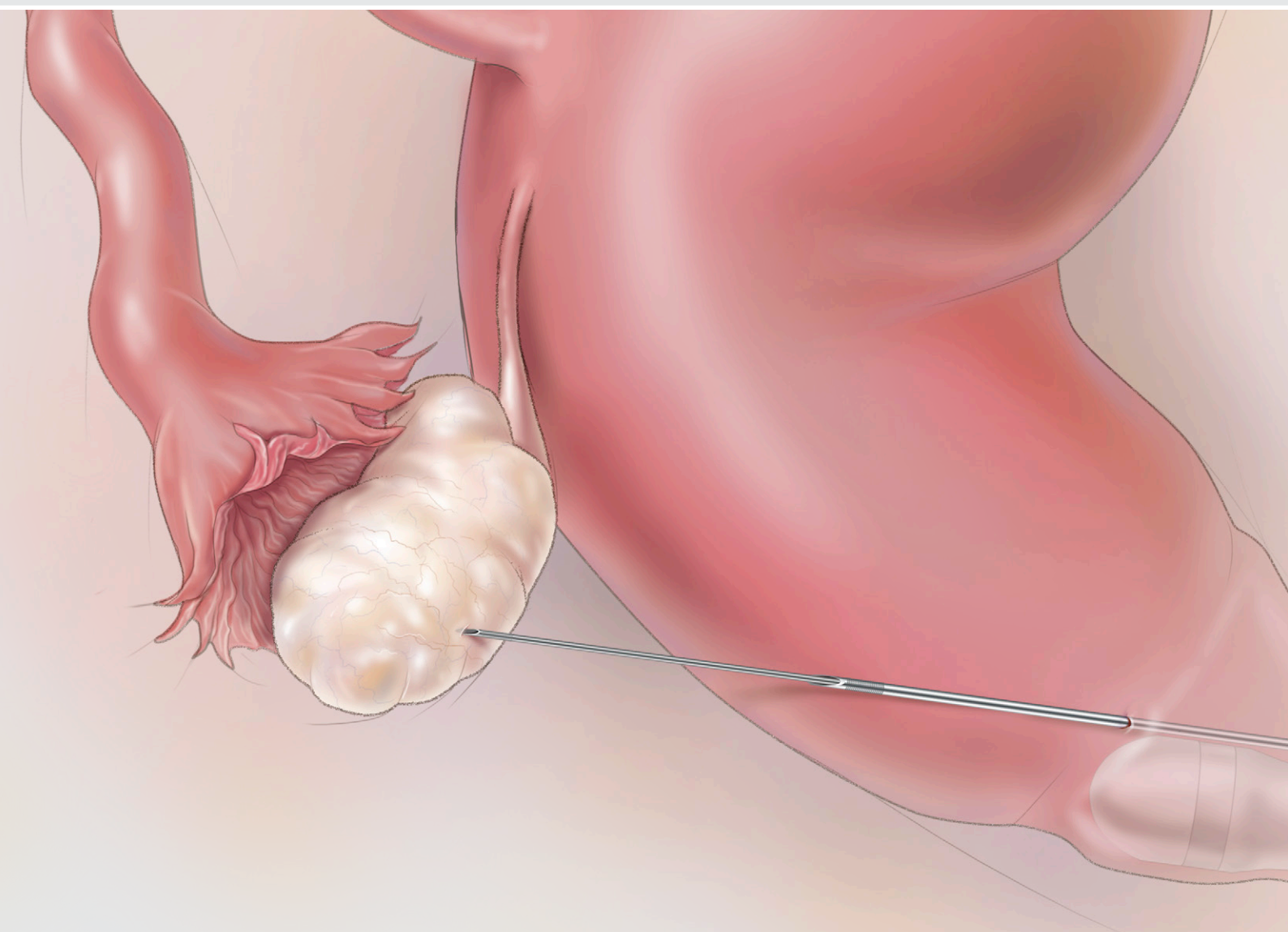


## SMALL GAGE OVUM ASPIRATION NEEDLES

### Features and benefits of ovum aspiration needles.

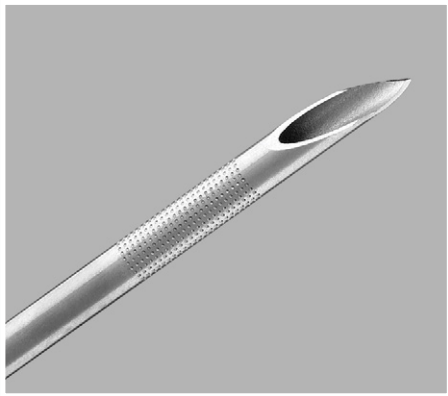

#### Small Gage Ova-Stiff™

- Extra-stiff cannula for precise placement of tip in follicle.
- Standard RWH noncoring bevel minimizes tissue entering lumen.
- Translucent FEP vacuum line. Flushing is achieved by attaching the flushing cannula to a syringe and inserting the cannula into the flared end of the aspiration line at the silicone stopper.
- EchoTip® enhances visualization of needle tip when used with ultrasonic imaging equipment.

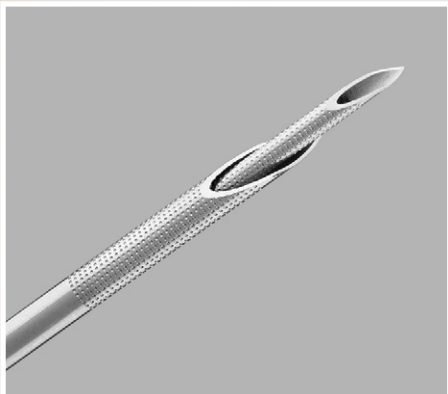

#### Immature Ovum Aspiration Needle

- EchoTip extends around the bevel, improving visibility and placement of the needle.
- 17 gage, 30 cm guide needle used to puncture vaginal wall, adding rigidity to the 20 gage, 35 cm aspiration needle.
- New stopper design accepts the tip of a Luer slip syringe, facilitating follicle flushing.



## Small Gage Ova-Stiff™ OVUM ASPIRATION NEEDLE

Order Number	Reference Part Number	gage	Length cm	Aspiration Line cm	Bevel	Handle Color
G20943	K-OPS-6035-RWH-B-ET	18	35	60	B	pink
G19715	K-OPS-6035-RWH-ET	18	35	60	RWH	pink
G26182	K-OPS-7030-RWH-ET	19	30	60	RWH	brown
G26638	K-OPS-7035-RWH-ET	19	35	60	RWH	brown
G38845	K-OPS-7035-RWH-ET-US	19	35	60	RWH	brown
G26631	K-OPS-8035-RWH-ET	20	35	60	RWH	yellow
G26056	K-OPS-8035-RWH-B-ET	20	35	60	B	yellow
G44359	K-OPS-8035-RWH-B-ET-US	20	35	60	RWH	yellow
G50735	K-OPS-2130	21	30	60	A	green
G44360	K-OPS-2130-US	21	30	60	A	green
G50736	K-OPS-2135	21	35	60	A	green

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.



## IMMATURE OVUM ASPIRATION NEEDLE

Order Number	Reference Part Number	Guide Needle gage	Guide Needle Length cm	Aspiration Needle gage	Aspiration Needle Length cm	Aspiration Tubing Length cm
G38609	K-IOPS-2035-1730	17	30	20	35	60
G44358	K-IOPS-2035-1730-US	17	30	20	35	60

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

### Customer Service

**EMEA:** EDI - [www.cookmedical.com/edi.do](http://www.cookmedical.com/edi.do)  
 Distributors: +353 61239240, [ssc.distributors@cookmedical.com](mailto:ssc.distributors@cookmedical.com)  
 Austria: +43 179567121, [oe.orders@cookmedical.com](mailto:oe.orders@cookmedical.com)  
 Belgium: +32 27001633, [be.orders@cookmedical.com](mailto:be.orders@cookmedical.com)  
 Denmark: +45 38487607, [da.orders@cookmedical.com](mailto:da.orders@cookmedical.com)  
 Finland: +358 972519996, [fi.orders@cookmedical.com](mailto:fi.orders@cookmedical.com)  
 France: +33 171230269, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
 Germany: +49 6950072804, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
 Hungary: +36 17779199, [hu.orders@cookmedical.com](mailto:hu.orders@cookmedical.com)  
 Ireland: +353 61239252, [ie.orders@cookmedical.com](mailto:ie.orders@cookmedical.com)  
 Italy: +39 0269682853, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
 Netherlands: +31 202013367, [nl.orders@cookmedical.com](mailto:nl.orders@cookmedical.com)  
 Norway: +47 23162968, [no.orders@cookmedical.com](mailto:no.orders@cookmedical.com)  
 Spain: +34 912702691, [es.orders@cookmedical.com](mailto:es.orders@cookmedical.com)  
 Sweden: +46 858769468, [se.orders@cookmedical.com](mailto:se.orders@cookmedical.com)  
 Switzerland - French: +41 448009609, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
 Switzerland - Italian: +41 448009609, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
 Switzerland - German: +41 448009609, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
 United Kingdom: +44 2073654183, [uk.orders@cookmedical.com](mailto:uk.orders@cookmedical.com)

[www.cookmedical.com](http://www.cookmedical.com)

**Americas:** EDI - [www.cookmedical.com/edi.do](http://www.cookmedical.com/edi.do)  
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335  
 E-mail: [orders@cookmedical.com](mailto:orders@cookmedical.com)

**Australia:**  
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283  
 E-mail: [cau.custserv@cookmedical.com](mailto:cau.custserv@cookmedical.com)

