



# Embryo transfer catheters




< Touch to  
learn more

# Embryo transfer catheters




**Straight**

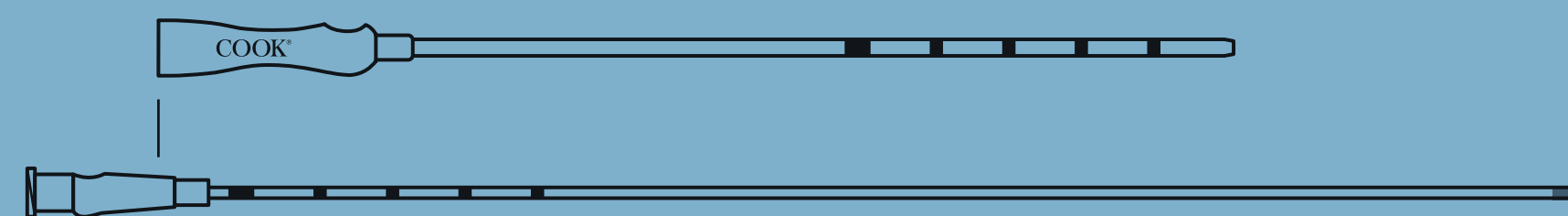

**Precurved with bulb tip**




**Specialty**

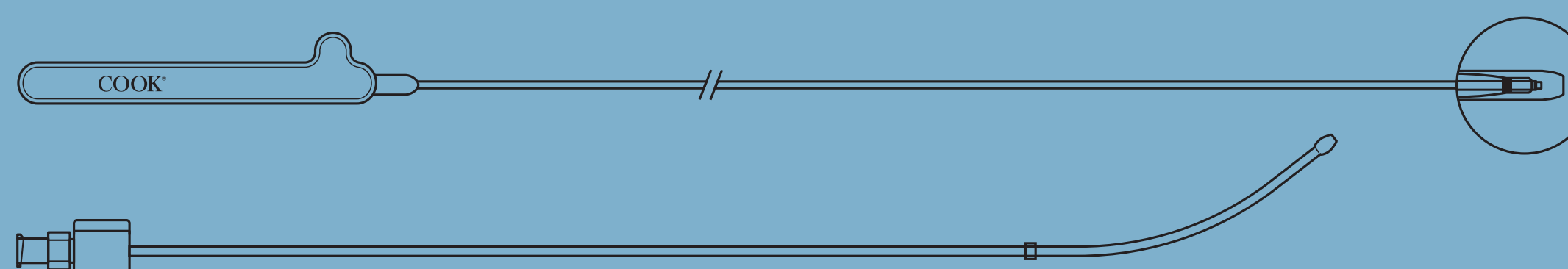

**Obturator**



# Embryo transfer catheters

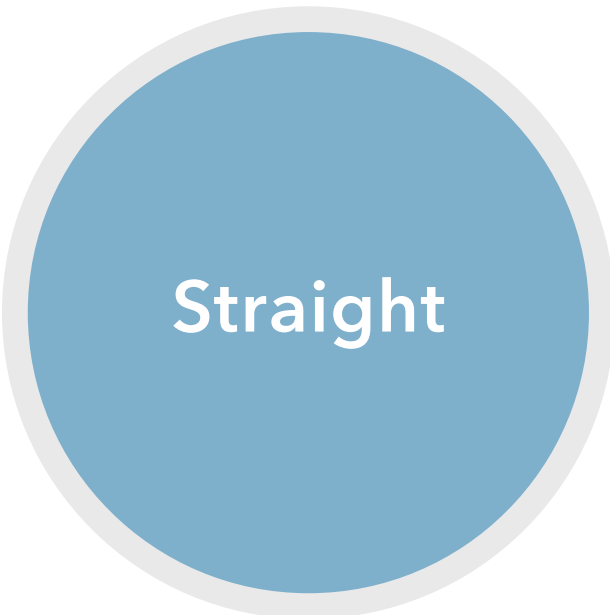

## Flexible guide catheter

### Soft-Pass™

GC = 6.8 Fr, 12 or 17 cm  
TC = 4.4 Fr, 19.7 or 24.7 cm

### Soft-Pass™ Embryo Transfer Catheter Set

With support cannula  
K-J-SPPE-681710  
Without support cannula  
K-J-SPPE-681200  
K-J-SPPE-681700




### Soft-Pass™ Embryo Transfer Catheter Set with EchoTip®

With support cannula  
K-J-SPPE-681210-ET  
K-J-SPPE-681710-ET  
Without support cannula  
K-J-SPPE-681200-ET  
K-J-SPPE-681700-ET



### Soft-Pass™ Mock Transfer Catheter Set

J-SPPE-681200-MC  
J-SPPE-681700-MC  
J-SPPE-681710-ET-MC

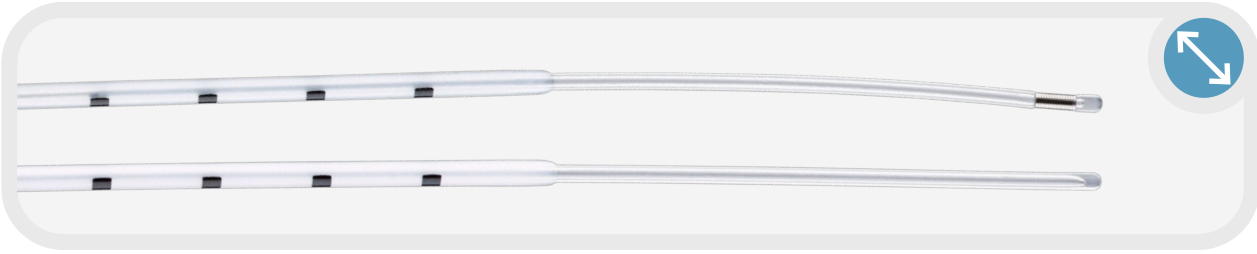

## Semirigid guide catheter

### Soft-Pass™

GC = 6.8 Fr, 12 or 17 cm  
TC = 4.4 Fr, 19.7 or 24.7 cm


### Soft-Pass™ Embryo Transfer Catheter

With support cannula  
K-J-SP-681710  
Without support cannula  
K-J-SP-681700



### Soft-Pass™ Embryo Transfer Catheter With EchoTip®

K-J-SP-681700-ET

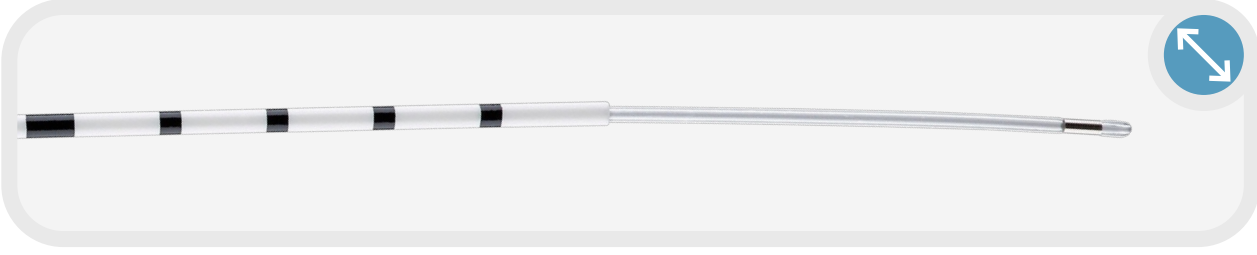

## Flexible guide catheter

### Pivot Embryo Transfer Catheter Set


K-PETS-5000  
GC = 7.0 Fr, 17.3 cm  
TC = 2.8 Fr, 24 cm



GC = Guide catheter, TC = Transfer catheter


Some products or part numbers may not be available in all markets. Speak with our Cook representative for details. Please see product risk information in the IFU at [cookmedical.eu](http://cookmedical.eu).

# Embryo transfer catheters

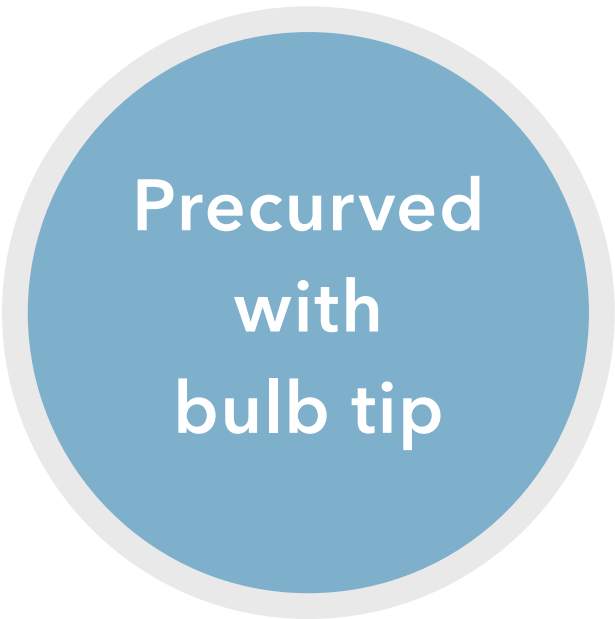


## Guardia™ Access Embryo Transfer Catheter

K-JETS-6019  
K-JETS-7019  
GC = 6.6 Fr, 16.7 or 17.3 cm  
TC = 2.8 Fr, 23 or 24 cm

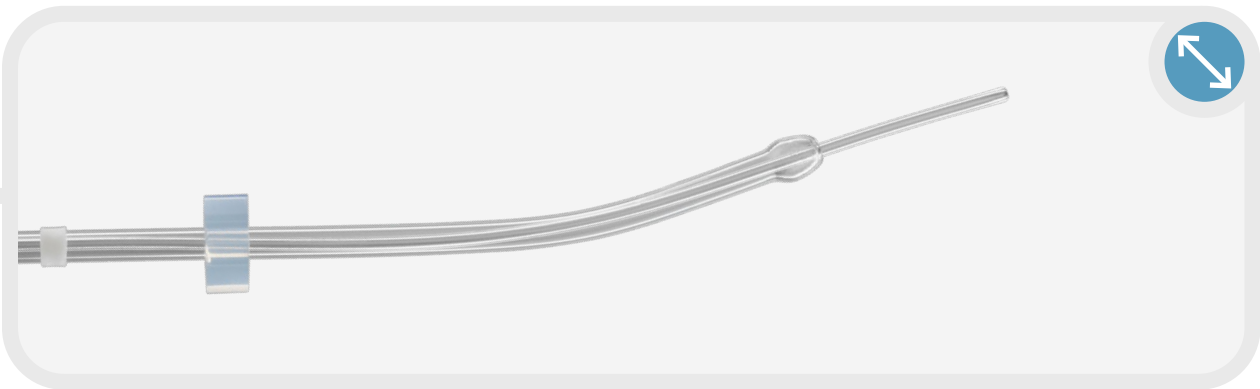

## Guardia™ AccessET Curved Embryo Transfer Catheter

K-JETS-6019-ET  
K-JETS-7019-ET  
GC = 6.6 Fr, 16.7 or 17.3 cm  
TC = 2.8 Fr, 23 or 24 cm

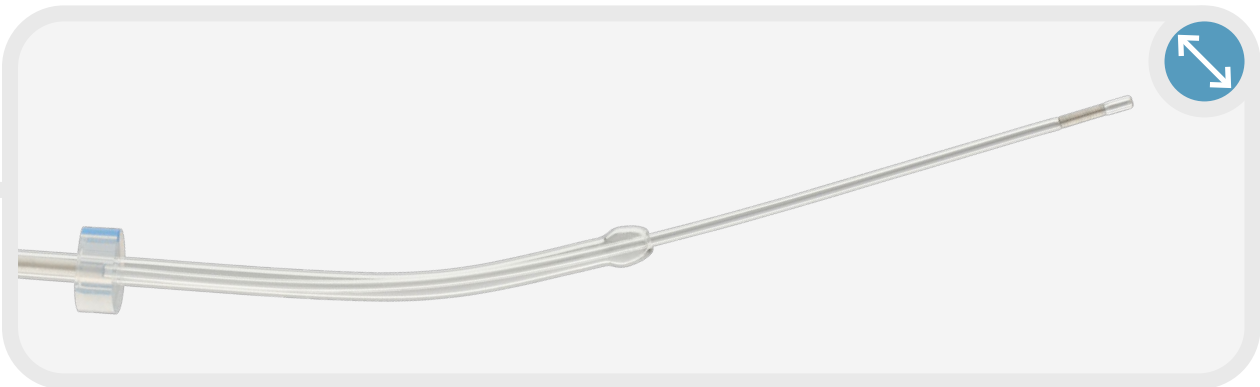

## Guardia™ Access Embryo Transfer Catheter with Internal Support Cannula

K-JETS-7019-INT  
GC = 6.6 Fr, 17.3 cm  
TC = 2.8 Fr, 25 cm



## Guardia™ Access Nano Embryo Transfer Catheter

K-JETS-551910-S  
GC = 5.5 Fr, 17.3 cm  
TC = 2.8 Fr, 24 cm



GC = Guide catheter, TC = Transfer catheter

Some products or part numbers may not be available in all markets. Speak with our Cook representative for details. Please see product risk information in the IFU at [cookmedical.eu](http://cookmedical.eu).





**Echosight® Jansen-Anderson  
Intratubal Transfer Set**

K-J-JITS-572900  
GC = 5.7 Fr, 29 cm  
TC = 3 Fr, 40 cm




# Embryo transfer catheters




## Malleable

J-UOB-2828  
2.8 Fr, 28 cm  
J-UOB-4028  
4 Fr, 28 cm

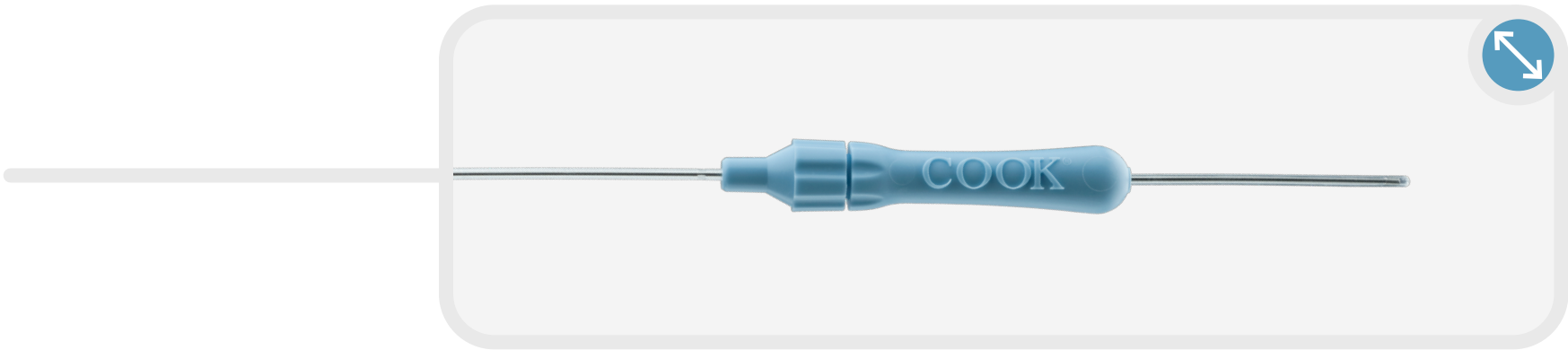


## Stiff

J-UOB-4028-ST  
4 Fr, 28 cm

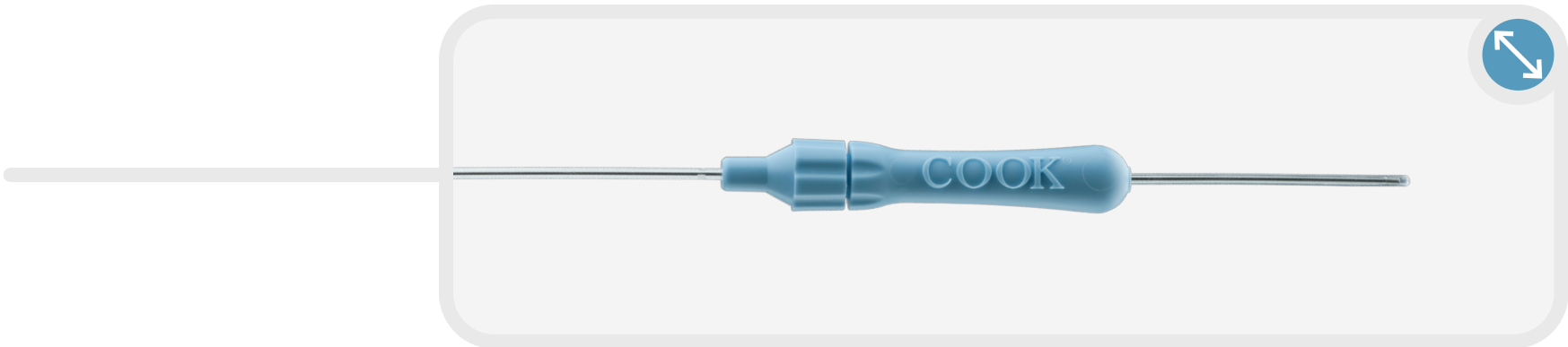


## Guardia™ Access Malleable

K-JET-7002  
4 Fr, 19 cm  
For use with K-JETS-7019

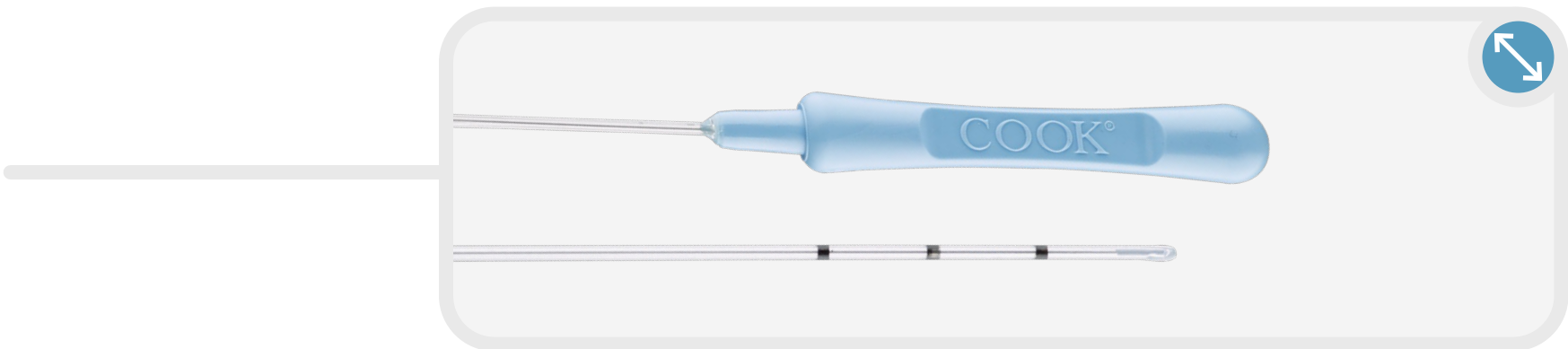

## Soft-Pass™ Malleable


J-SP-4420  
4.4 Fr, 20.1 cm

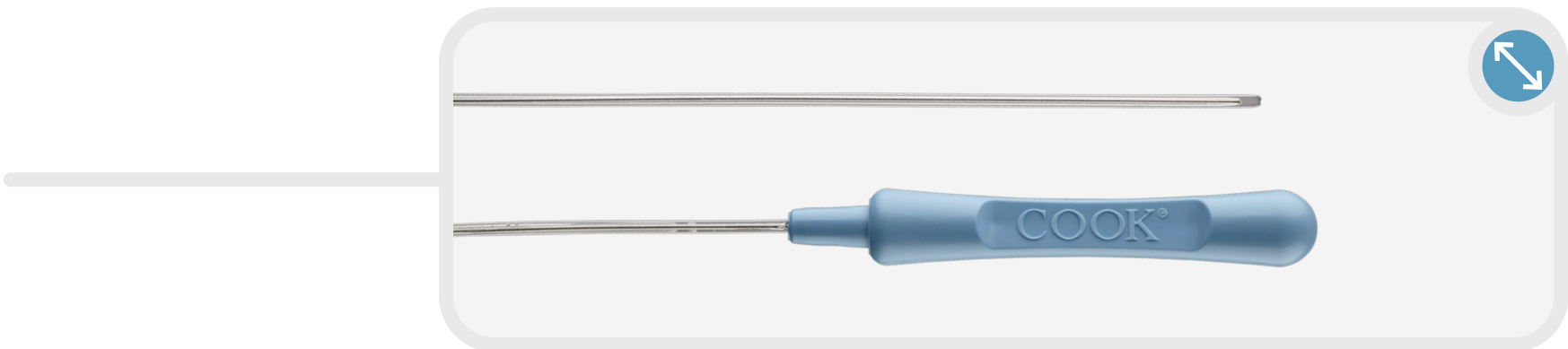

Some products or part numbers may not be available in all markets. Speak with our Cook representative for details.  
Please see product risk information in the IFU at [cookmedical.eu](http://cookmedical.eu).



## Pivot Embryo Transfer Catheter Set

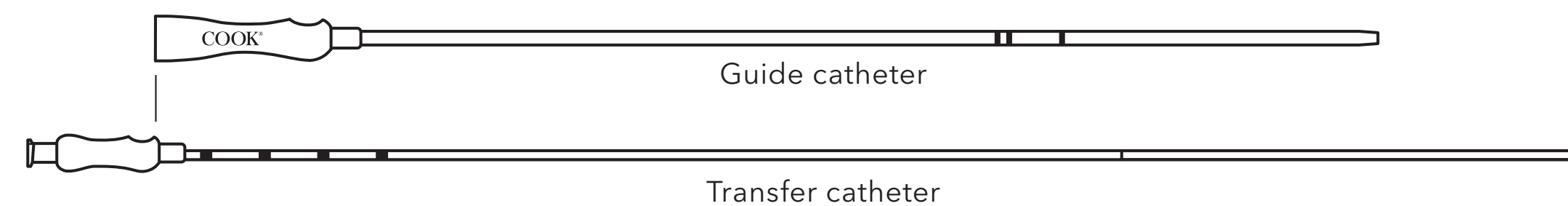


Straight

- The soft, flexible guide catheter is designed to follow the curves in the cervix in order to facilitate placement of the transfer catheter.
- The depth markings on the guide catheter aid in correct catheter placement.
- The soft 2.8 Fr transfer catheter features Microvol™ technology, which reduces the volume of medium that is required for embryo transfer.
- The integrated support cannula aids the placement of the transfer catheter.



## Soft-Pass™ Embryo Transfer Catheter Set

Straight

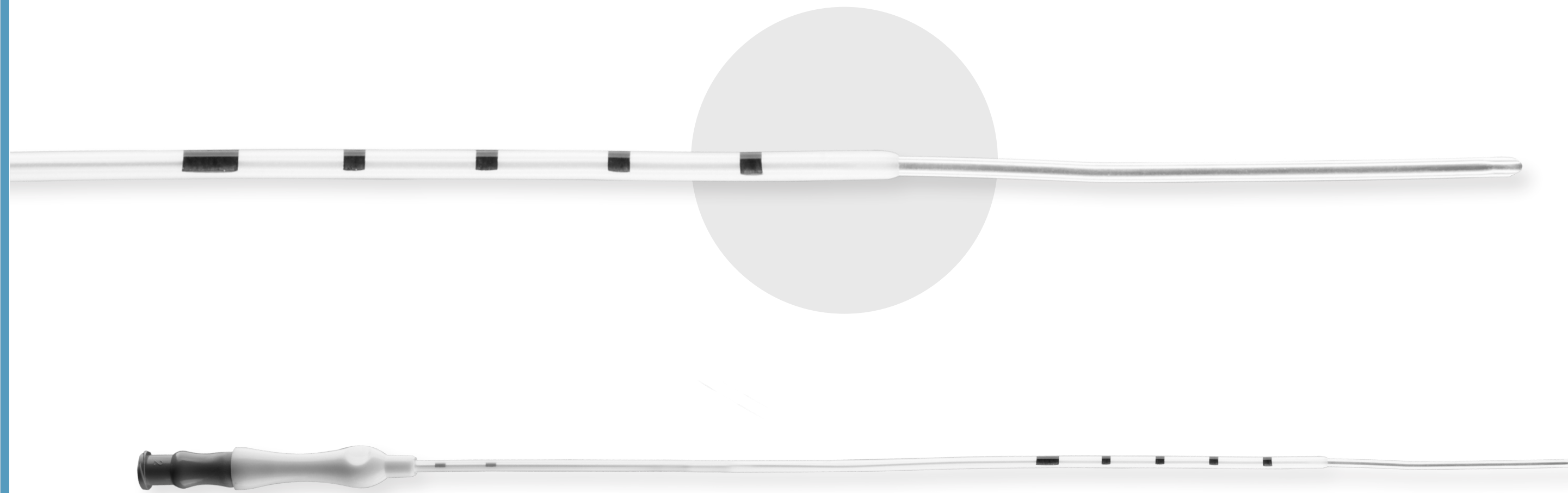

- The soft, flexible guide catheter is designed to facilitate placement of the set.
- Sets are available with transfer catheters that have an integrated stainless steel support cannula to aid catheter placement.





## Soft-Pass™ Embryo Transfer Catheter Set with EchoTip®

Straight




- The soft, flexible guide catheter is designed to facilitate placement of the set.
- Sets are available with transfer catheters that have an integrated stainless steel support cannula to aid catheter placement.
- The transfer catheter's EchoTip echogenic tip improves visibility under ultrasound.



## Soft-Pass™ Mock Transfer Catheter Set



Straight

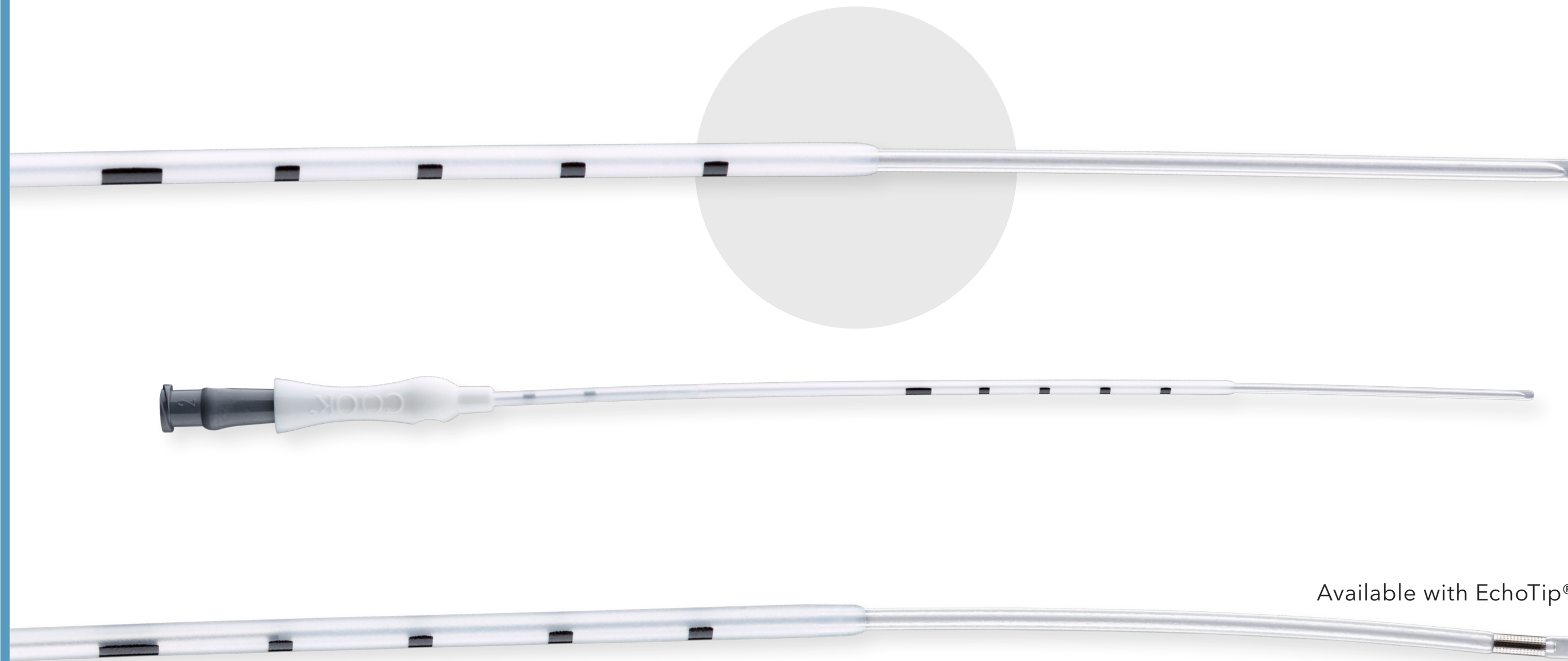


- The soft, flexible guide catheter is designed to facilitate placement of the set.
- Sets are available with transfer catheters that have an integrated stainless steel support cannula to aid catheter placement.
- The transfer catheter has a closed tip.



## Soft-Pass™ Embryo Transfer Catheter

Straight

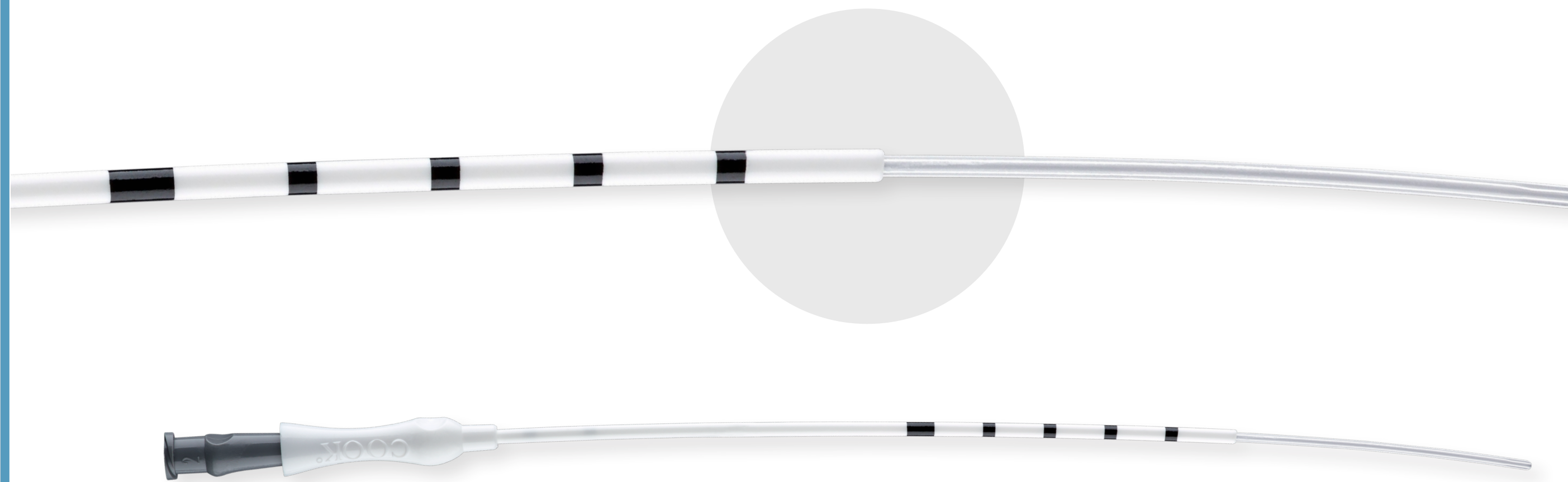

Used to place in vitro fertilized (IVF) embryos into the uterine cavity.



## Soft-Pass™ Embryo Transfer Catheter with EchoTip®

Straight



Used to place in vitro fertilized (IVF) embryos into the uterine cavity.

- The transfer catheter's echogenic tip improves visibility under ultrasound.






## Guardia™ Access Embryo Transfer Catheter



Precurved  
with  
bulb tip

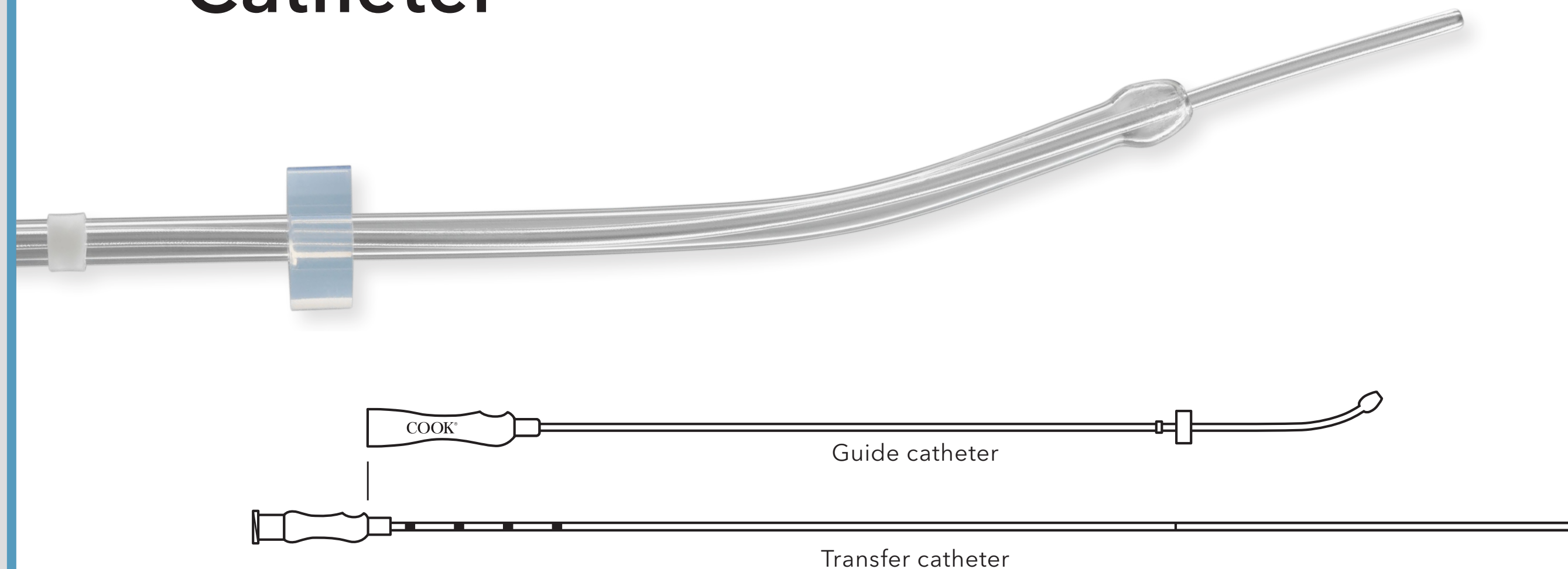


- The guide catheter is curved to facilitate insertion.
- The guide catheter's rounded bulb tip eases passage through the cervix.
- The transfer catheter is soft and flexible.
- The Microvol™ technology minimizes the volume of media that is required.
- The cervical stop is set at 4 cm and can be adjusted to 5 cm.
- The catheters are packaged in separate sterile pouches.



## Guardia™ AccessET Curved Embryo Transfer Catheter



Precurved  
with  
bulb tip

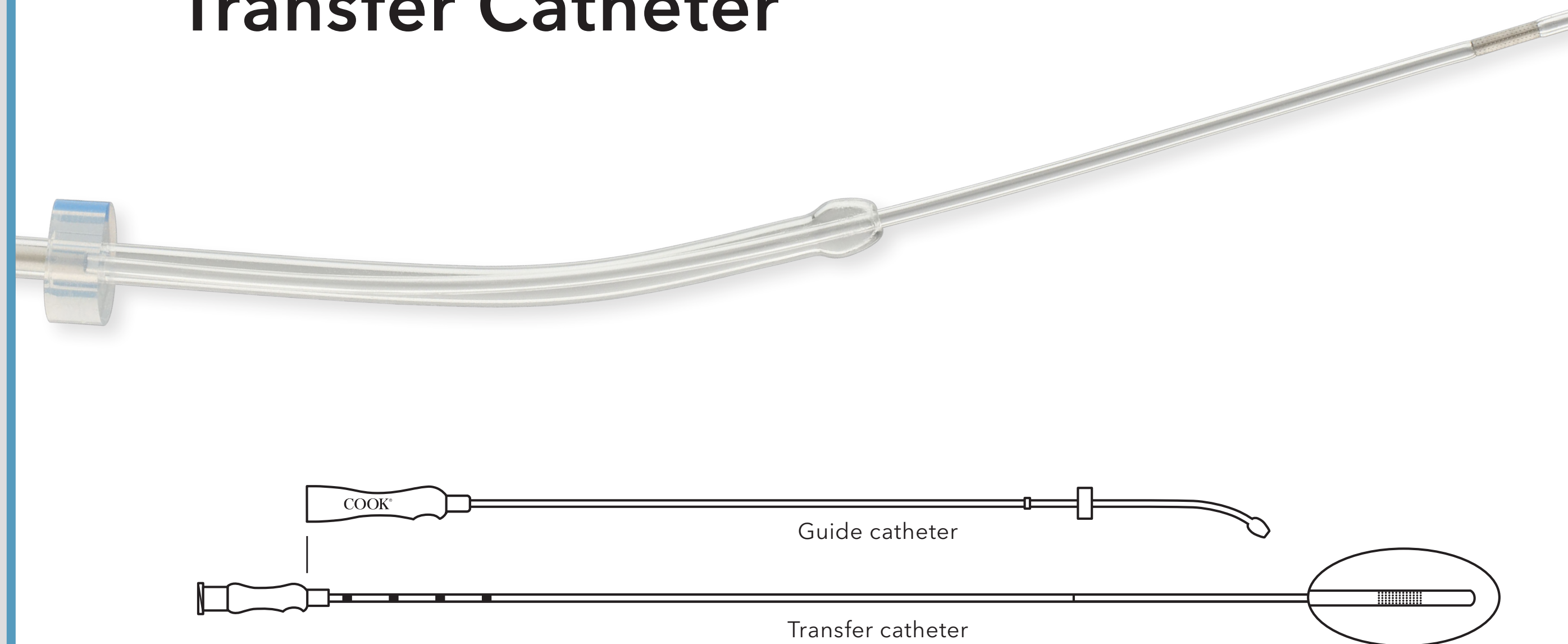


- The transfer catheter's echogenic tip improves visibility under ultrasound.
- The guide catheter is curved to facilitate insertion.
- The guide catheter has a rounded bulb tip to ease passage through the cervix.



## Guardia™ Access Embryo Transfer Catheter with Internal Support Cannula



Precurved  
with  
bulb tip

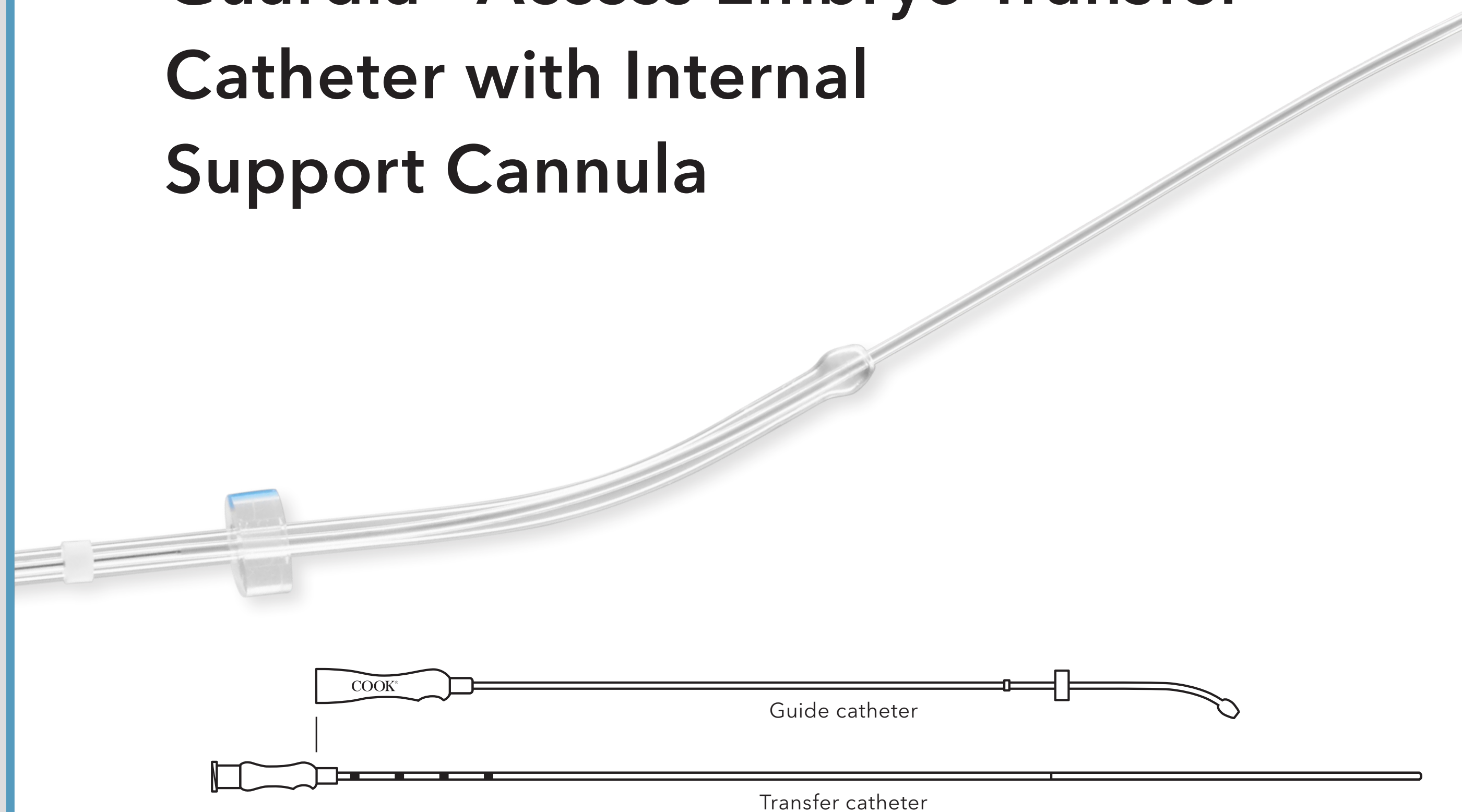

Used to place in vitro fertilized (IVF) embryos into the uterine cavity.


- The integrated support cannula aids in placing the transfer catheter.



## Guardia™ Access Nano Embryo Transfer Catheter



Precurved  
with  
bulb tip

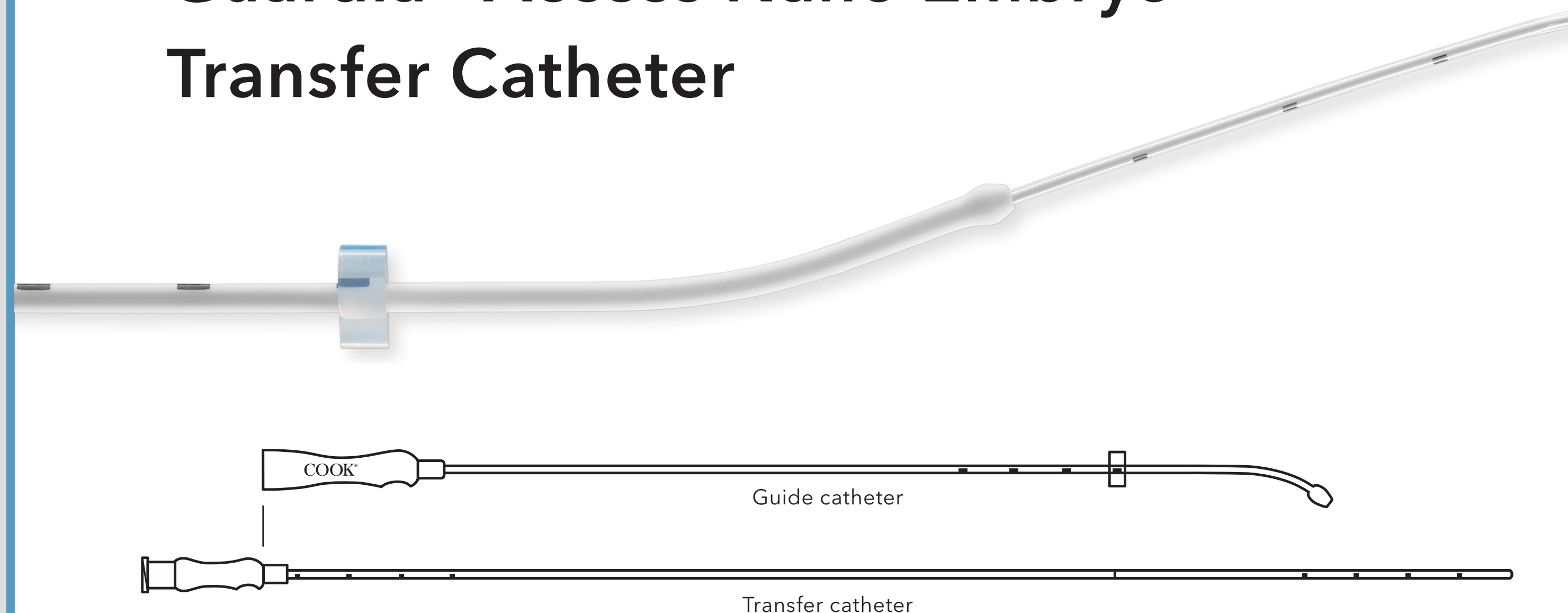


- The outer diameter of the guide catheter's bulb tip is reduced for smaller anatomy.
- The transfer catheter has an internal support cannula.
- The guide catheter is curved to facilitate insertion.
- The guide catheter has a rounded bulb tip to ease passage through the cervix.





## Echosight® Jansen-Anderson Intratubal Transfer Set

Specialty

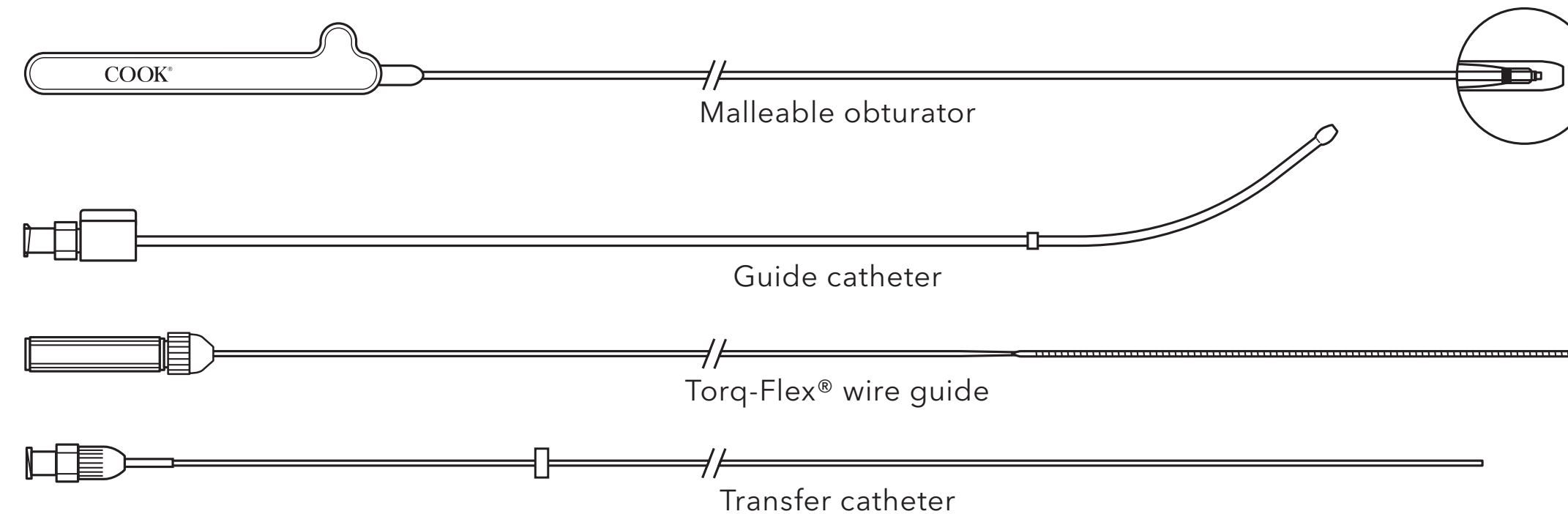


Used to inject either sperm, gametes, or embryos into the uterine ostium of the fallopian tube via ultrasound guidance.

- The catheters' polyethylene Echosight material enhances visualization under ultrasound.
- The transfer catheter has a female Luer-lock adapter.



## Malleable: Guardia™ Obturator

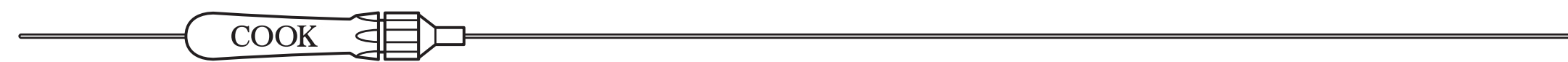

Obturator


- Works with all Cook Medical embryo transfer catheters and many other embryo transfer catheters.
- The hub locks in place.



## Stiff: Guardia™ Obturator



Obturator

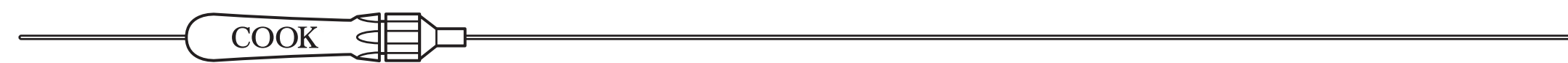

- The adjustable length works with a variety of Cook Medical and other embryo transfer catheters.
- The hub locks in place.



## Guardia™ Access Malleable Obturator



Obturator


- Made to fit the standard Guardia Access Embryo Transfer Catheter (K-JETS-7019) and fits the AccessET.
- Malleable for maneuverability.





## Soft-Pass™ Malleable Obturator

Obturator

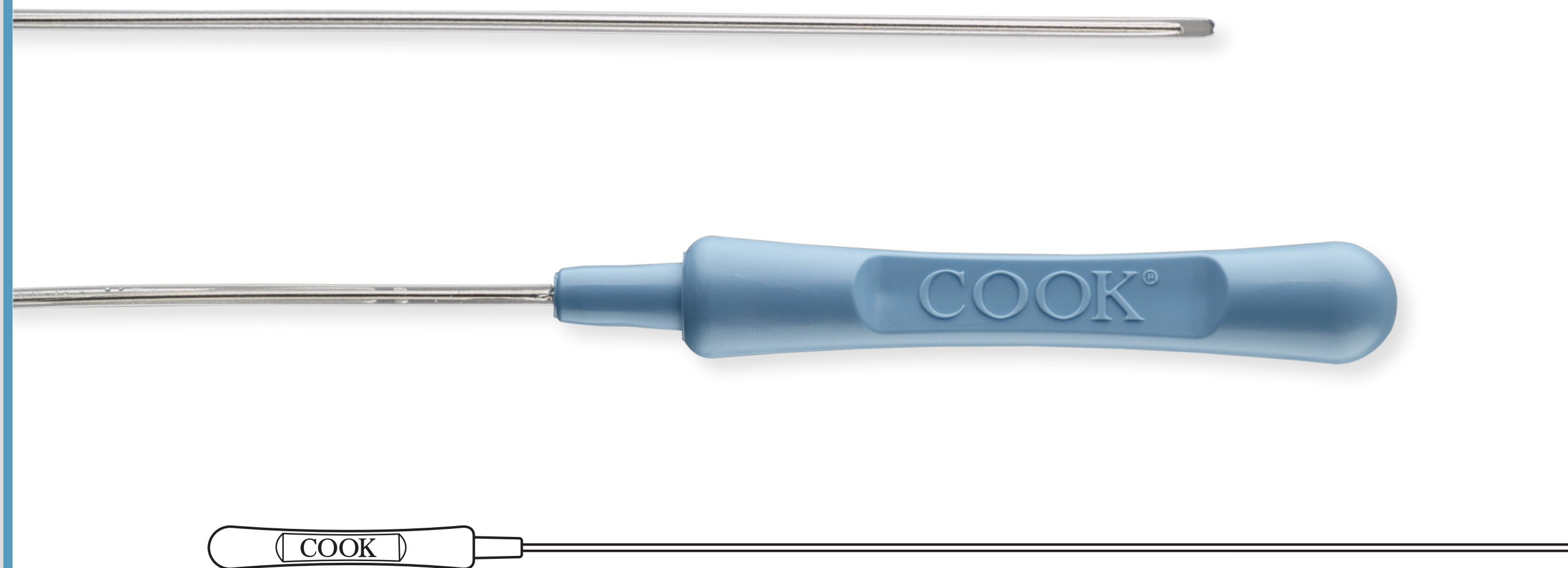

Used to supplement and assist uterine access of an approved embryo transfer device for placement of in vitro fertilized (IVF) embryos into the uterine cavity.

- Designed to fit Soft-Pass catheters.