

# Maintain catheter placement and reduce leakage during sonohysterography.

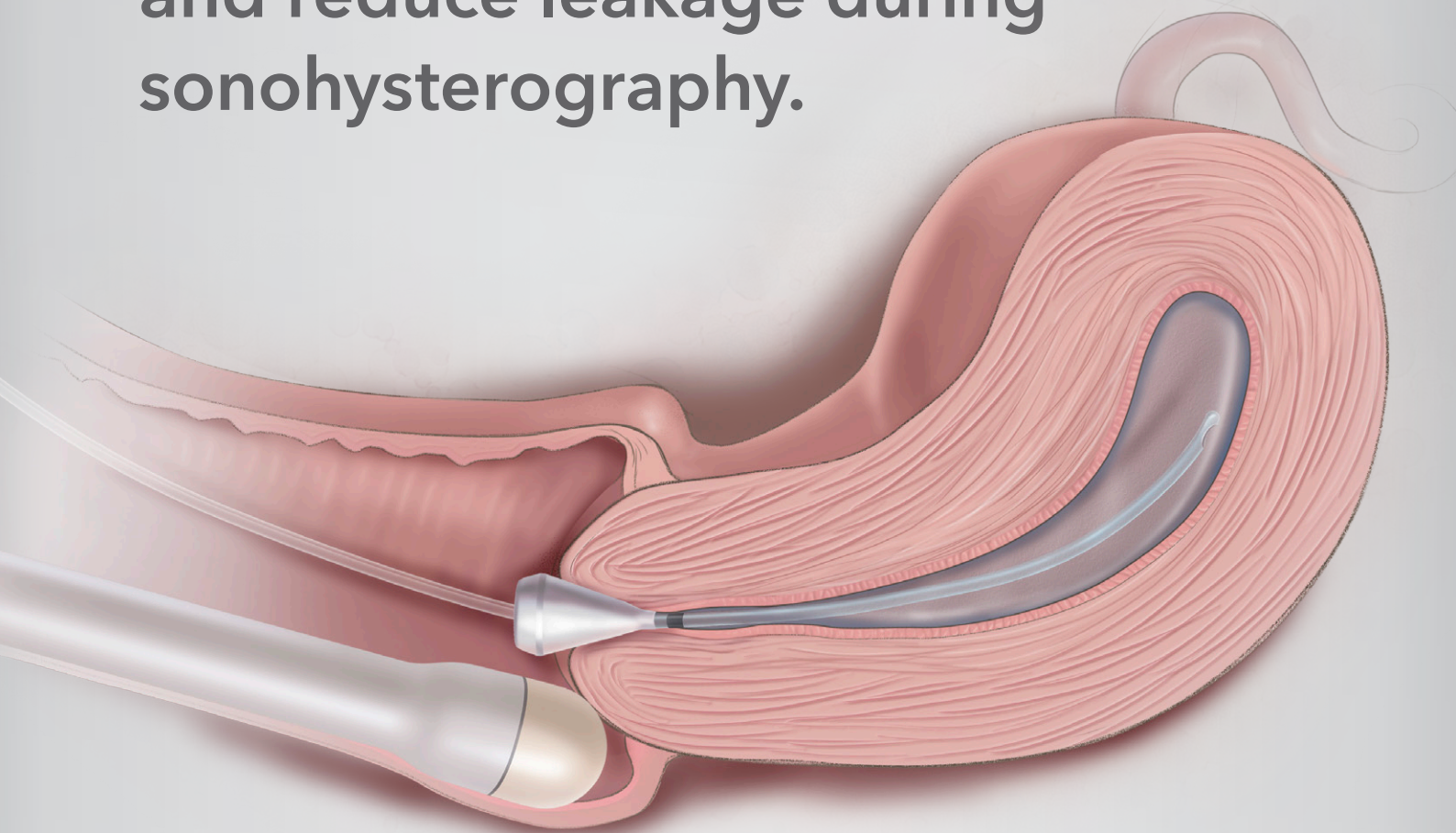


Illustration by Lisa Clark

## Goldstein

SONOHYSTEROGRAPHY CATHETER

### Versatility

Formable material maintains the catheter shape for ease of placement.

### Efficiency

The rounded tip eases passage of the catheter through the cervix. A large sideport helps with the instillation of saline.

### Control

The signature acorn-shaped positioner occludes the cervix to reduce the loss of instilled saline during the procedure.



[www.cookmedical.com](http://www.cookmedical.com)



## Goldstein

### SONOHYSTEROGRAPHY CATHETER

The Goldstein Sonohysterography Catheter is used to access the uterine cavity for sonohysterography. It is available in either PTFE or polyethylene material.

- The adjustable acorn-shaped positioner maintains catheter placement while reducing the leakage of instilled saline at the external cervical os.
- An ink band provides a reference mark to indicate that proper catheter positioning has been achieved.
- The round, closed catheter tip eases the introduction of the catheter into the uterus.
- The large oval sideport facilitates the instillation of saline solution.
- The formable polytetrafluoroethylene (PTFE) versions maintain their shape to aid in placement.

Order Number	Reference Part Number	Fr	Length cm
PTFE			
G18664	J-GSHC-522620	5.2	26
G18665	J-GSHC-523620	5.2	36
Polyethylene			
G17149	J-GSHC-532600	5.3	26
Polyethylene with Outer Support Sheath			
G19012	J-GSHC-542600-SV	5.4	26

*Sold in boxes of 10.*

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

### Customer Service

**EMEA:** EDI - [www.cookmedical.com/edi.do](http://www.cookmedical.com/edi.do)  
 Distributors: +353 61239240, [ssc.distributors@cookmedical.com](mailto:ssc.distributors@cookmedical.com)  
 Austria: +43 179567121, [oe.orders@cookmedical.com](mailto:oe.orders@cookmedical.com)  
 Belgium: +32 27001633, [be.orders@cookmedical.com](mailto:be.orders@cookmedical.com)  
 Denmark: +45 38487607, [da.orders@cookmedical.com](mailto:da.orders@cookmedical.com)  
 Finland: +358 972519996, [fi.orders@cookmedical.com](mailto:fi.orders@cookmedical.com)  
 France: +33 171230269, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
 Germany: +49 6950072804, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
 Hungary: +36 17779199, [hu.orders@cookmedical.com](mailto:hu.orders@cookmedical.com)  
 Ireland: +353 61239252, [ie.orders@cookmedical.com](mailto:ie.orders@cookmedical.com)  
 Italy: +39 0269682853, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
 Netherlands: +31 202013367, [nl.orders@cookmedical.com](mailto:nl.orders@cookmedical.com)  
 Poland: +48 223060159, [pl.orders@cookmedical.com](mailto:pl.orders@cookmedical.com)  
 Spain: +34 912702691, [es.orders@cookmedical.com](mailto:es.orders@cookmedical.com)  
 Sweden: +46 858769468, [se.orders@cookmedical.com](mailto:se.orders@cookmedical.com)  
 Switzerland - French: +41 448009609, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
 Switzerland - Italian: +41 448009609, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
 Switzerland - German: +41 448009609, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
 United Kingdom: +44 2073654183, [uk.orders@cookmedical.com](mailto:uk.orders@cookmedical.com)

**Americas:** EDI - [www.cookmedical.com/edi.do](http://www.cookmedical.com/edi.do)  
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335  
 E-mail: [orders@cookmedical.com](mailto:orders@cookmedical.com)

**Australia:**  
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283  
 E-mail: [cau.custserv@cookmedical.com](mailto:cau.custserv@cookmedical.com)



[www.cookmedical.com](http://www.cookmedical.com)