




# Extract kidney stones with the nitinol advantage.

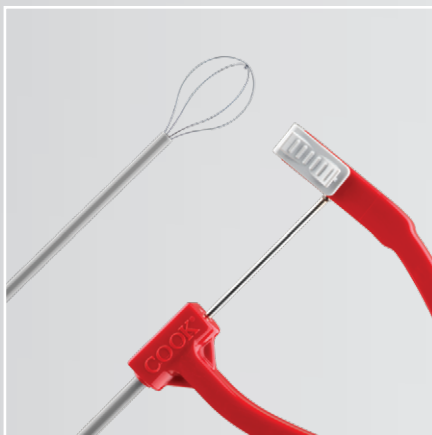

**Perc NCircle®**  
Nitinol Tipless Stone Extractor

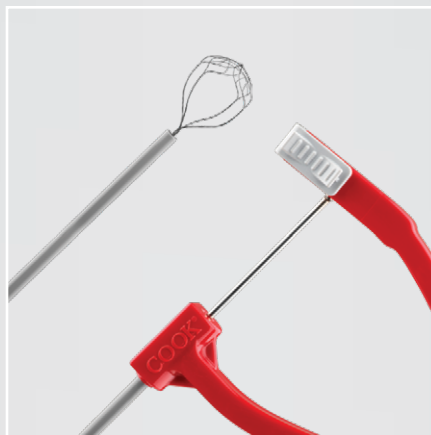

**Perc NCompass®**  
Nitinol Stone Extractor

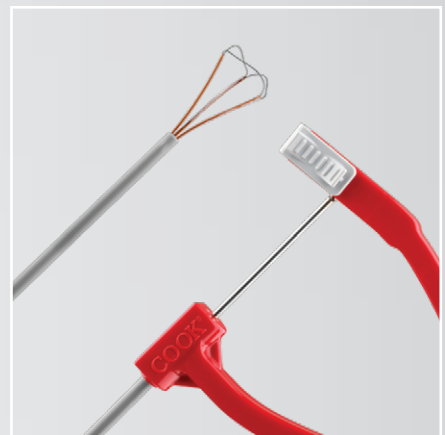

**Perc NGage®**  
Nitinol Stone Extractor

Cook Medical offers the industry's most comprehensive line of nitinol stone extractors.

Three of our extractors are specially designed to remove kidney calculi: Perc NCircle, Perc NCompass, and Perc NGage. With three extractor options, you can easily select the right device for the procedure.

With our percutaneous extractors, you can quickly and confidently remove large stones and fragments that might be too big for traditional forceps and graspers.

Each of these devices is constructed of nitinol wire and has a stainless steel cannula.

## Features

### Nitinol wire construction


Minimizes trauma and allows you to remove larger kidney stones and fragments.

### Stainless steel cannula

Provides strength and reliability during kidney stone extraction.

### FEP coating


Reduces friction, protects the working channel of the nephroscope, and minimizes the risk of trauma to tissue.

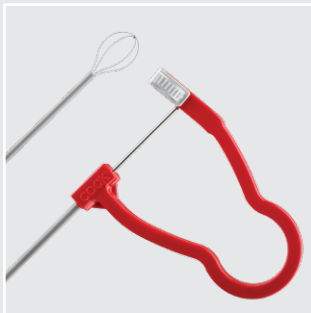

## Perc NCircle<sup>®</sup>


NITINOL TIPLESS STONE EXTRACTOR

The tipless design of the Perc NCircle allows you to position it directly against the mucosal lining with minimal trauma. Reducing trauma helps to reduce bleeding and maintain a clear working field. The small, flexible nitinol wires have an advantage over traditional PCNL forceps or graspers, because the wires allow you to remove large stones or fragments.

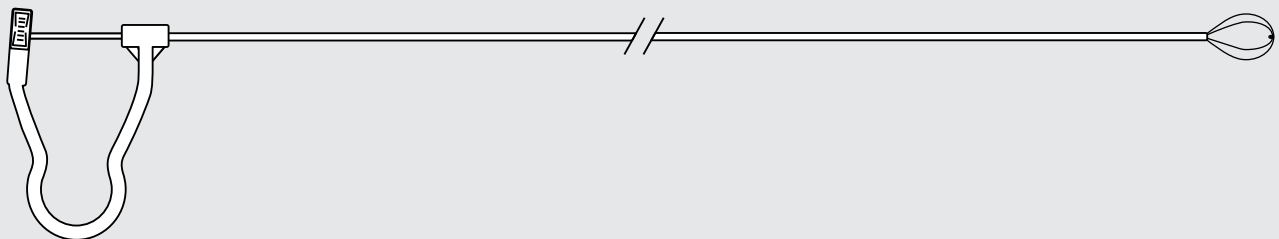

Order Number	Reference Part Number	Fr	Length cm	Basket Diameter cm
G32862	NTRSE-100038	10.0	38	2


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

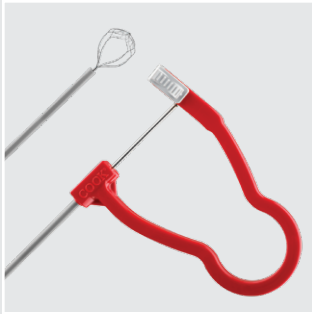

## Perc NCompass®


NITINOL STONE EXTRACTOR

The multiwire geometry of the Perc NCompass creates two basket designs: a standard 4-wire design at full deployment and an intricate 16-wire design during closure. At full deployment, the distal portion of the basket produces a netting effect to help prevent stone migration. The tiptless basket reduces the risk of trauma to surrounding tissue.

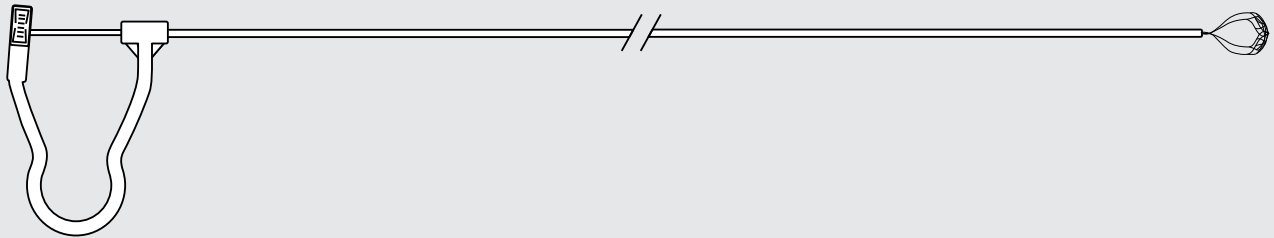

Order Number	Reference Part Number	Fr	Length cm	Basket Diameter cm
G34931	NCTR-100038	10.0	38	1.5


Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

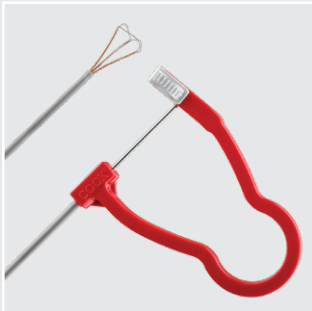

## Perc NGage®

NITINOL STONE EXTRACTOR

The Perc NGage is like two extractors in one. Use it to catch, reposition, and release stones as you would with graspers, or use it as a basket to secure stones for extraction. The resilient nitinol wire allows you to completely and consistently open the extractor to remove large stone and fragments and, when necessary, to disengage stones.

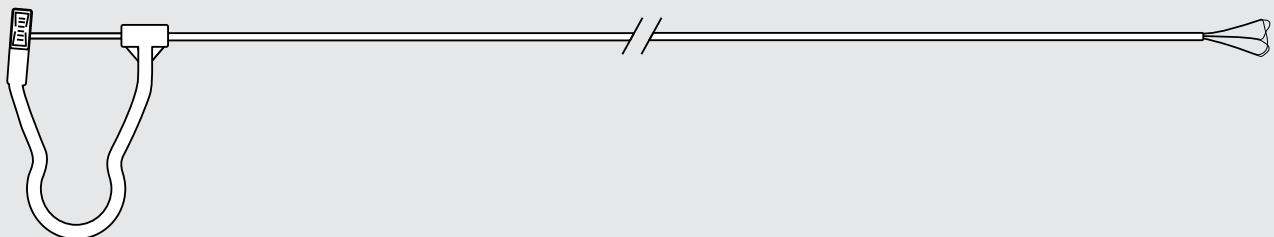

Order Number	Reference Part Number	Fr	Length cm	Basket Diameter cm
G56401	NGRE-100038	10.0	38	1.5

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

### Customer Service - Urology

**EMEA:** EDI - [www.cookmedical.com/edi.do](http://www.cookmedical.com/edi.do)  
Distributors: +353 61239240, [ssc.distributors@cookmedical.com](mailto:ssc.distributors@cookmedical.com)  
Austria: +43 179567121, [oe.orders@cookmedical.com](mailto:oe.orders@cookmedical.com)  
Belgium: +32 27001633, [be.orders@cookmedical.com](mailto:be.orders@cookmedical.com)  
Denmark: +45 38487607, [da.orders@cookmedical.com](mailto:da.orders@cookmedical.com)  
Finland: +358 972519996, [fi.orders@cookmedical.com](mailto:fi.orders@cookmedical.com)  
France: +33 171230269, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
Germany: +49 6950072804, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
Hungary: +36 17779199, [hu.orders@cookmedical.com](mailto:hu.orders@cookmedical.com)  
Ireland: +353 61239252, [ie.orders@cookmedical.com](mailto:ie.orders@cookmedical.com)  
Italy: +39 0269682853, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
Netherlands: +31 202013367, [nl.orders@cookmedical.com](mailto:nl.orders@cookmedical.com)  
Norway: +47 23162968, [no.orders@cookmedical.com](mailto:no.orders@cookmedical.com)  
Poland: +48 223060159, [pl.orders@cookmedical.com](mailto:pl.orders@cookmedical.com)  
Spain: +34 912702691, [es.orders@cookmedical.com](mailto:es.orders@cookmedical.com)  
Sweden: +46 858769468, [se.orders@cookmedical.com](mailto:se.orders@cookmedical.com)  
Switzerland - French: +41 448009609, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
Switzerland - Italian: +41 448009609, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
Switzerland - German: +41 448009609, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
United Kingdom: +44 2073654183, [uk.orders@cookmedical.com](mailto:uk.orders@cookmedical.com)  
**[www.cookmedical.com](http://www.cookmedical.com)**

**Americas:** EDI - [www.cookmedical.com/edi.do](http://www.cookmedical.com/edi.do)  
Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335  
E-mail: [orders@cookmedical.com](mailto:orders@cookmedical.com)

**Australia:**  
Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283  
E-mail: [cau.custserv@cookmedical.com](mailto:cau.custserv@cookmedical.com)

