Full stream ahead.








Upgrade to the UroStream Hydrophilic Wire Guide and see how we've enhanced your access as you look to navigate through impacted stones, ureteral obstructions, and tortuous anatomy.


Designed with a unique radiopaque polymer (tungsten-impregnated) jacket to aid visualization under fluoroscopy, the UroStream wire guide was shown to have superior radiopacity to other wire guides available on the market.*

It features a reinforced, kink-resistant nitinol core, and it has a flexible distal tip to reduce the risk of undue trauma.

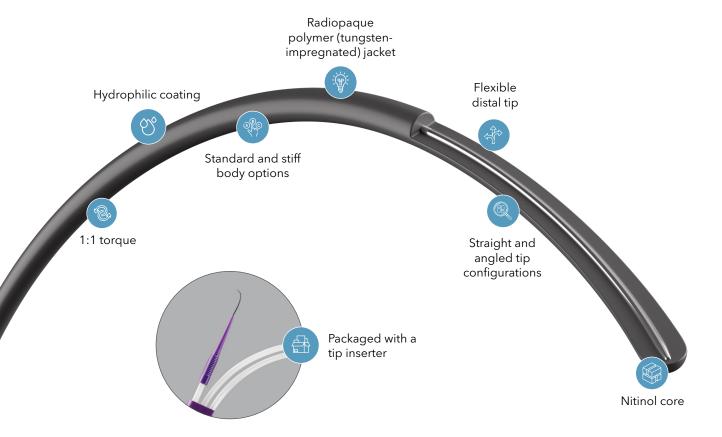



Get started with a UroStream wire guide and explore the benefits of our complete product portfolio to streamline your procedures.

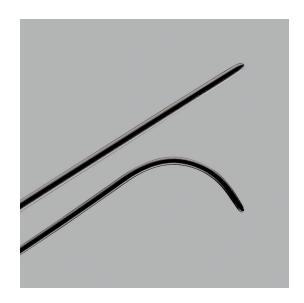

UroStream[™] Hydrophilic Wire Guide

Used to gain access, to establish a tract, and to assist in the placement, replacement, and exchange of medical devices during urological procedures.

Order Number	Referene Part Number	Diameter inch	Length cm	Taper Length cm	Flexible Tip Length cm
Standard Sha	ft with a Straight Tip: Nitin	ol with a Tungsten-	Impregnated Po	lyurethane Jacket	
G59149	UROS-025150	0.025	150	14	3
G59151	UROS-035150	0.035	150	14	3
G59155	UROS-038150	0.038	150	14	3
Standard Sha	ft with an Angled Tip: Nitin	ol with a Tungsten	-Impregnated P	olyurethane Jacket	
G59150	UROS-025150-A	0.025	150	14	3
G59152	UROS-035150-A	0.035	150	14	3
G59158	UROS-038150-A	0.038	150	14	3
Stiff Shaft wit	h a Straight Tip: Nitinol wit	h a Tungsten-Impre	egnated Polyure	thane Jacket	
G59153	UROS-035150-S	0.035	150	20	3
G59156	UROS-038150-S	0.038	150	20	3
Stiff Shaft wit	h an Angled Tip: Nitinol wi	th a Tungsten-Impr	egnated Polyur	ethane Jacket	
G59154	UROS-035150-AS	0.035	150	20	3
G59159	UROS-038150-AS	0.038	150	20	3

^{*}Radiopacity testing of Cook UroStream wire guide and currently marketed wire guides. TS210234-R (2021).

Customer Support

EU Website: cookmedical.eu EDI: cookmedical.eu/edi

Distributors: +353 61239240, ssc.distributors@cookmedical.com

Austria: +43 179567121, oe.orders@cookmedical.com Belgium: +32 27001702, be.orders@cookmedical.com Denmark: +45 38487607, da.orders@cookmedical.com Finland: +358 972519996, fi.orders@cookmedical.com France: +33 171230269, fr.orders@cookmedical.com Germany: +49 6950072804, de.orders@cookmedical.com Hungary: +36 17779199, hu.orders@cookmedical.com Ireland: +353 61239252, ie.orders@cookmedical.com Italy: +39 0269682853, it.orders@cookmedical.com Netherlands: +31 202013367, nl.orders@cookmedical.com

Poland: +48 223060159, pl.orders@cookmedical.com Spain: +34 912702691, es.orders@cookmedical.com

Sweden: +46 858769468, se.orders@cookmedical.com Switzerland - French: +41 448009609, fr.orders@cookmedical.com Switzerland - Italian: +41 448009609, it.orders@cookmedical.com Switzerland - German: +41 448009609, de.orders@cookmedical.com

United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com EDI: cookmedical.com/edi.do Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335 E-mail: customersupport@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283 E-mail: cau.custserv@cookmedical.com

