

Simplify  
percutaneous  
dilation.



Ultraxx™  
NEPHROSTOMY BALLOON

# Ultraxx™

## NEPHROSTOMY BALLOON

The Ultraxx Nephrostomy Balloon is intended to dilate the musculofascial tract, renal capsule, and parenchyma to establish and maintain a percutaneous tract. The balloon has a rated burst pressure of 20 atm, and the functional dilating portion of the balloon is 15 cm long. The combination of high pressure and balloon length helps a physician establish a nephrostomy tract.

The Ultraxx is available with either a standard PTFE or a clear vinyl access sheath and can be purchased as a set with the Cook Inflation Device.



## Features

---

Balloon port grip	Provides added control and aids handling during the inflation process
-------------------	-----------------------------------------------------------------------

---

Sheath lock	Holds the sheath in place and out of the way during dilation
-------------	--------------------------------------------------------------

---

Inflation-deflation port	Provides uniform pressure during inflation; when the balloon is deflated, the port provides complete fluid evacuation for a less complicated removal process.
--------------------------	---------------------------------------------------------------------------------------------------------------------------------------------------------------

---

Platinum band	Provides radiopacity to ensure precise balloon dilation positioning
---------------	---------------------------------------------------------------------

---


Tapered tip	Eases insertion and allows placement over a 0.038 inch diameter wire guide
-------------	----------------------------------------------------------------------------

---

Clear sheath*	Enables visualization within the sheath
---------------	-----------------------------------------

---

\* The clear sheath is only for specific part numbers.

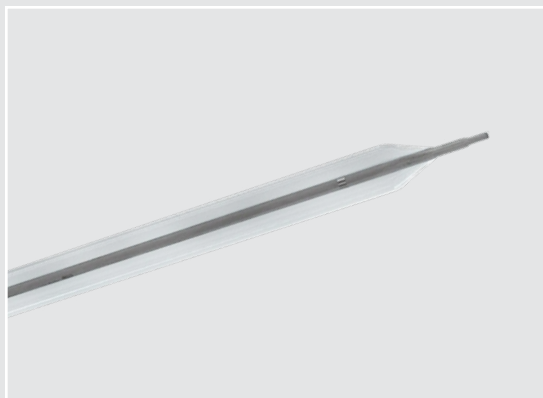

## Ultraxx™

NEPHROSTOMY BALLOON

Ultraxx features a high-pressure balloon for establishing and maintaining a percutaneous tract.

Order Number	Reference Part Number	Inflated Balloon Diameter mm	Balloon Length cm
<b>Set</b> With Clear Vinyl Sheath			
G47646	UNBS-6-15-CS-01	6	15
G47648	UNBS-8-15-CS-01	8	15
G47658	UNBS-10-15-CS-01	10	15
With PTFE Sheath			
G47645	UNBS-6-15-01	6	15
G47647	UNBS-8-15-01	8	15
G47649	UNBS-10-15-01	10	15

Order Number	Reference Part Number	Inflated Balloon Diameter mm	Balloon Length cm
<b>Balloon Catheter Only</b> With Clear Vinyl Sheath			
G57618	UNB-6-15-CS-01	6	15
G57619	UNB-8-15-CS-01	8	15
G57620	UNB-10-15-CS-01	10	15
With PTFE Sheath			
G57615	UNB-6-15-01	6	15
G57616	UNB-8-15-01	8	15
G57617	UNB-10-15-01	10	15

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Support & Distribution for details.

### Customer Service

EU Website: [cookmedical.eu](http://cookmedical.eu)  
 EDI: [cookmedical.eu/edi](http://cookmedical.eu/edi)  
 Distributors: +353 61239240, [ssc.distributors@cookmedical.com](mailto:ssc.distributors@cookmedical.com)  
 Austria: +43 179567121, [oe.orders@cookmedical.com](mailto:oe.orders@cookmedical.com)  
 Belgium: +32 27001702, [be.orders@cookmedical.com](mailto:be.orders@cookmedical.com)  
 Denmark: +45 38487607, [da.orders@cookmedical.com](mailto:da.orders@cookmedical.com)  
 Finland: +358 972519996, [fi.orders@cookmedical.com](mailto:fi.orders@cookmedical.com)  
 France: +33 171230269, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
 Germany: +49 6950072804, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
 Hungary: +36 17779199, [hu.orders@cookmedical.com](mailto:hu.orders@cookmedical.com)  
 Iceland: +354 8007615, [is.orders@cookmedical.com](mailto:is.orders@cookmedical.com)  
 Ireland: +353 61239252, [ie.orders@cookmedical.com](mailto:ie.orders@cookmedical.com)  
 Italy: +39 0269682853, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
 Netherlands: +31 202013367, [nl.orders@cookmedical.com](mailto:nl.orders@cookmedical.com)  
 Poland: +48 223060159, [pl.orders@cookmedical.com](mailto:pl.orders@cookmedical.com)  
 Spain: +34 912702691, [es.orders@cookmedical.com](mailto:es.orders@cookmedical.com)  
 Sweden: +46 858769468, [se.orders@cookmedical.com](mailto:se.orders@cookmedical.com)  
 Switzerland - French: +41 448009609, [fr.orders@cookmedical.com](mailto:fr.orders@cookmedical.com)  
 Switzerland - Italian: +41 448009609, [it.orders@cookmedical.com](mailto:it.orders@cookmedical.com)  
 Switzerland - German: +41 448009609, [de.orders@cookmedical.com](mailto:de.orders@cookmedical.com)  
 United Kingdom: +44 2073654183, [uk.orders@cookmedical.com](mailto:uk.orders@cookmedical.com)

USA Website: [cookmedical.com](http://cookmedical.com)  
 EDI: [cookmedical.com/edi.do](http://cookmedical.com/edi.do)  
 Americas:  
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335  
 E-mail: [customersupport@cookmedical.com](mailto:customersupport@cookmedical.com)  
 Australia:  
 Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283  
 E-mail: [cau.custserv@cookmedical.com](mailto:cau.custserv@cookmedical.com)



CC, URO-8.5X11