


Full stream ahead




UroStream™
HYDROPHILIC WIRE GUIDE

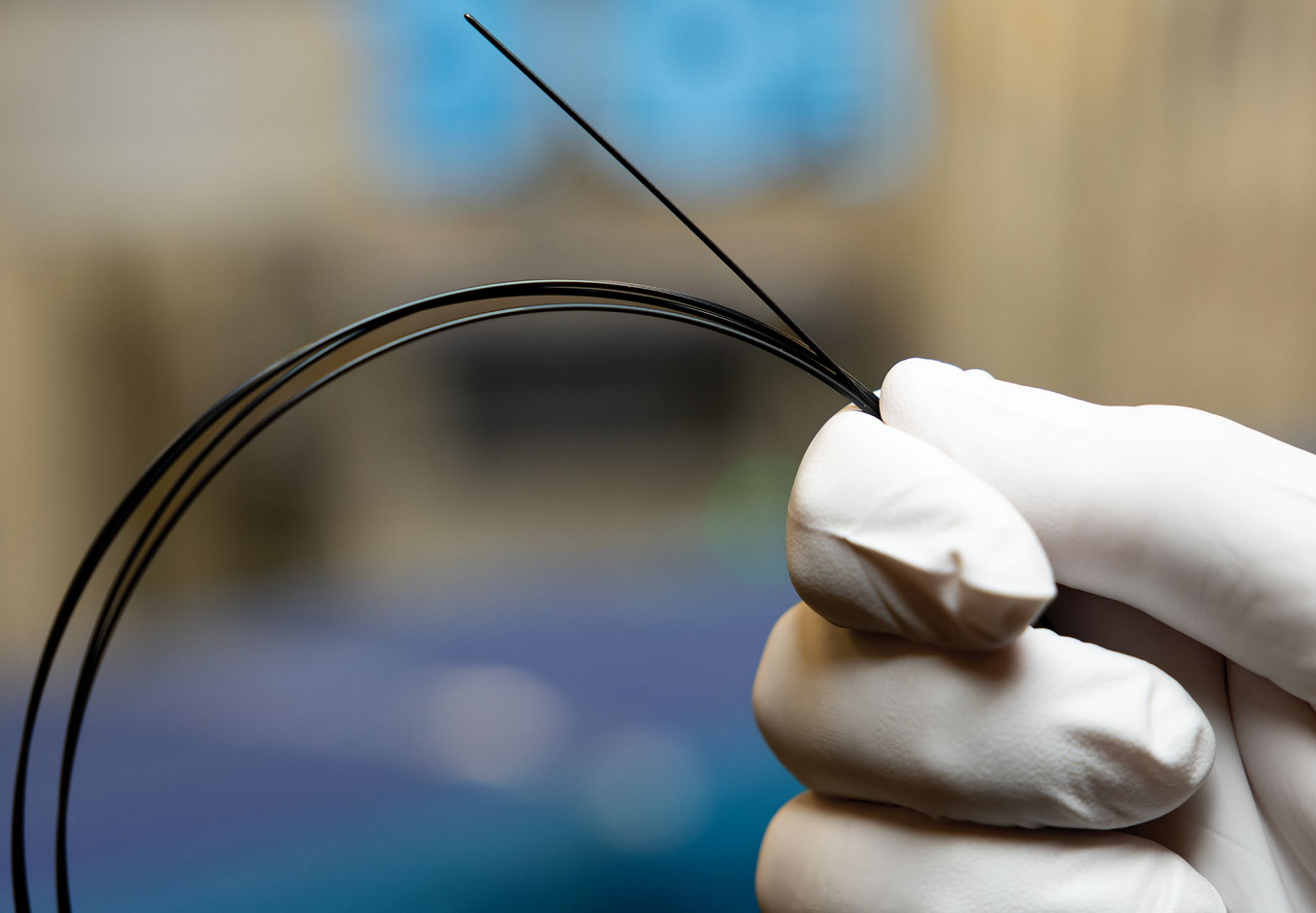

Upgrade to the UroStream Hydrophilic Wire Guide and see how we've enhanced your navigation through impacted stones, ureteral obstructions, and tortuous anatomy.

Designed with a unique radiopaque polymer (tungsten-impregnated) jacket to aid visualisation under fluoroscopy, the UroStream wire guide was shown to have superior radiopacity to other wire guides available on the market.*

It features a reinforced, kink-resistant nitinol core, and it has a flexible distal tip to reduce the risk of undue trauma.

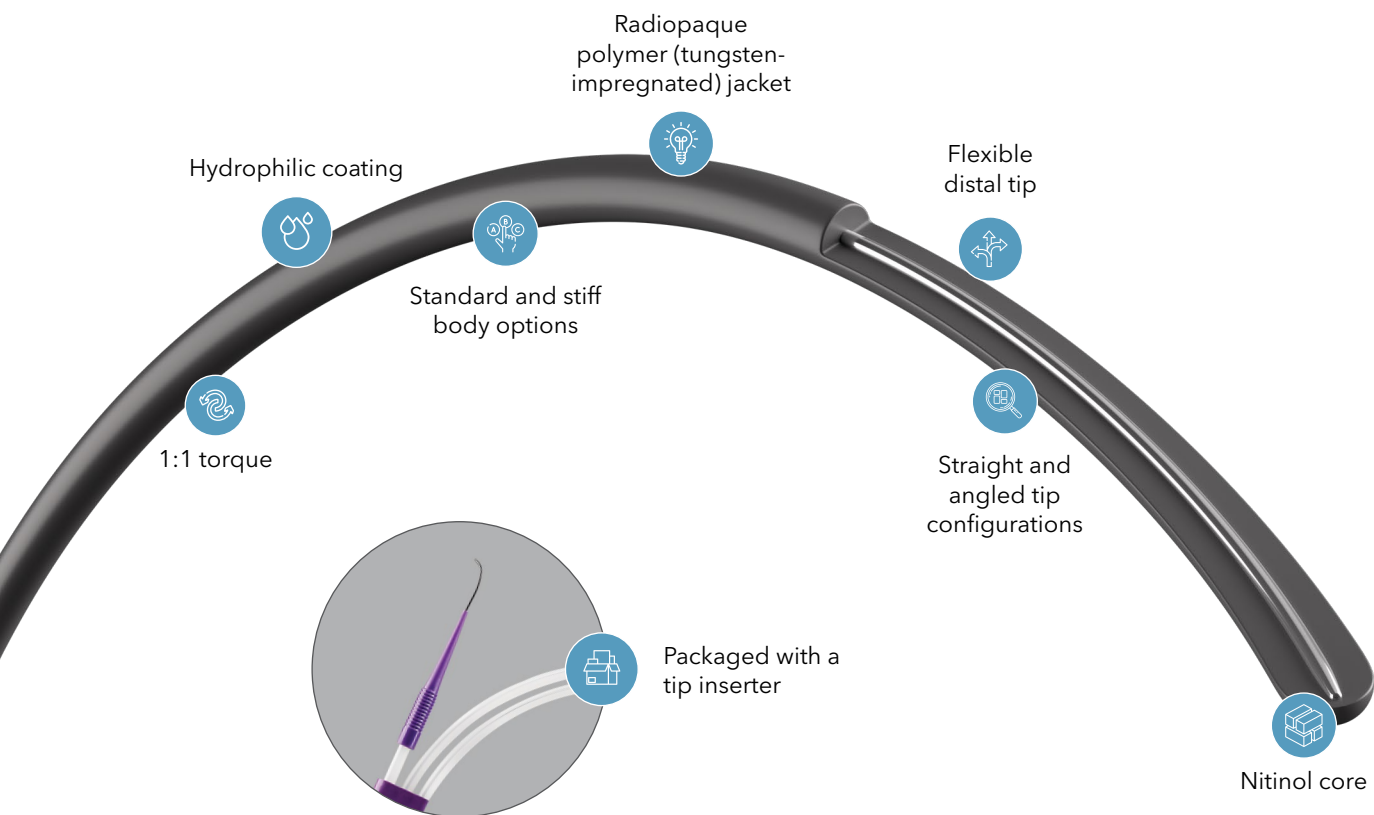


*Radiopacity testing of Cook UroStream wire guide and currently marketed wire guides. TS210234-R (2021).



Crafted by Cook—from start to finish!

Get started with a UroStream wire guide and explore the benefits of our complete product portfolio to streamline your procedures.

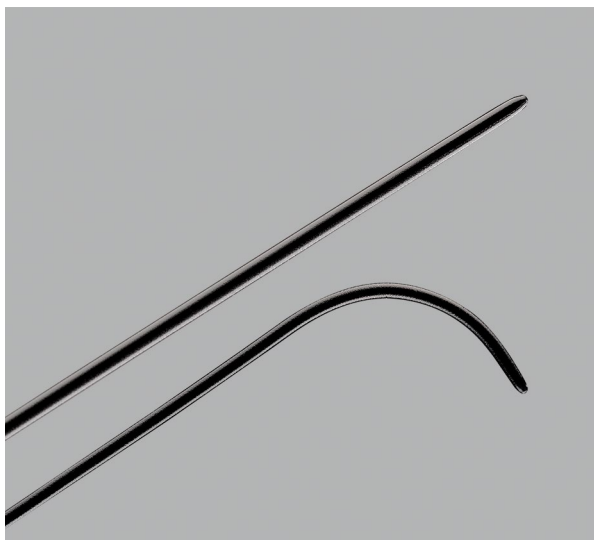

UroStream™ Hydrophilic Wire Guide

Used to gain access, to establish a tract, and to assist in the placement, replacement, and exchange of medical devices during urological procedures.

Order Number	Referene Part Number	Diameter inch	Length cm	Taper Length cm	Flexible Tip Length cm
Standard Shaft with a Straight Tip: Nitinol with a Tungsten-Impregnated Polyurethane Jacket					
G59149	UROS-025150	0.025	150	14	3
G59151	UROS-035150	0.035	150	14	3
G59155	UROS-038150	0.038	150	14	3
Standard Shaft with an Angled Tip: Nitinol with a Tungsten-Impregnated Polyurethane Jacket					
G59150	UROS-025150-A	0.025	150	14	3
G59152	UROS-035150-A	0.035	150	14	3
G59158	UROS-038150-A	0.038	150	14	3
Stiff Shaft with a Straight Tip: Nitinol with a Tungsten-Impregnated Polyurethane Jacket					
G59153	UROS-035150-S	0.035	150	20	3
G59156	UROS-038150-S	0.038	150	20	3
Stiff Shaft with an Angled Tip: Nitinol with a Tungsten-Impregnated Polyurethane Jacket					
G59154	UROS-035150-AS	0.035	150	20	3
G59159	UROS-038150-AS	0.038	150	20	3

Customer Support

EU Website: cookmedical.eu
EDI: cookmedical.eu/edi
Distributors: +353 61239240, ssc.distributors@cookmedical.com

Austria: +43 179567121, oe.orders@cookmedical.com
Belgium: +32 27001702, be.orders@cookmedical.com
Denmark: +45 38487607, da.orders@cookmedical.com
Finland: +358 972519996, fi.orders@cookmedical.com
France: +33 171230269, fr.orders@cookmedical.com
Germany: +49 6950072804, de.orders@cookmedical.com
Hungary: +36 17779199, hu.orders@cookmedical.com
Ireland: +353 61239252, ie.orders@cookmedical.com
Italy: +39 0269682853, it.orders@cookmedical.com
Netherlands: +31 202013367, nl.orders@cookmedical.com
Poland: +48 223060159, pl.orders@cookmedical.com
Spain: +34 912702691, es.orders@cookmedical.com
Sweden: +46 858769468, se.orders@cookmedical.com
Switzerland - French: +41 448009609, fr.orders@cookmedical.com
Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
Switzerland - German: +41 448009609, de.orders@cookmedical.com
United Kingdom: +44 2073654183, uk.orders@cookmedical.com

USA Website: cookmedical.com

EDI: cookmedical.com/edi.do

Americas:

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: customersupport@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283

E-mail: cau.custserv@cookmedical.com



CC, URO-A4